



User Manual



www.Evertoneprosage.com

Thank you for your purchase of Prosage Thermo. Please read the user manual before use.

The **Prosage Thermo** is the only professional Hot and Cold massage, percussive therapy gun with warm up technology for everyday pains.

- Prepares your muscles for the deep tissue massage
- Helps prevent future injuries
- Helps Improve range of motion
- Helps activate your muscles to perform at peak levels.

Safety Instruction

- For adult use only.
- Use only on a dry surface of the body by lightly pressing and moving across the body for approximately 60 seconds per region.
- Use the **Prosage Thermo** only on the body's soft muscle tissue as desired without producing pain or discomfort. Do not use on the head or any hard, bony area of the body.
- Use only the attachments provided.
- Bruising can occur regardless of control setting or pressure applied. Check treated areas frequently and immediately stop at first sign of pain or discomfort.
- Keep fingers, hair, or any other body part away from the shaft and backside of the applicator as pinching may occur.
- Do not place any objects into the ventilation holes of the **Prosage Thermo**.
- Do not immerse in water or allow water to enter the ventilation holes of the **Prosage Thermo**.
- Do not drop or misuse the **Prosage Thermo**.
- Only recharge with supplied 9VDC **Prosage Thermo** Charger
- Never leave the **Prosage Thermo** operating or charging unattended.

Product Features

This Set Includes:

- 1 x Massager Device
- 8 x Interchangeable Attachments
- 1 x Charging Base/ Holding Stand
- 1 x DC Adaptor

Attachments:


1. Silicone Ball Head
2. Silicone Cone Head
3. Silicone Accupressure Head
4. Heat Attachment
5. Cold Gel Attachment
6. Silicone Fork Head
7. Silicone Fingers Head
8. Scalp Massage Head

MicrocurrentTek pads
Extension handle



Charge the device

- Remove the adaptor cover at the bottom of the charging base. (See Fig. 1)
- Place the 9VDC adaptor into compartment (See. Fig 2) and secure the cover.
- Then the adaptor jack will be fixed firmly on the charging base (See Fig. 3)
- Plug the AC adaptor cable into a wall outlet and then place the device on the charging base (See Fig. 4)
- The charging indicator light on the device (above the control button) **will blink Red**. Upon full charge the indicator **light will turn Blue**. (*Please Note: The device cannot be used whilst charging.)
- Charge for 3-4 Hours for full charge. This will charge the device for use of 70 mins.

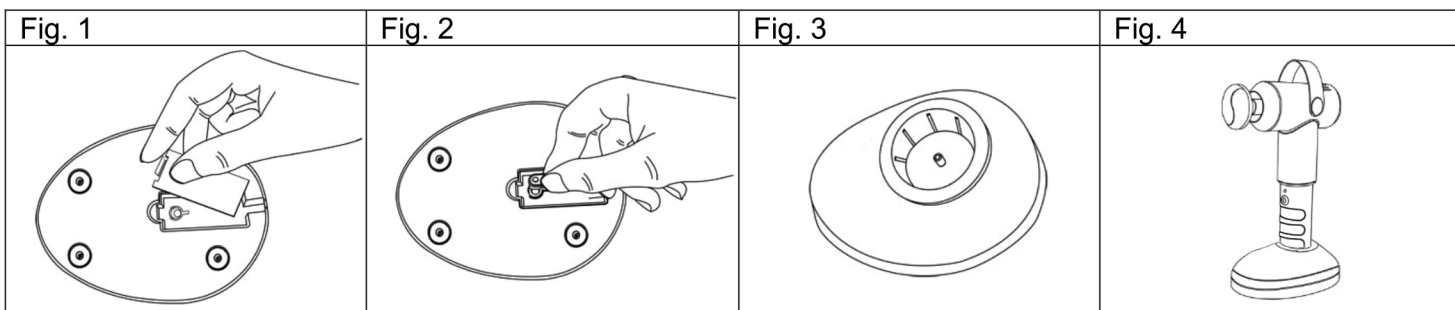


Operating Instructions

- Select the appropriate interchangeable attachment for your needs, then screw the selected attachment (clockwise) the into the massager head.
 - *Please note: The Cold Gel attachment needs to be placed in the freezer for 1-2 hours before use.
- Long press the motor button on the control board for 3 seconds to start the motor; Default motor rotation (1), blue light panel on; Press the motor button on the main control board again to make the motor faster (in second gear); Press again, the motor becomes faster (in third gear); Press the cycle again to (1) in turn; Press the button for another 3 seconds to stop the motor and put out the blue light.
- The blue light on the shaft will indicate the device is ON and will auto start at LO level. Press the power button again to select MID level, and once more for HI level.
 - *Please note: The Theralite function is not intended for use with the Hot and Cold attachment heads. The Theralite will be blocked during use with the Hot & Cold attachments.
- Keep the heating switch on the control panel for 3 seconds, the PTC heating plate is powered on and the red plate is on, then keep the heating switch on for 3 seconds and the red light is off.
- EMS functions: massage, hammering, kneading, firming, acupuncture, intensity 1-5; Specific operation: long button on the control panel open EMS function mode button 3 seconds, the default massage, press the mode button again to switch to the beat, to press the mode button again to switch to the kneading, press the mode button again, switching to compact file to press the mode button again to switch to the acupuncture and moxibustion, massage, press the mode button once to switch in circulation, longer toggle mode button on the control panel for 3 seconds to shut down the EMS functions; When each gear is on, the mode intensity is 1. The "+" and "-" buttons on the control panel can be pressed to adjust the intensity.



Press the MOTOR button for 3 seconds to Turn ON

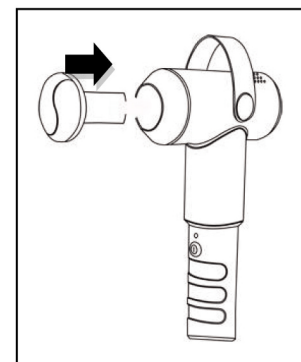

- Charging status: charging with a blue light on and off to prompt; Full of blue lights always on; All functions of pressing the switch in charging state do not work;
- Blocking protection tip: When the red light flashes continuously for 5 times in a row and stops, the motor button can be started normally.
- Low voltage alarm, the red light flashes for 5 times in a row, and all switches cannot work when pressed;
- Massage the desired body part at the desired speed, applying pressure as needed.
- A blocking protection is set up, when you continue to press your body with device too hard, the device will blink and automatically shuts off to protect internal electronic components from burning. Press START to resume operating the device.
- Use the silicone strap to better assist your grip and angle whilst in use.
- The massager is preset with a 10 minute AUTO-OFF timer.
- When the power/intensity begins to lessen (slow down), please place the device on the charging base to recharge.

Maintenance and Caring

- For Maintenance use a slightly damp towel to clean the device and dry with a soft dry cloth
- Do not use the product for over 30 minutes continuously. Doing so may activate the over heating protection switch and stop operation automatically. If this happens, please wait at least 30 minutes in order to let the device cool down before re using again.

Technical Specification

Battery Capacity	2200 mAH
Power by	AC adaptor input 100-240V/50-60Hz
Output Voltage	9V DC 1A
Max. Current	6000 mA
Power Consumption	40 W
1st mode Speed	3100 Rpm 51 Hz
2nd mode Speed	3600 Rpm 60 Hz
3rd mode Speed	4200 Rpm 70 Hz
Charging Time	3-4 hours
Work for when fully charged	70 mins
Auto Timer	10 minutes
Weight only the massager	700 gms / 1.54 lbs
weight - full set	1582 gms / 3.48 lbs
Overall size of massager	23.7 x 8.5 x 14.4 cm

Warranty

1. Damage caused by failure to follow the instructions, human factors, force majeure, etc. is not covered by warranty.
2. Maintenance cost will be charged for product out of warranty.
3. You can register for your free warranty on www.EVERTONEPROSAGE.com.

In order for this free warranty to be valid and to allow us to meet your requirements, please register your purchase(s) with evertoneprosage.com. Please visit www.evertoneprosage.com and complete all details. You must have a valid email address to register. Please retain your original receipt to validate your purchase.

Disclaimer

In the process of the preparation of this manual, we have sought to make the content correct and complete. However, there is no guarantee that there are no errors or omissions in the manual.

We reserve the right to change specifications of the products mentioned in this manual at any time without prior declaration. No part of this manual shall be reproduced, transmitted, transcribed or stored in a searchable system or translated into any other languages without the written authorization of our company. The company is not responsible for the loss caused by non-standard or wrong operation.

The **Prosage Thermo** System does not guarantee specific results for individuals who use **Prosage Thermo** due to many factors that can influence health & weight loss. **Prosage Thermo** is not intended to diagnose, treat, cure or prevent any disease, illness or medical condition.

This is not a weight loss system. Individual results may vary. We make no warranties or guarantees. No individual results should be considered typical or guaranteed. The results of **Prosage Thermo** may be temporary. This product has not been evaluated by the FDA. Patent Pending Applied.

To Purchase

To replace or repurchase attachments, or to purchase new attachments and ad-ons, please visit us at www.evertoneprosage.com

Customer Service:

Phone: 516-653-0666

Email: CustomerService@Evertoneusa.com

www.evertoneprosage.com



MC0522 Eng. 8.1.2019



Register now for preferred Savings Club and save 20% on every order
www.igia.com



Join the Club & Save 20% Off Every Order

- As an IGIA™ Savings Club Member, You will receive 20% OFF any item advertised on our website.
- "No Worry" Membership Renewal: A one time fee of \$29.99 for a years worth of savings
- Save BIG on every purchase you make, potentially saving you hundreds of dollars
- We Guarantee our Club will pay for itself, time after time.

Join now & start saving today!

Visit us at **www.IGIA.com**

CONTACT US

Call us 1(516) 653-0666
Email: Sales.IgiaNY@gmail.com

Register for your FREE warranty

Thank you for shopping with us! We are confident that you will be very satisfied with your purchase. To protect your purchase and register for your warranty, please visit us at www.IGIA.com. Located at the top right hand of the website you will find an icon listed "Warranty Registration"; This is where you will fill in the necessary fields to register and receive your warranty.

We appreciate your support and look forward to your future savings! If you have any questions, please don't hesitate to contact us at: Sales.IgiaNY@gmail.com

Visit us at www.IGIA.com

0716