

## 2-IN-1 CELLULITE REMOVER

Item No. 207634

User Guide



Thank you for purchasing the 2-in-1 Cellulite Remover [Fat Breaker]. Please take a moment to read this guide and store it for future reference.


## INTRODUCTION

Stubborn fat cells are high in alpha receptors, which prevents them from releasing their fat content regardless of your diet and exercise. This device works by assisting in breaking down subcutaneous fat using vibration and heat. This device also helps to stimulate collagen by penetrating deep into the skin tissue.

## FEATURES

- Cordless and rechargeable
- Three-speed vibration massage relaxes muscles, increases circulation and helps break up stubborn fat cells
- Optional heat mode helps destroy fat cells and promote the production of collagen
- Charges via USB
- Includes USB charging cable and carrying case

## IDENTIFICATION OF PARTS

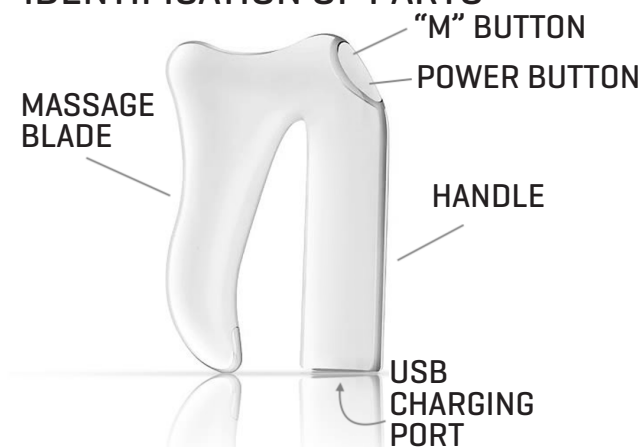

## OPERATING INSTRUCTIONS

- Before use, remove the device from the packaging and charge it using the included USB cable.
- While the device is charging, the BLUE and GREEN LEDs inside the POWER/M button will flash. When the device is fully charged, the BLUE / GREEN LEDs will remain steadily on. Unplug the charging cable when charging is complete.
- Press and hold the POWER button for 2 to 3 seconds to turn the device ON. The BLUE light indicates that the device is on in LOW speed mode.
- Press the POWER button a second time for MEDIUM speed. The light will be GREEN.
- Press the POWER button a third time for HIGH speed. The light will be YELLOW.
- Press the POWER button a fourth time to go back to LOW speed. The light will be BLUE.
- To activate the heat, press the M button. [Please allow 4 to 5 minutes for the device to warm up.] **NOTE:** *The heat function can only be used while the vibration massage is on.*
- Place the device on the target area of skin and apply gentle pressure, slowly and steadily moving the massage blade across the length of the target area. Recommended use: twice a day, 20 minutes per session.
- Press and hold the POWER button for 2 to 3 seconds to turn the device OFF. **NOTE:** *This device is recommended for torso, abdomen, arms and legs only. Use on dry skin only.*

## SPECIFICATIONS

Rated Power	1.2 W
Rated Voltage	5V
Battery Capacity	2000 MAH
Charging Time	5 Hours
Heating Temperature	113°F
Massage Modes	Low / Medium / High with Heat ON or OFF
Product Size	4.2" L x 1.3" W x 5.2" H

## WARNING

- Consult your physician before using any massage device.
- Consult your physician if you have a history of heart problems, have a pacemaker, medical implant or you are currently pregnant or plan to become pregnant. Use of this device is not recommended for those persons.
- Keep this device out of reach of children. It is not a toy.
- Do not use this device on areas of the body other than those specified above.
- Turn this device OFF after each use.
- Do not attempt to open, modify or repair this device. There are no user-serviceable parts.

**PLEASE NOTE: THIS DEVICE IS NOT A WEIGHT LOSS DEVICE.** *We make no warranties or guarantees that you will lose a specific amount of weight for a particular period of time. This product has not been evaluated by the FDA. This device is not intended to diagnose, treat, cure or prevent any disease, illness or medical condition.*