

IGIÁ MICROTONE™

MICROCURRENT FACIAL & BODY MASSAGER



INSTRUCTION MANUAL

1.



SATISFACTION GUARANTEED

Congratulations on your purchase of the Microtone™ Microcurrent Facial & Body Massager. Please read this manual thoroughly before use and keep this manual in a safe place for future reference.

In the process of the preparation of this manual, we have sought to make the content correct and complete. However, there is no guarantee that there are no errors or omissions in the manual. We reserve the right to change specifications of the products mentioned in this manual at any time without prior declaration. No part of this manual shall be reproduced, transmitted, transcribed or stored in a searchable system or translated into any other languages without the written authorization of our company. The company is not responsible for the loss caused by non-standard or wrong operation.

WARNINGS:

1. Do not use this device if you have heart problems, high blood pressure, diabetes, are pregnant, or are a minor.
2. Only use this device as instructed.
3. If you experience pain or discomfort in your skin while using this product, stop using it.
4. Do not let this product fall into water as it can short circuit.
5. Please only use the USB charger that came packaged with this product.

Customer Service:

Phone: (516) 653 0666

Email: OnlineOrders2017@gmail.com


www.IGIA.com



EFFECT

MICROTONE™
MICROCURRENT FACIAL & BODY MASSAGER

1. Using high frequency vibration mode, promote muscle movement, saving time and effort.
2. High efficiency current pulse modulation technology.
3. Stimulation of deep subcutaneous tissue to help with elasticity.
4. The Mechanical vibration and stimulation accelerate and helps target muscle tissue, promote formation of different aerobic exercise, and helps to enhance skin.
5. The use of Blue light on the skin promotes collagen.
6. The use of Y type platinum massage wheel design, allows the device to contour to the face, hips and legs.
7. Help target bags under the eyes.

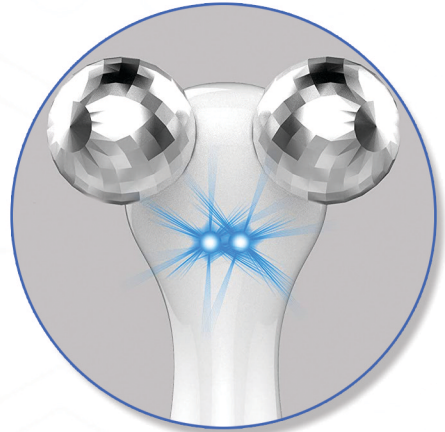

3.

MICROTONE™
MICROCURRENT FACIAL & BODY MASSAGER

SPECS

Name: Microtone™ Massager

Model: BK3473

Product appearance material: ABS

Weight: 190g

Product specifications: 160x87x61mm

Polymer lithium battery: 200mA voltage: 3.7V

Output power: 1W

Charging time: 5-6 hours.

PRODUCT FEATURES

This platinum-coated microcurrent face and body massager uses high frequency vibration to increase the skin's natural elasticity, smoothing fine lines and wrinkles while promoting collagen production. A range of 5 microcurrent modes and 5 vibration intensity levels work together to circulate blood within the skin, bringing fresh oxygen and nutrients to the deepest layers and toning sub-dermal muscles for a firmer, younger looking face and body. Gentle Blue light fluorescent technology further penetrates the cells, discouraging bacterial growth and encouraging collagen production.

The unique Y-shape design with sphered ends targets acupressure points, releasing tension and promoting lymphatic drainage. Roll the device over the forehead, cheeks, arms, legs and torso for full therapeutic effect. It can also be used under the eyes to reduce the appearance of inflammation and dark circles.



MICROTONE MICROCURRENT FACE + BODY MASSAGER WITH BLUELIGHT TECHNOLOGY

MICROTONE™
MICROCURRENT FACIAL & BODY MASSAGER

HOW TO USE

- 1.) Make sure the device is fully charged.
- 2.) Press the ON/OFF button to turn on the device. The LCD display screen will light up. This is where you can monitor your microcurrent mode, microcurrent speed and vibration intensity throughout the therapy session.
- 3.) Press the MODE button to set your microcurrent vibration mode. Choose from 5 settings to find a pulsation interval that feels best.
- 4.) Use the MC button to choose a speed of microcurrent pulsation that feels best. You may wish to start slower and work your way to a faster speed throughout the duration of the session by pressing the MC button to cycle through the speeds.
- 5.) Use the VM button to choose the vibration intensity that feels best. Again, you may choose to change this at any point in the session by pressing the VM button to cycle through the modes.
- 6.) Press the LED button to activate the Bluelight LED fluorescent technology.
- 7.) Hold the device by the ergonomic handle and move the two spheres over the skin's surface using gentle sweeping downward motions. Apply pressure as desired and circular motions to stimulate deeper layers of skin and muscle. Use lightly on the face and around the eyes, as these areas tend to be more delicate. Use the device on each area of the body for as long as desired.
- 8.) When you are ready to end the therapy session, press the ON/OFF button once again. The LCD display screen will turn off. Recharge the device and store until next use. Note that the device has a safety shut-off feature after 20 minutes without operation.

1. Pulse ball
2. LCD display
3. Vibration switch/ Gear key
4. Pulse current button.
5. Bactericidal switch/Fluorescence detection mode.
6. Mode key
7. Power main switch.

