



## Mega Sport Deluxe Shoulder and Waist Massager



Instruction manual

KH254

This user manual provides both operation information of this product and its detailed specifications. Read this manual carefully and understand it thoroughly before operation. This product is designed for home use only.

## **Safety Precautions**

- Persons with following condition must consult the doctor before using this unit:
- Acute disease
- Malignant tumor
- Infectious disease
- Pregnancy
- Cardiac dysfunction
- High fever
- Abnormal blood pressure
- Skin sensory disorders or skin problems
- Receiving medical treatment, especially those feeling discomfort
- Should not be used by persons in the first trimester of pregnancy, fitted with a pacemakers or AICD, or being treated for an existing deep vein thrombosis

## **WARNING**

- If you feel any problems with your body or skin, consult your medical practitioner and follow their instructions
- Do not use this unit while sleeping
- Do not use on your broken skin or damaged sensitive nerve area
- Always starting with the lower level then increase the level gradually
- For safety, do not use this device in high temperature and high humidity environment.

## **Electricity safety attention**

- The product should never be left unattended when plugged in. Unplug power from outlet when it is not in use and before putting on or taking off parts.
- Close supervision is necessary when this product is used by, or near

children, invalids, or disabled persons.


- Please check the voltage and make sure it is complied with the product specification.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Please do not put in and pull out the power with wet hands, to avoid electronic shock or error
- Please do not bind the product with the power cord, to avoid electric shock or fire breaking out because of broken power cord.
- Please do not use the product with broken power line and plug.
- Do not operate this massager under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.

#### **NOTE:**


- Preset working time is 7 minutes for each massage cycle. It'll turn off automatically after 7 minutes. You can turn it on again for more cycles.
- To avoid overheat, it is recommended that you should not use more than 3 cycles continuously. Always let it cool down before you use it for more cycles.

**After use: When the product is not in use, the power adaptor should be unplugged from electrical outlet.**

**PRODUCT FEATURES**

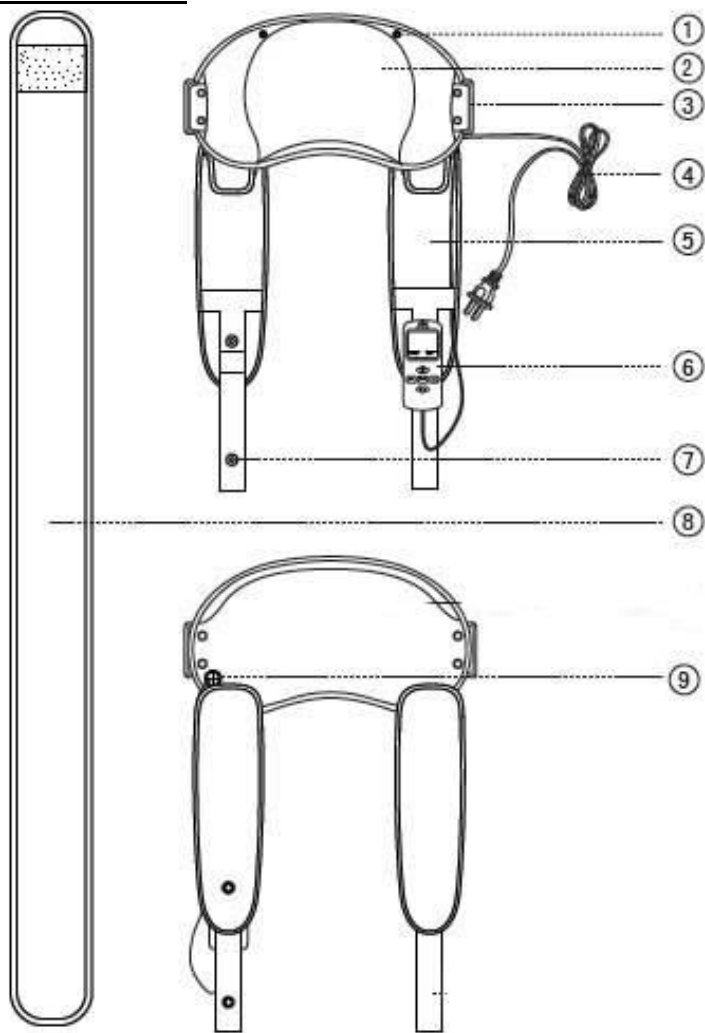

<p>① Ventilation hole</p> <p>② The front of mainframe</p> <p>③ Buckle loop</p> <p>④ Power cord</p> <p>⑤ Shoulder belt</p>	<p>⑥ Handy Controller</p> <p>⑦ Secure strap</p> <p>⑧ Waist Belt</p> <p>⑨ The back of mainframe</p>
---	--

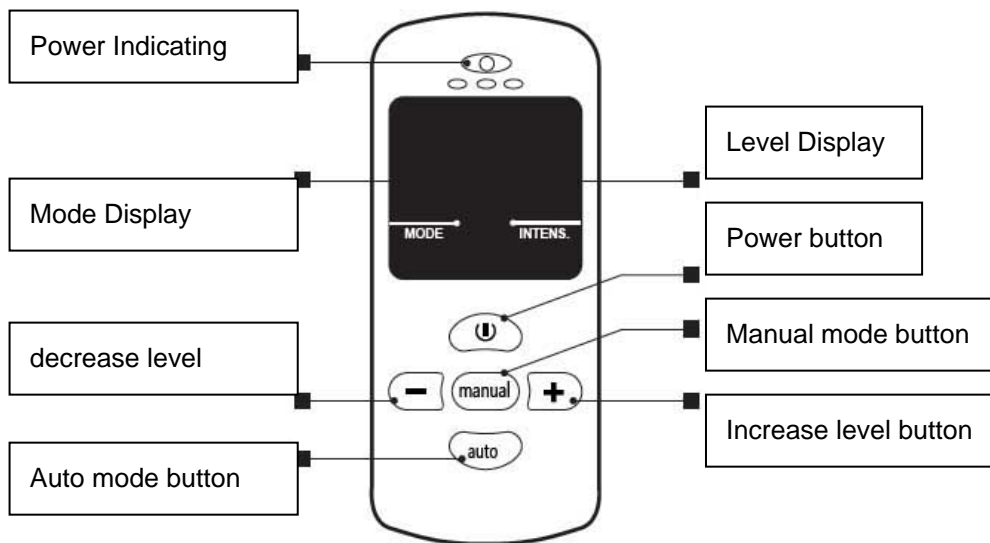


NO. OF MANUAL MODE AVAILABLE	10 (under mode 10, the display will show 0)
NO. OF AUTO MODE AVAILABLE	1 (mode display show E)
LEVEL AVAIL OF EACH MODE	10 (under level 10, the display will show 0)

### **How to operate the unit**


1. Connect the power plug to the wall outlet AC110V/ 60Hz, the power indicator light will be on in red
2. Press "power" button once, the massager starts working at mode 1 and level 1, the LCD displays MODE 1 and INTENS.1,
3. For selecting the mode, press "Manual" button or "Auto" button.
4. To increase or decrease the power level (intensities), press "+" or "-" button to adjust to meet the require level.
5. Press "auto" button, the unit would be under auto mode, and the mode is showing "E", the power level could be adjusted. And the massage pattern will be changed automatically.
6. To avoid the highest power level tapping, it'll always go back to default power level 1 whenever the user changes the mode.
7. Preset working time is 7 minutes for each massage cycle.

### **For using the shoulder belt**

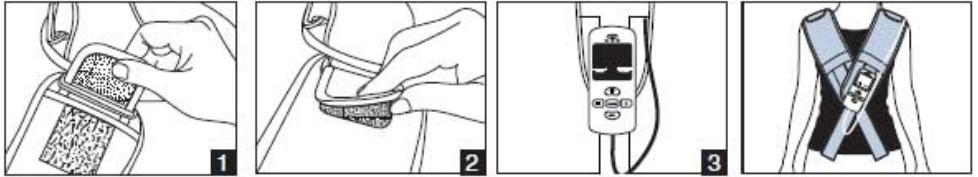

1. The buckles loop is designed on the mainframe, it can be connected with the shoulder belt.
2. Connect the shoulder strap with the mainframe with the Velcro on the shoulder strap.
3. The control could be attached onto the shoulder Strap.
4. The design of secure strap is for user to hold so that it will not drop down easily during massaging.

### **For using the belt**

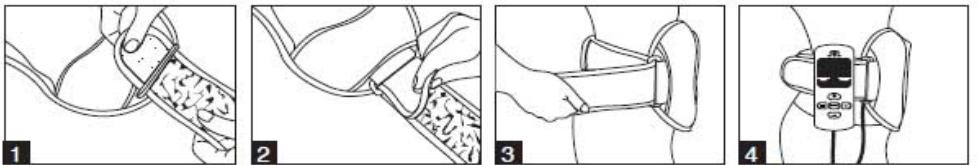

1. The mainframe with buckles loop could be connected to the belt
2. The user could adjust the length with the Velcro tape.
3. To place the massager around the position of body.
4. The control could be attached on the belt when the user massage.

## Using the product with different postures




### Shoulder:

Two kneading head with curve line design helps relieve the stiff muscles, fatigue , shoulder pain and promote the body circulation.




### Back:

Placing the massager on the back helps stimulate the skin and hypoderm, and the conduction to nervous system helps promote the blood circulation, and increase endocrine and immunity.

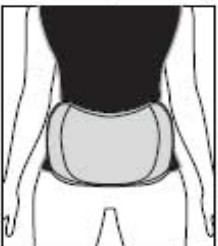


### Waist and belly:

To massage the waist and tummy, helps stretch the tendon and promote the circulation, reduce the fat through the vibration and kneading massage continuously.

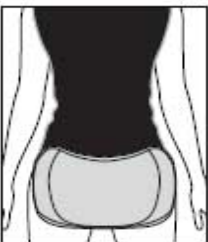

### Buttocks:

To massage the buttocks, it helps stretch the tendon and promote the circulation, reduce the fat through the vibration and kneading massage continuously.

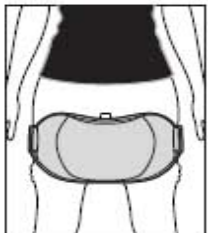

### Thigh and calf:

Placing on the thigh and calf, it helps stimulate the musculature of thigh and calf, reinstate the muscle strength, relieve fatigue and aching muscles.

## SPECIFICATION

Product Name	Carepeutic® MEGA SPROT
Model No	KH254
Power Input	AC110V/ 60Hz
Max. loading current	0.2A
Power Consumption	20 W (110v version)

### Warning :

Do not use the product for more than 30 minutes continuously each time, it may activate the over-heat protection device and stop operating automatically. Please wait approximately 30 minutes for the device to cool down before re-using.

### Usage and Storage Condition

- Store in a dry, dust free location in a temperature between 10°C to 40°C and 30% to 90% relative humidity.
- Do not spill any liquid into the device.
- Before cleaning the product, make sure to turn off the product and unplug the power from outlet.
- Use dry cloth to wipe the PU cover without using any chemicals.

### Limited Warranty:

Please contact seller for warranty. Limited warranty from manufacturer only for original purchaser against manufacturing or workmanship defect under normal use. Due to the hygiene or rehabilitation nature of this product, it is non-returnable.

Carepeutic® is a registered trademark, all rights reserved, © Q.L.M.

[www.ideatomorrow.com](http://www.ideatomorrow.com)

