




Foot massager manual

Please read this manual carefully before using and store in a safe place

● PRODUCTS AND ACCESSORIES

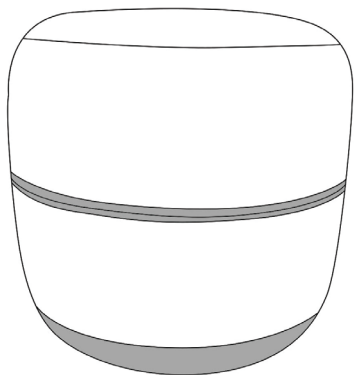

Main part
Upper and bottom cover



charger




Manual



Screw nuts 5


● NAME AND FUNCTION OF EACH COMPONENT

EXTERNAL STRUCTURE

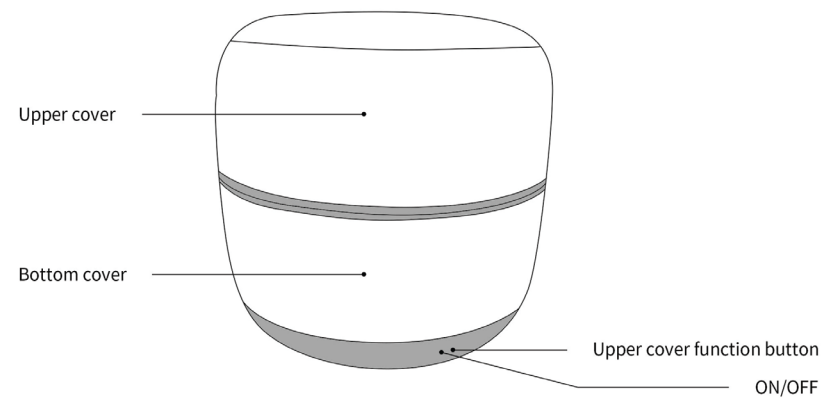


● NAME AND FUNCTION OF EACH COMPONENT

INTERNAL STRUCTURE

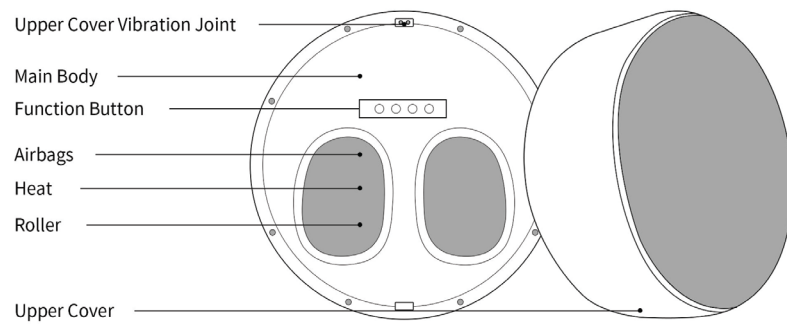

● PREPARATION BEFORE USING

INSTALLATION POSITION · In order to prevent signal interference, please Keep at least 1 m away from the IV, radio and other video and audio devices.

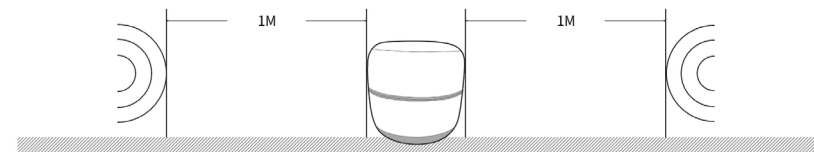

- When charging, straighten the adapter so that it does not get entangled and damaged
- Charging time shall not exceed 3 hours
- Do not charge with a non-original adapter

GROUNDING INSTRUCTIONS




- This product must be grounded to Avoid electric shock
- Do not place the product on uneven ground for use

SURROUNDING ENVIRONMENT

- Keep the foot massager away from damp places (for example: bathroom) to avoid electric shock
- Keep the foot massager away from high temperatures sources (for example: stove) to avoid fire or damage material



● USAGE RECOMMENDATIONS

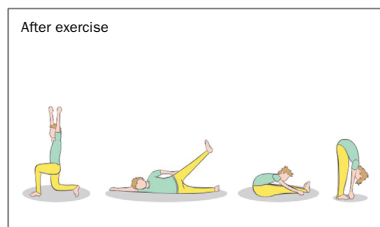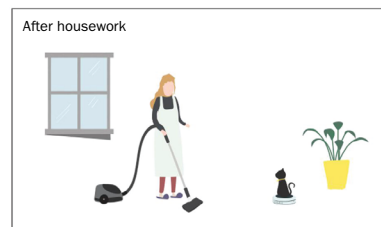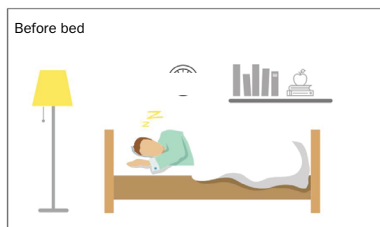

● OPERATION GUIDE OF MASSAGER




POWER CONNECTION

- Lithium battery design, in the case of sufficient power, no need to connect the adapter for use; If the power is low, it can be used after charging



Upper Body Function


- Under the condition that the upper cover is closed, the upper cover cushion can be used for massage

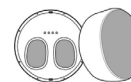


Foot massager

- After opening the top cover, you can use the foot massager


● OPERATION GUIDE OF FOOT MASSAGER

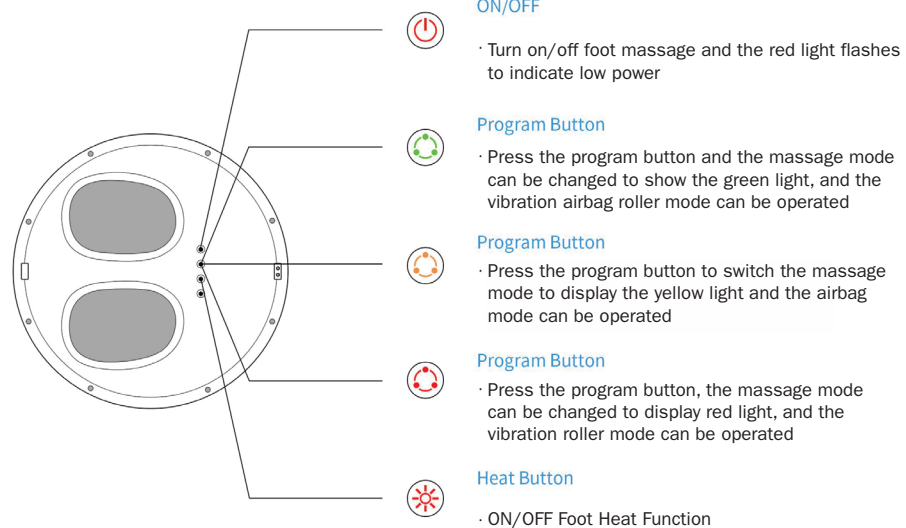

● OPERATION GUIDE OF FOOT MASSAGER

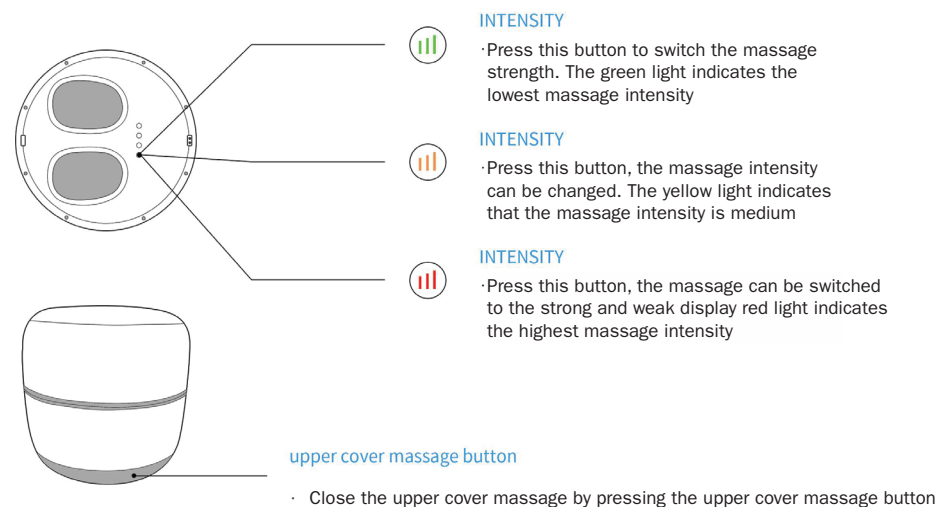

● AFTER USE INSTRUCTION

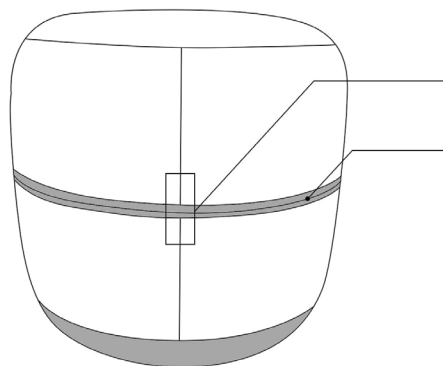

Caution

- Align the joint of the top cover with the joint of the main engine base
- Make sure top cover and main engine base are fully closed

* the upper and lower cover is not fully closed, the machine and the finish will be damaged

● SAFETY PRECAUTIONS



CAUTION

- Before using, please read the "SAFETY PRECAUTIONS" carefully to ensure proper operation

CONTROLLER

- Do not use wet hands to operate the controller or pull the controller wire
- Do not splash water on the controller to avoid function failure
- Do not place heavy objects on the controller

CHARGER

- Do not use adapters outside of the original factory to avoid fire or damage to the product
- The charging time of the adapter shall not exceed 3 hours
- Do not plug the adapter with wet hands to avoid electric shock
- Remove the adapter by holding the plug of the adapter rather than the wire
- Do not squeeze adapter wires by products or other heavy objects
- Do not damage or damage excessive bending, pulling, twisting and knotting adapter wires
- When the adapter wire is damaged, do not disassemble or replace it without permission. A qualified professional should be contacted for replacement and repair

Usage

environment

- Do not use in higher than 40°C

● SAFETY PRECAUTIONS

| | |
|------------------------|--|
| Usage environment | <ul style="list-style-type: none">· Do not use in lower than 0 °C· Do not use this product in damp or dusty environment to avoid electric shock or machine failure· Do not use during a lightning strike and unplug the adapter from the power socket· Do not place the product on uneven ground for use· Do not place this product near a stove or in direct sunlight |
| Prohibited Disassembly | <ul style="list-style-type: none">· Do not disassemble the internal components of this product without permission to avoid the possibility of machine failure· In addition to the repair technician designated by the company, the product must not be disassembled, repaired or modified. If there is any problem, please contact the dealer in time |
| Remove plug | <ul style="list-style-type: none">· Remove the power plug from the socket during maintenance. In addition, do not pull the plug with wet hands to avoid electric shock and injury· Remove the power plug from the socket after use or when not in use· In case of power failure, please immediately cut off the power switch and unplug the power plug from the socket |
| Attention | <ul style="list-style-type: none">· Only people less than 265lb are acceptable |

● SAFETY PRECAUTIONS

| | |
|------------------|---|
| Attention | <ul style="list-style-type: none">· Do not use other therapeutic equipment in the same time of using this unit· This product can only be used by one person. Do not use two or more people at the same time· When using this product, please take off the ankle decorations· Before use, please confirm whether the fabric of all contact parts is damaged. If there is any damage, please stop using immediately and entrust to repair in time· Please stop using it immediately and consult your doctor if you feel any physical abnormalities during use |
| CAUTION | <ul style="list-style-type: none">· It is recommended to use this product for 10 minutes· Make sure there are no animals or other items in the main machine before use· If you do not use this product for a long time, please charge it every other month, otherwise it will damage the battery· Do not stand on this product to avoid injury or malfunction· Do not push or pull on floors and carpets that are easily damaged· This product must be kept in a well-ventilated place and not exposed to strong sunlight |
| MEDICAL WARNINGS | <p>Not recommended for the following people.</p> <ul style="list-style-type: none">· Patients with serious hyperostosis· Swollen lump and patient thrombosis· Patients with serious heart conditions |

● SAFETY PRECAUTIONS

MEDICAL WARNINGS

- Pregnant Women
- People that consume large amounts of alcohol
- People with diabetes
- People suffering from high fever
- Patients with malignant tumors or malignant abscesses
- People with distorted joints

● TROUBLE SHOOTING GUIDE

| Problems/Issues | Possible Cause | Trouble Shooting Method |
|-----------------------------------|--|-------------------------|
| Sounds like tapping or hum sounds | The normal sound of the motor working normally or the friction between the massage wheel and the cloth cover | Normal |
| Unable to turn on | No Battery | Use after charging |

If problems are not from above reasons, please turn off power supply and contact the retailer

Note : Do not dismantle the product by yourself, and our company shall shoulder no responsibility to the damage of the product or the inquiry of people caused by improper use. Thanks for your cooperation

● Kindly Remind

Be sure to check for the following symptoms during year-round use

- Burnt smell
- The electricity is abnormal
- The power cord heats up
- Other abnormal phenomenon
- If the above phenomenon occurs, in order to prevent the occurrence of faults and accidents, please stop using foot massager - cut off the power switch, unplug the power switch, and contact the dealer you bought this product

● SPECIFICATION

| | | | |
|--------------|------------------------------|-----------------|----------|
| Name | Foot Massager | Rated Voltage | 110-220v |
| Package Size | 455X455X455mm | Rated Power | 50/60Hz |
| Color choice | Grey/pink/blue | Noise | ≤55dB |
| G. W | 12Kg | Safety Design | Class I |
| N. W | 9.96Kg | Weight Capacity | 120Kg |
| Standard | GB4706.1-2005 GB4706.10-2008 | | |