

How to Use Hot/Cold Compress Function

- 1 / Ensure the skin is clean and dry and free of makeup or other products.
- 2 / Hold the **power/suction strength button** for 4 seconds to turn the device on.
- 3 / Select the either compression treatment method with the **mode selection button**. We recommend using the **hot compress mode** before extraction to open up the pores, and the **cold compress mode** afterwards to close them.
- 4 / Gently glide the device's compression plate across your entire face and neck. Be careful not to leave the device while on the hot treatment option in one position for too long on your skin in order to avoid a rash or burn. Note that the device itself will get hot in compression mode, and is recommended to not run on this mode for longer than 30 minutes each session.
- 5 / Press and hold the **power/suction strength button** to shut off the device when done.

Hot/Cold Compress

Extractions can be tough on your skin if you're not careful. That's why Clara has a built-in hot/cold compress to make them easy, and painless.

Hot Compress Mode


Why? Use before extractions to help soften blackheads / debris for easier removal.

How? Press the **mode selection button** and select **hot compress mode**.

Cold Compress Mode

Why? Use after extractions to tighten the pores and reduce redness.


How? Press the **mode selection button** and select **cold compress mode**.

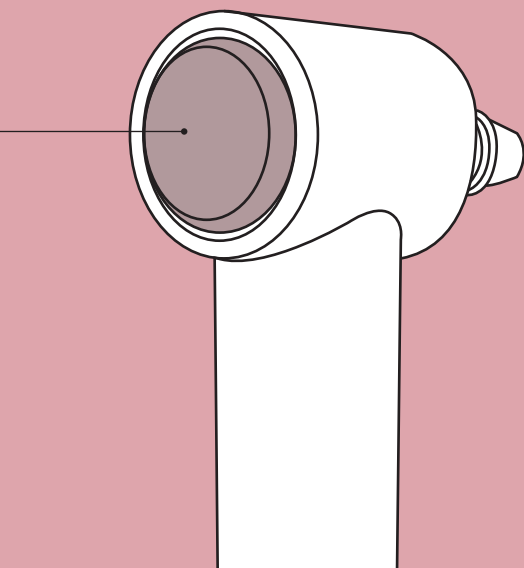

How to Remove Blackhead with Clara Pore Vacuum?


- 1 / Ensure the skin is clean and dry and free of makeup or other products.
- 2 / Hold the **power/suction strength button** for 4 seconds to turn the device on. You will see the spiral icon appear on the screen.
- 3 / We recommend using the HOT compress before extraction to open up the pores, and the COLD compress afterward to close them. Read how to use hot/cold compress section for further instructions.
- 4 / Select one of the eight suction tips. You can add the silicone tip to the small suction tips if you have sensitive skin for a more gentle experience.
- 5 / Press the **power/suction strength button** again to turn on the vacuum suction, and each time you want to increase the suction strength. (Level 1 is the lowest strength, and Level 5 is the strongest.) We recommend starting at Level 1 and increasing after a few sessions as your skin becomes more used to the therapy.
- 6 / You can combine your suction treatment with any of the three light therapy options (R for Red Light, G for Green Light, and B for Blue Light) for optimal results. Press **Mode Selection** to cycle through the different light options. Read about the benefits of each under Light Therapy Section.
- 7 / Press and hold the **power/suction strength button** to shut off the device when done.

Treat each of the areas for up to 5 minutes per session. We recommend that you use a lower suction intensity and work your way up after a few uses.

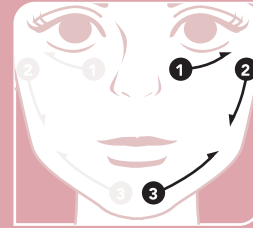


Cheeks and Chin

Gently glide the device across the face following arrows in order.

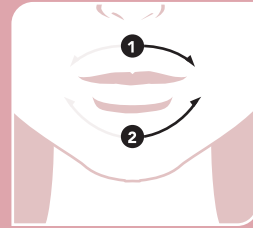


Lip Area

Following arrows and gently glide the device across the face, making sure NOT to let it touch the lips.

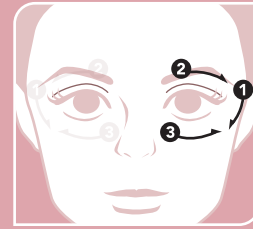


Eye Area

Glide the device downwards following the arrows on each side, ensuring to avoid the eyelid area.

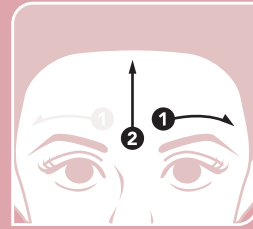


Forehead

Following the arrows in each direction, glide the device outwards and upwards.

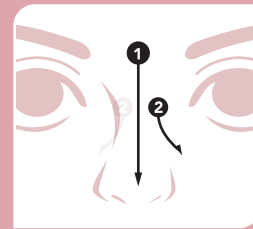

Nose

In a downwards motion, gently glide the device following the arrows.

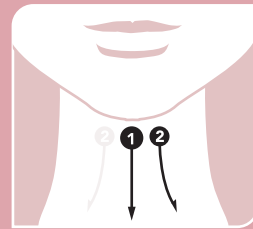

Neck

Use the device in gentle, downwards motions following the arrow in a horizontal direction.

What to Expect with Clara?

Mild Suction Noise

Expect a very low humming noise during use.

Minor Redness

Both during and after your treatments with Clara, you are likely to experience a very mild reddening of the skin. If you experience signs of hypersensitivity after 48 hours, please stop the use of the device immediately and seek the advice of a physician's if needed.

Minor Dryness

Since you will be extracting impurities from the skin, you can anticipate mild dryness immediately after treatment. We recommend applying your favourite moisturizer after each session.

Storage and Maintenance

After use, unplug the device and store it in a dry, safe place. Ensure that the device is never exposed to direct sunlight or submerged underwater.

Discard and replace the suction filter as needed and clean the suction tips after each use.

How to Clean Clara

Ensure that the device is unplugged and powered off. Use a clean, dry cloth to wipe down the device. The suction tips can be cleaned with warm water and gentle soap. Allow tips to air dry thoroughly before storage.

Troubleshooting

The device won't power on

- Confirm that the power cord is appropriately connected to the device, USB charger and to a working power outlet.

The suction is weak

- Confirm that the device has adequate, solid contact with your skin.
- Confirm that the device is on an appropriate suction power level (there are 5 options, which can be switched by pressing the POWER button while the device is turned on.)
- Confirm that the filter is positioned correctly inside of the suction base, and not clogged with debris. Replace with a new one as needed.

Warnings

- 1 / This device is suitable for those ages 16 and above and anyone with reduced physical, sensory or mental capabilities or lack of experience/knowledge if they have been given supervision or instruction concerned the use of this device in a safe way and understanding the potential risks.
- 2 / Always unplug the device when not in use. Do not leave the device unsupervised.
- 3 / Keep the device away from wet surfaces. Do not use near a bathtub, sink, or any household furniture containing water.
- 4 / If the item has fallen into water or other liquids, do not reach for the machine. Unplug immediately.
- 5 / Keep the product and the power cord away from hot surfaces.

- 6 / The device may get hot after extended period of use, particularly with the compression therapy, and is recommended to be powered off after 30 minutes of use.

- 7 / Use this product only for its intended use as described in this manual. Do not use it for anything other than skincare.

- 8 / Use the product with caution. Do not hold the device while powered on one point of the face for an extended period of time to avoid damage or bruising. Ensure that the device is in constant motion across the face throughout use.

- 9 / The following persons should not use this appliance unless approved by a medical professional:

- Women who are pregnant or currently having their menstrual period (usage may cause skin problems due to hormonal imbalance.)
- Persons with heat-sensitive skin or skin diseases such as rashes, atopic dermatitis, etc.
- Persons with skin injuries such as open wounds, swelling, severe sunburns, etc.

- 10 / Keep the product away from children.

- 11 / Clean the device regularly to avoid the accumulation of grease and other residues.

Battery Disposal

Note: Opening and disassembling the device will void its warranty. The process is not reversible.

This action must only be undertaken when the device is ready to be disposed of.

All batteries should be recycled or disposed of as per state or local guidelines, and should not be thrown away in household waste. To remove the battery, cut open the silicone plastic shell and carefully remove the lithium battery. Wear gloves during this process for your safety. Do not dispose of in a fire as the battery will likely explode and cause injury.



Children must not play with this device

Cleaning and user maintenance may not be carried out by children without supervision

Here to Help

We pride ourselves on providing premium beauty tools which is why all of our products are backed by a 24-month warranty. If you have any questions or need further assistance, you can send us an email, connect with us on Live Chat or give us a call! Our friendly support team is here to help.

1 • 800 • 285 • 0183
support@fancii.com




LIVE ONLINE CHAT AT:
fancii.com

CLARA

Tell Me More About Clara

Clara is a multi-functional beauty tool that uses microdermabrasion technology to exfoliate away dead skin cells and suction tips to easily remove blackheads and excess oil from pores. Clara has a built-in hot and cold compress and three light therapy options, to stimulate collagen production, reduce the appearance of wrinkles and reveal healthy-looking skin.

Contents

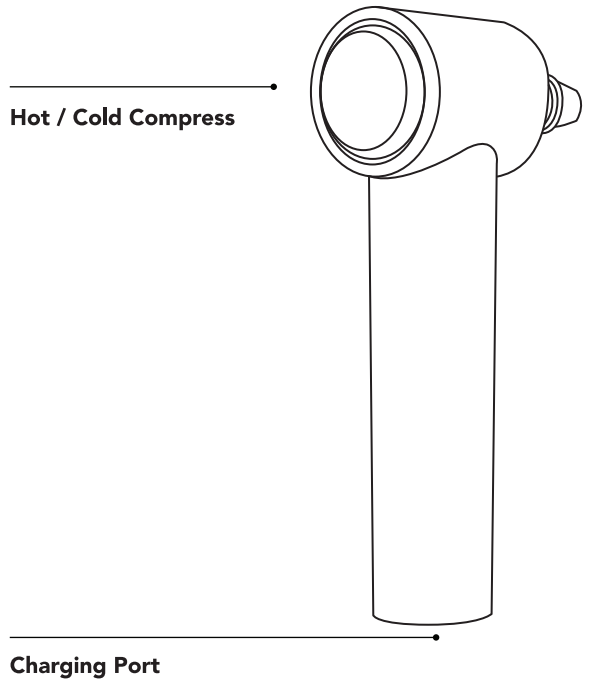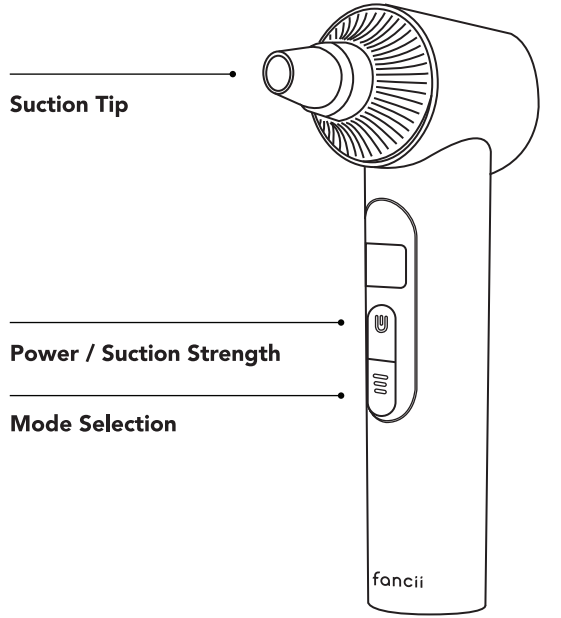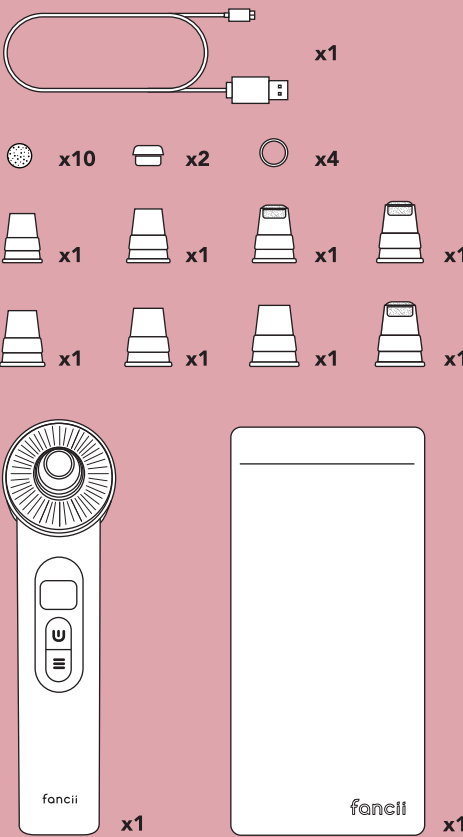

Suction Mode

Suction Mode
This symbol appears when Clara is ready to use.

Suction Strength Level
Set the suction strength level between 1 (weakest) to 5 (strongest) at your desire.

Hot Compress Mode
Use before blackhead removal to help open the pores for easier extractions.

Cold Compress Mode
Used after blackhead removal to shrink and tighten the pores.

Light Therapy Mode

R Red Light Mode
Softens the skin and restores elasticity.


G Green Light Mode
Reduces redness and inflammation.

B Blue Light Mode
Shrinks the pores and rejuvenates the skin.

Setting-Up Clara

- 1 / Remove the Clara device and the other components from the box.
- 2 / Insert in the USB cord (included) to the charging port at the bottom of the device. Plug the USB cord into a reliable power source such as a computer, or USB DC-5V adapter.
- 3 / The battery takes approximately 3.5 hours to fully charge. The device will glow green when it's charged and will be usable for up to 120 minutes before needing to be plugged in again.
- 4 / Attach the rubber ring around the middle of the suction base and insert the foam filter inside.

- 5 / Attach any suction tip to the suction base.

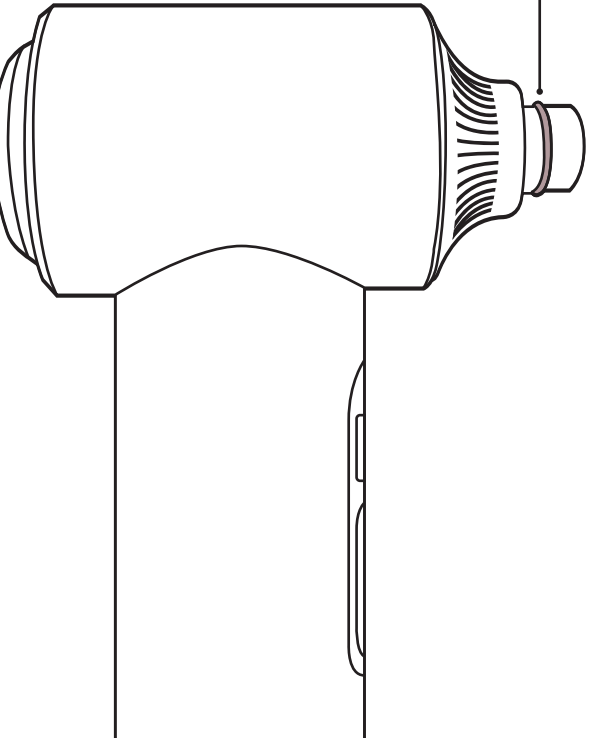

How to Choose the Right Tip?

Clara has a few different size tips, how do you know which to choose? The rule is, the wider the suction mouth, the stronger the suction will be. We recommend starting small and working your way up with repeated treatments. You can also choose the size based on the area you're using it on.

Microdermabrasion Tips

These tips are designed to exfoliate and suction away dead skin cells (and other impurities) that can lead to blocked pores, uneven skin texture, and dull-looking skin.

Microdermabrasion Tip (S)
Ideal for delicate areas around the eyes and nose.



Microdermabrasion Tip (M)
Ideal for the cheeks and jawline.



Microdermabrasion Tip (L)
Ideal for larger areas such as the forehead and neck.



Suction Tips

The suction tips are perfect for removing excess oils and other impurities that lead to clogged pores, inflammation, spots, blackheads, and acne breakouts. Clara includes a variety of sizes to accommodate different areas of the skin.

Suction Tip (XS)
Lowest suction strength that's ideal to start with for delicate areas around the eyes and nose.



Suction Tip (S)
Low suction strength that's ideal for delicate areas around the eyes and nose.



Suction Tip (M)
Medium suction strength that's ideal for the cheeks and jawline.



Suction Tip (L)
Strong suction that's best for larger areas of the face.



Suction Tip (XL)
Strongest suction ideal for areas such as the forehead and neck.



Silicone Tips

Clara comes with two reusable silicone tips that attach to the smaller suction tips and helps reduce suction and provide a smoother surface for those with delicate or sensitive skin.

How to Use LED Therapy

While exfoliating and decongesting your pores choose one of the three LED settings for added benefits.


Red Light
Softens the skin and restores elasticity.



Green Light
Reduces redness and inflammation.



Blue Light
Shrinks the pores and rejuvenates the skin.

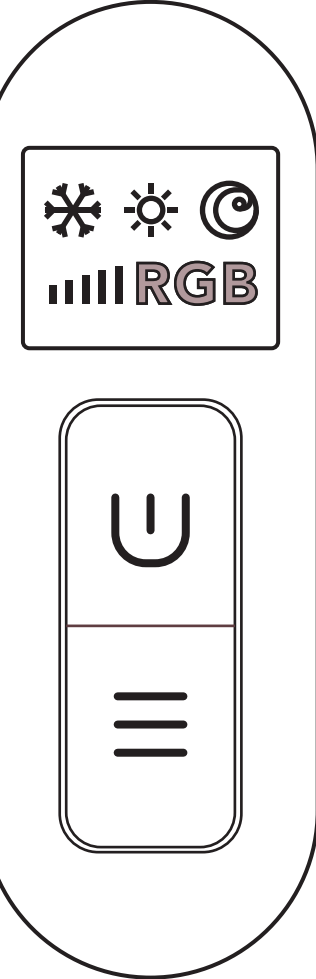

COME FIND US ON
#FanciiAF @FanciiCo



fancii