

IGIA

SONIC PLAQUE REMOVER

USER MANUAL

Whiter
Brighter
Healthier
Teeth!




- Removes Stains
- Whitens Teeth
- Promotes Healthier Gums
- Recommended By Dental Professionals

Recommended
by Dental
Professionals



Gentle Sonic Vibrations

Product structure

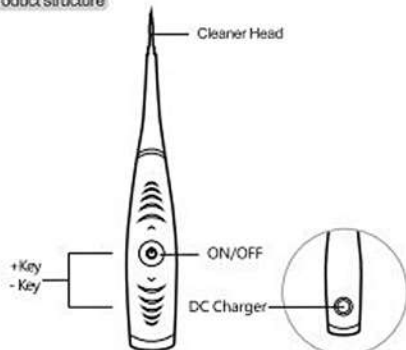

This product is small in shape and convenient to carry. It allows you to clean your teeth at home. This product uses a clean head with the same structure as dentist. It can effectively remove obvious dental calculus, and also have certain cleaning effect on dental plaque and tooth stains.

Product features

- >> Professional grade alloy metal cleaning head, cleaning effect very nice;
- >> Five adjustable, suitable for different needs;
- >> Intelligent high frequency vibration, effective removal of dental calculus and so on;
- >> Whole body design, beautiful and generous;
- >> The battery is durable and practical.

Attentions

- >> Metal head of the tooth cleaner, not replaceable
- >> Do not place products in direct sunlight. Put at -10 C to 40 C indoor environment


- >> Do not wash or soak the product for a long time
- >> Products are restricted to adults only. Please keep them properly after use, avoid children's contact, so that children will not be injured by mistake.
- >> The product contains batteries. Do not place products in a fire or high temperature container to avoid fire or explosion
- >> In case of product quality problem, please do not disassemble yourself. Please return to original plant for repair

Product SPEC

Rated voltage : 3.7V	Rated current : 200mA
Battery capacity : 300mAh	Standby current : ≤20uA
Product size : 183.6x30x22.3mm	Product weight : 53.6g
Charging voltage : 5V	Motor noise : ≤60dB

Usage method

The first use, will be recharged for more than 2 hours until completely discharged, Reference to charging description. (Pic1)

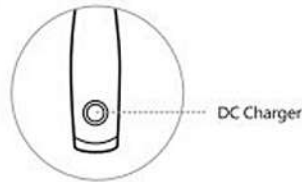

Pic 1

Press the boot button to turn on, adjust the key position through the upper and lower adjustment buttons, and start cleaning the parts that need cleaning. (Pic 2)
After cleaning the teeth, wash the metal head with flowing clean water. Clean the Department to wipe the body of the tooth cleaner.



Pic 2

Charger

1. This product uses the DC interface charging mode, the user takes out the distribution DC interface charging line, the DC interface is inserted into the machine back socket, and the USB interface is inserted at the end of the charger.
2. After entering the charging mode, the tooth cleaner enters the closed state. At this time, the machine will stop working
3. In charging, the light at the bottom is red and full of is green light. (Pic 3)



Pic 3