


Small massager leading brand

Skillful Scalp Massager

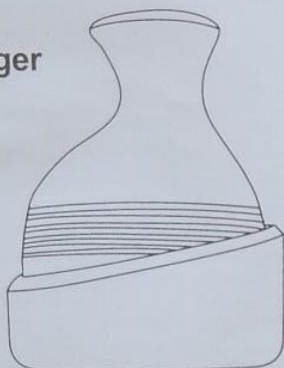

Product Manual

Introduction

Scalp massager, the original 3D meridian massage cluster, with 6 kinds of Chinese medicine massage effect which are opening, closing, thrusting, drawing, rubbing, scratching, can promote blood circulation, improve scalp troubles, pass through the meridians, eliminate fatigue, eliminate dandruff, stop itching, improve sleep ; it can provide the most comfortable meridian massage experience

Scope of application

For hair loss in the head, less hair, poor hair quality, the massager has a significant effect and can significantly improve fatigue.

Applicable object

1. Mental worker
2. people who often feel tired
3. Hair loss in the head, less hair, poor hair
4. College students, middle school students and primary school students
5. People who often have headaches and dizziness
6. people who suffer from insomnia and dreaminess.

Note: People who have the following conditions should not use:

1. Patients with hemorrhagic head.
2. Patient undergoing head surgery within three months

3

Precautions

To ensure that you use the scalp massager correctly and safely, please read the instructions carefully.


1. Do not use the massager at high temperatures;
2. If the head has redness and inflammation, bruises and other discomfort symptoms, or if you have had a head surgery, do not use it;
3. Do not use if you have bleeding disorders such as hemorrhagic disease or hemophilia;
4. For skin problems and obstacles caused by silicone products, please consult a physician before purchasing;
5. Patients with severe stroke, hypertension, and hypotension should use the massager with caution. Please consult a professional physician. If there is any discomfort during use, stop using it immediately.
6. Please avoid using the massager after drinking or when you are not feeling well;
7. Please use the power adapter configured with this product and pay attention to electricity safety;
8. The massager is designed as a home electronic health care product and is not suitable for diagnostic and therapeutic purposes;
9. Do not place the charging base in a humid environment;
10. Do not charge in excessively humid conditions.

Common faults and troubleshooting

Fault conditions	Cause Analysis	Processing method
Power dead	Is the battery exhausted?	Please charge the battery
Massage strength weakened	Is the battery exhausted?	Please charge the battery

5


Product description

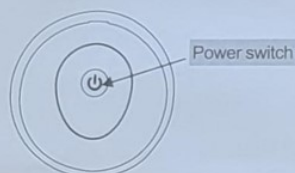

Product top view(Fig. 1)



Product main view(Fig. 2)

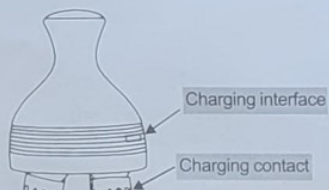

Main unit (Fig. 3)



Charging base (Fig. 4)

2

Method of operation

1. Press the button to turn it on first, the default is slow, the second is fast, the third is reversed quickly, each press is repeated for the first time, press and hold for 3 seconds to shut down.
2. Check if the machine voltage is too low. If the battery red indicator light flashes quickly, the machine should be charged in time.
3. The power adapter is a power adapter configured for the company.

Charging

1. During the charging process, the battery indicator light is always red. When the battery is fully charged, the indicator light turns yellow.
2. Do not use this product during charging (the charging is not supported in the power-on state. The product starts the self-protection program and automatically stops working);
3. The charging time is 2 hours.

Power on, mode selection, shutdown

1. In the shutdown state, short press the key, the massager starts working, and the default slow speed starts after switching the power on. The default time is 10 minutes.
2. After turning on the machine, press the switch key briefly, change the massage mode, and knead slowly. Fast kneading in the forward direction; slow kneading in forward direction; slow kneading in reverse direction; Fast kneading in the reverse direction;
3. In the power on state, long press the switch and the massager is turned off.

4

Warranty card (warranty record)

Name	Date of purchase	
Model		
Address	Color	
Warranty	Warranty content	Remarks

Product certification

Product Name: Scalp Massager
Product model: MZ-5568
Production Date:
Inspectors:
Inspection conclusion: qualified

Precautions: According to the scope and duration of the warranty, The warranty must be repaired by the local dealer store with the purchase invoice and warranty card.

This card is changed to void with no replacement card. (If the dealer store cannot repair the massager or has been relocated, please contact the manufacturer directly.)

We will not provide free warranty service for any malfunction caused by the following users:

1. When this warranty card or warranty certificate is not presented;
2. Man-made improper use of the wrong operation causes damage;
3. Unauthorized repair of the product that caused the product failure.

The company reserves the right of final interpretation of the above content

Please keep this card in a safe place

6