

OVERHEAT PROTECTION SYSTEM

The device will shut down upon operating for 15 mins due to overheating protection. We suggest restarting the device about 5 mins after the device shut down. However, you can restart the device immediately without any issues.

ERROR CODES AND SOLUTIONS

Error codes	Solutions
E1	The device run out of battery. Please charge the battery.
E2	Singal failure between PCB board abd battery. Please contact the service center.
E3	Motor failure. Please contact the sevice center.

CAUTION:

Do not disassemble the product by yourself. For any technical issues, please contact our service center.

WAYFLEX MASSAGER :

Evertone Inc.

Glen Cove NY 11542

Customer Service: cscare2021@gmail.com

www.avibaba.com

Scan for
More Products




Wayflex

Total Wave

Massage Device

ENGLISH

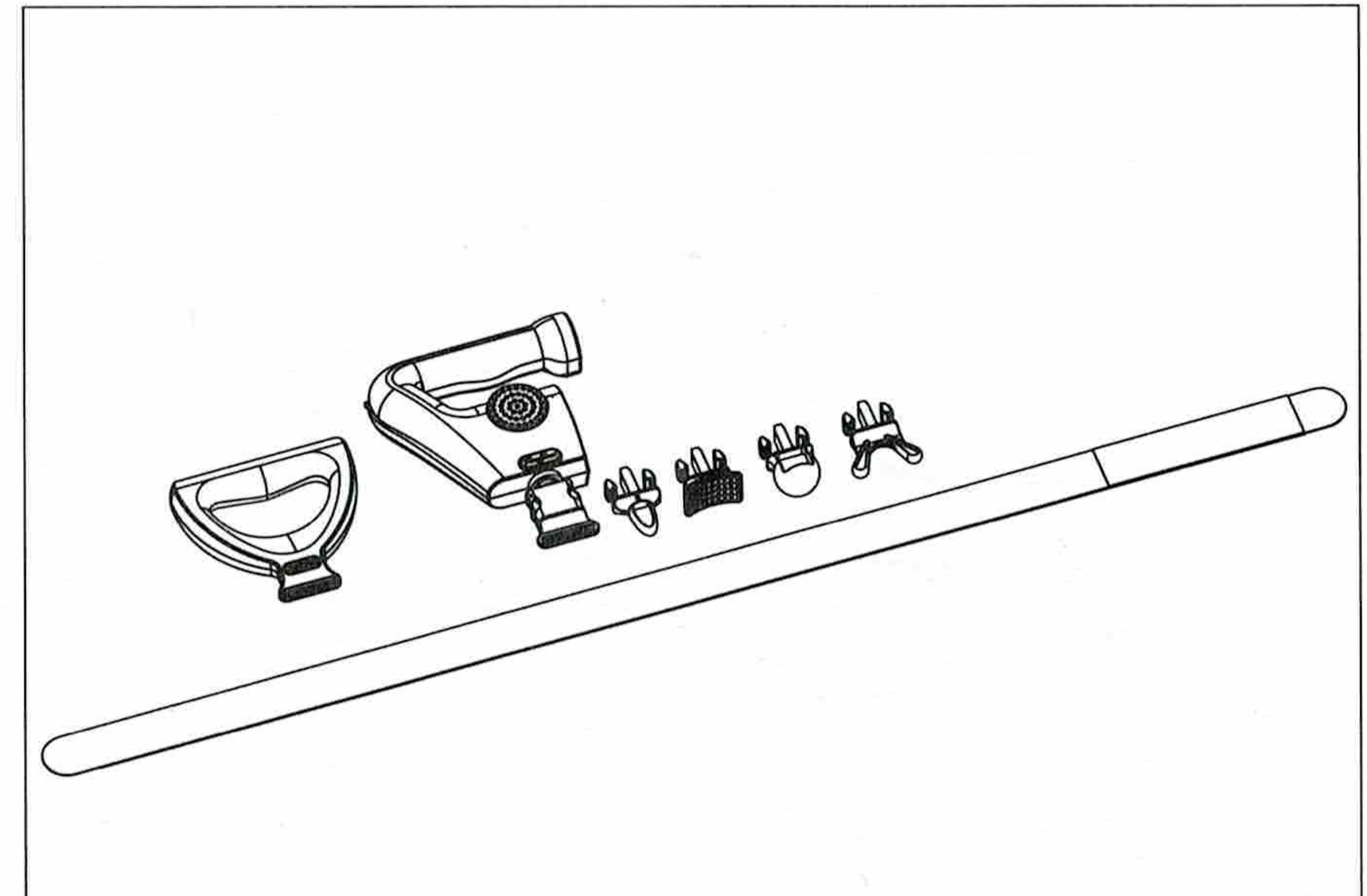

USER'S GUIDE

- ▲ Read all precautions and instructions before using this equipment.
- ▲ After reading, keep this user's guide in a safe place for future reference.

ATTACHMENT DEMONSTRATIONS



The ball head attachment works on large muscle groups for deep massage for areas like your shoulders, waist, back, and legs.

The triangle head is great for relieving stiff necks and typical shoulder tension and can help relieve headaches, and also it is perfectly shaped to relieve those muscles down your calves or your Achilles area, to help gain mobility and flexibility.



The curved head molds to your body and features stimulation nodules that help lactic acid break down in your stiff muscles, and stimulate blood flow to help increase circulation.

The pointed head is for targeting those really sore areas that need some extra attention. It also works on all those acupuncture points on your body for relief.








The adjustable belt suits whatever length you need for the whole body.

SPEED DATA

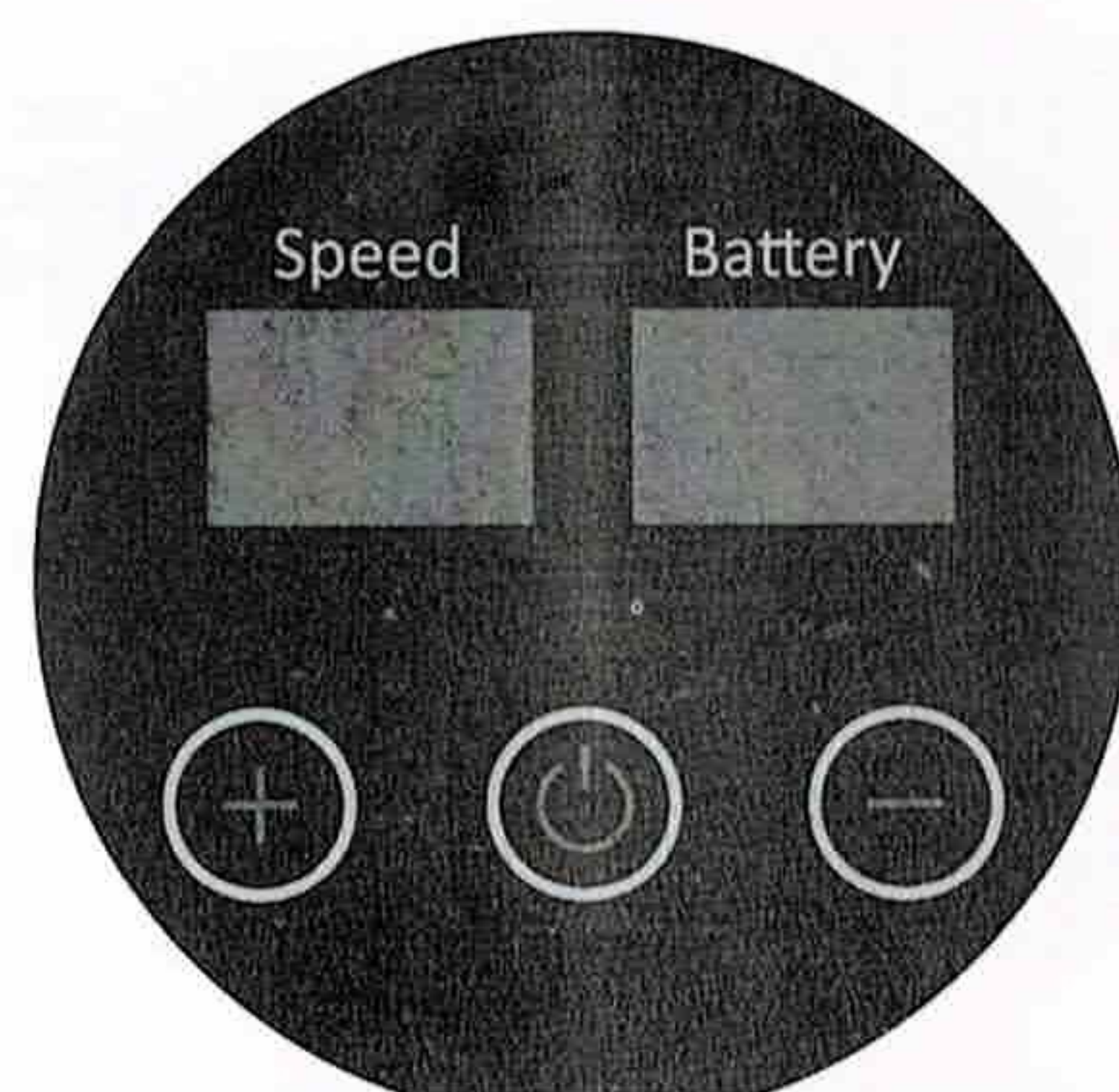
SPEED	VIBRATION PER SECOND
1	14
2	15
3	16
4	17
5	19
6	20
7	21
8	23
9	24
10	26
11	27
12	29
13	30
14	32
15	33
16	34
17	36
18	37
19	39
20	40
21	42
22	43
23	44
24	46
25	47
26	49
27	50
28	52
29	53
30	55

DISPLAY INSTRUCTIONS




	Press  to turn on/off the device
	Press + to increase the current speed
	Press - to decrease the current speed

POWER SWITCH

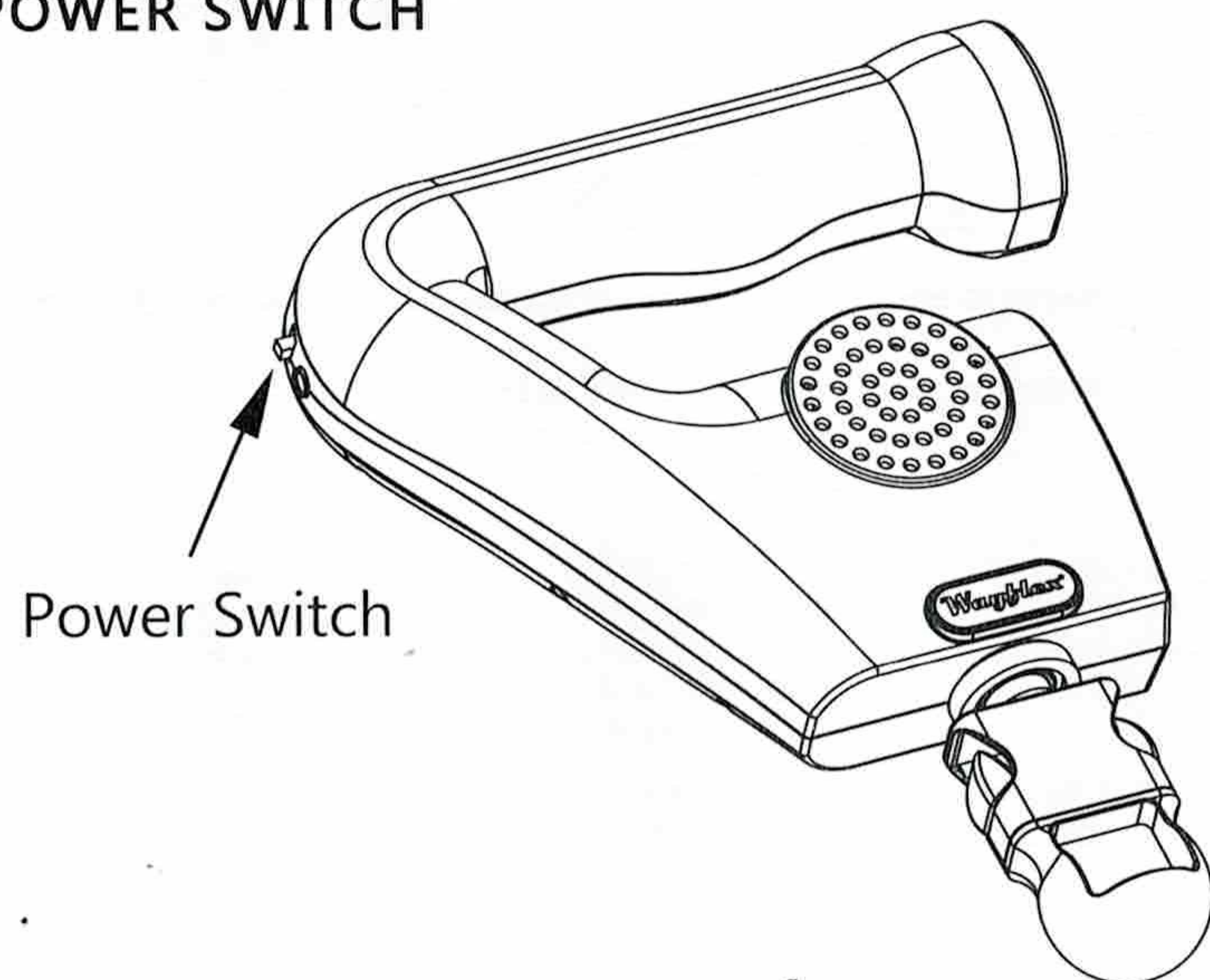


ATTACHMENT DEMONSTRATIONS


Belt attachment

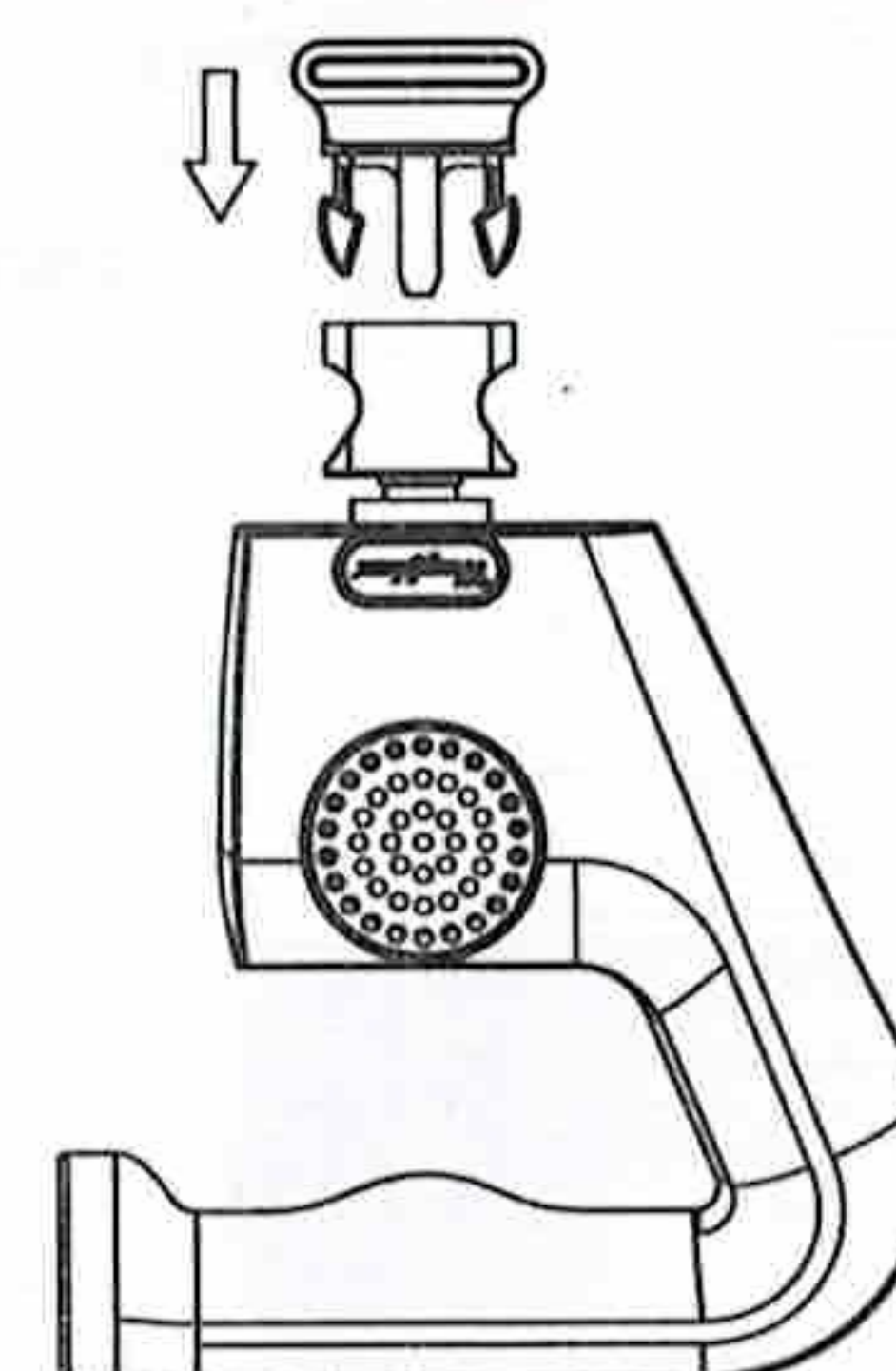


Belt attachment with adjustable belt and belt handle

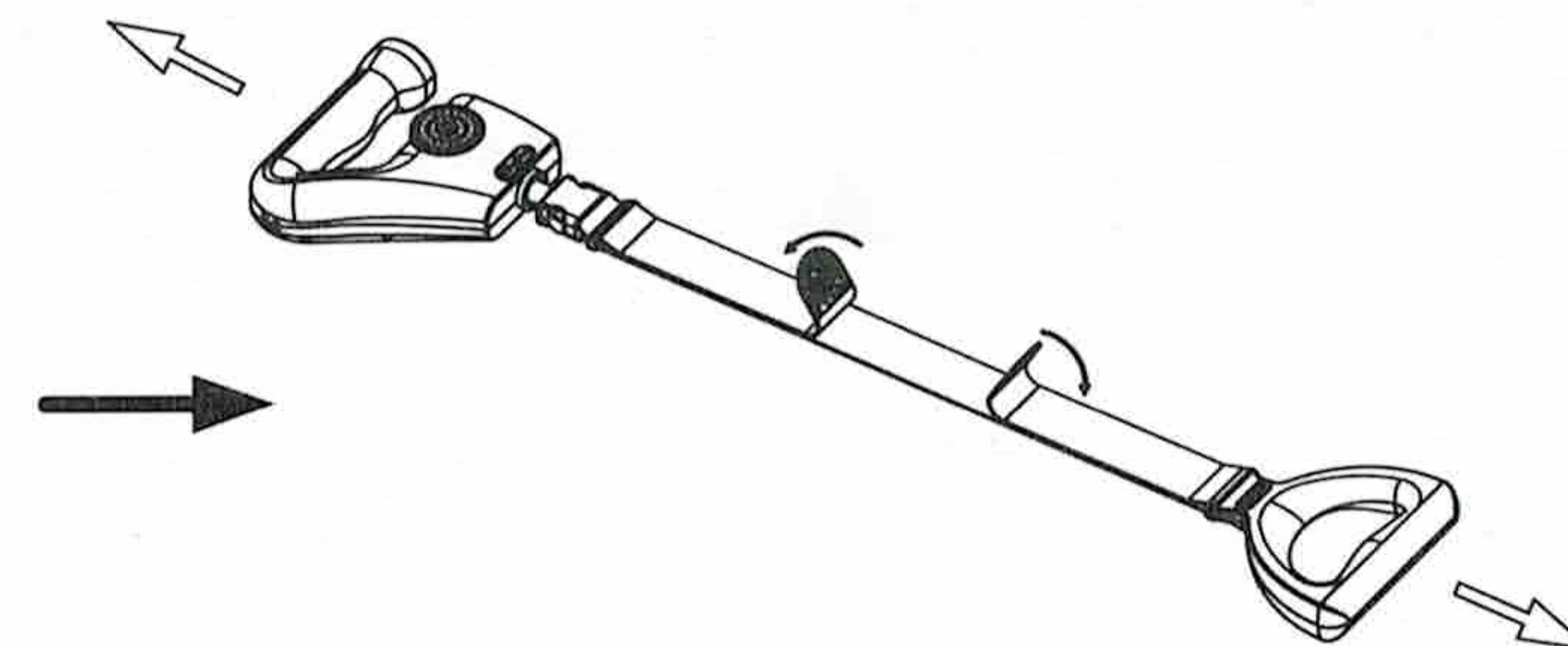


Curved head attachment

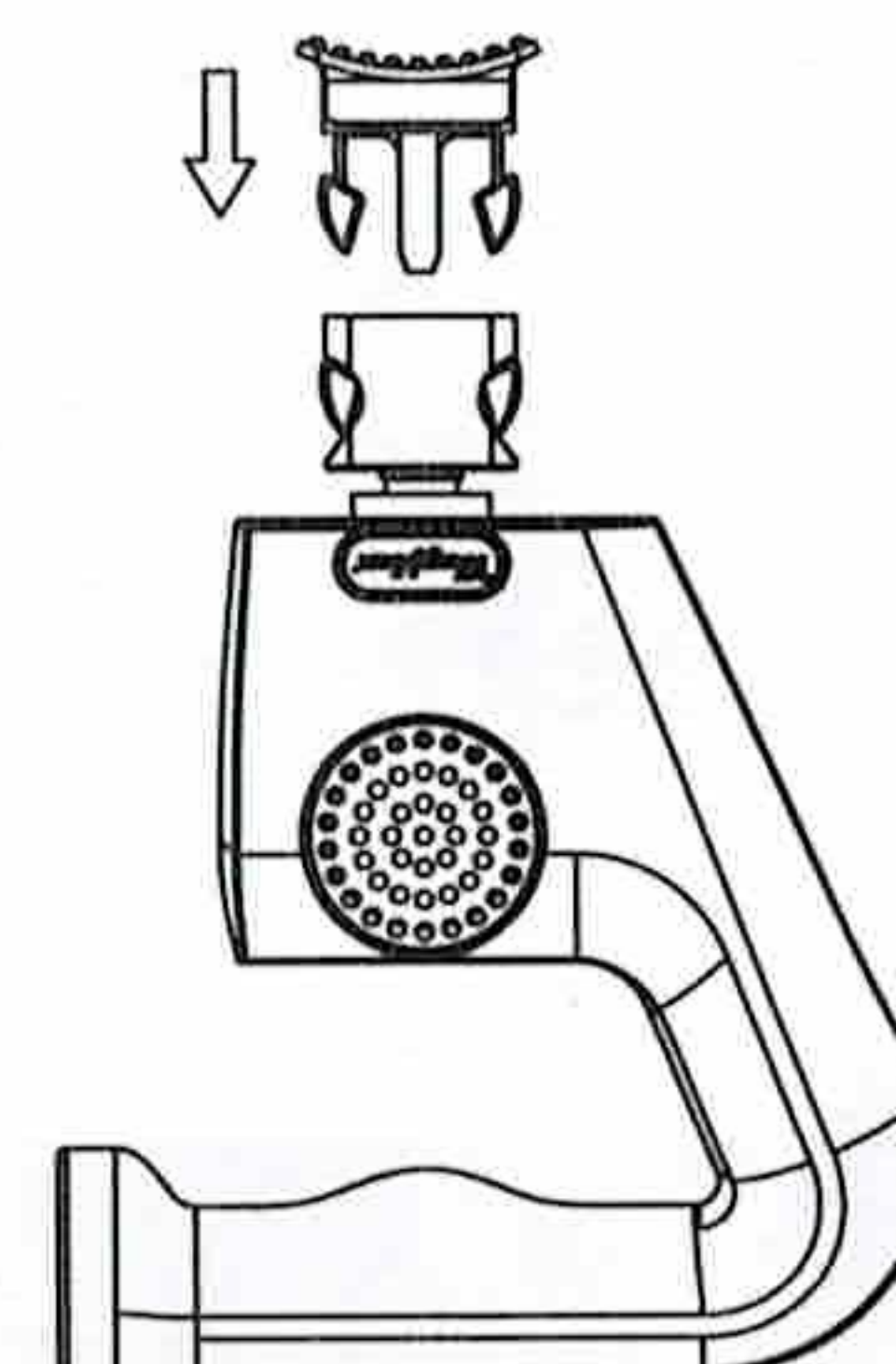


Pointed head attachment

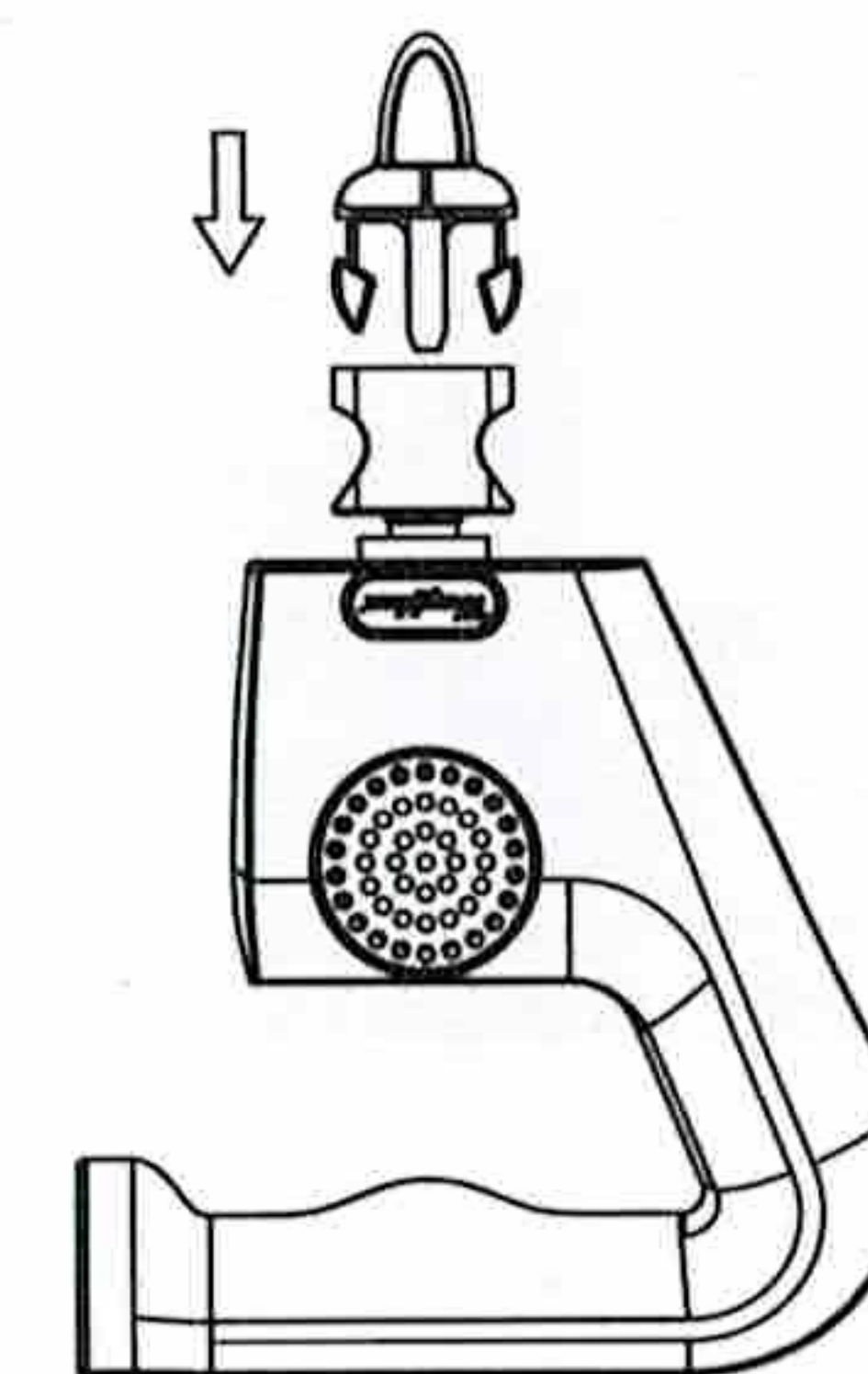

Triangle head attachment

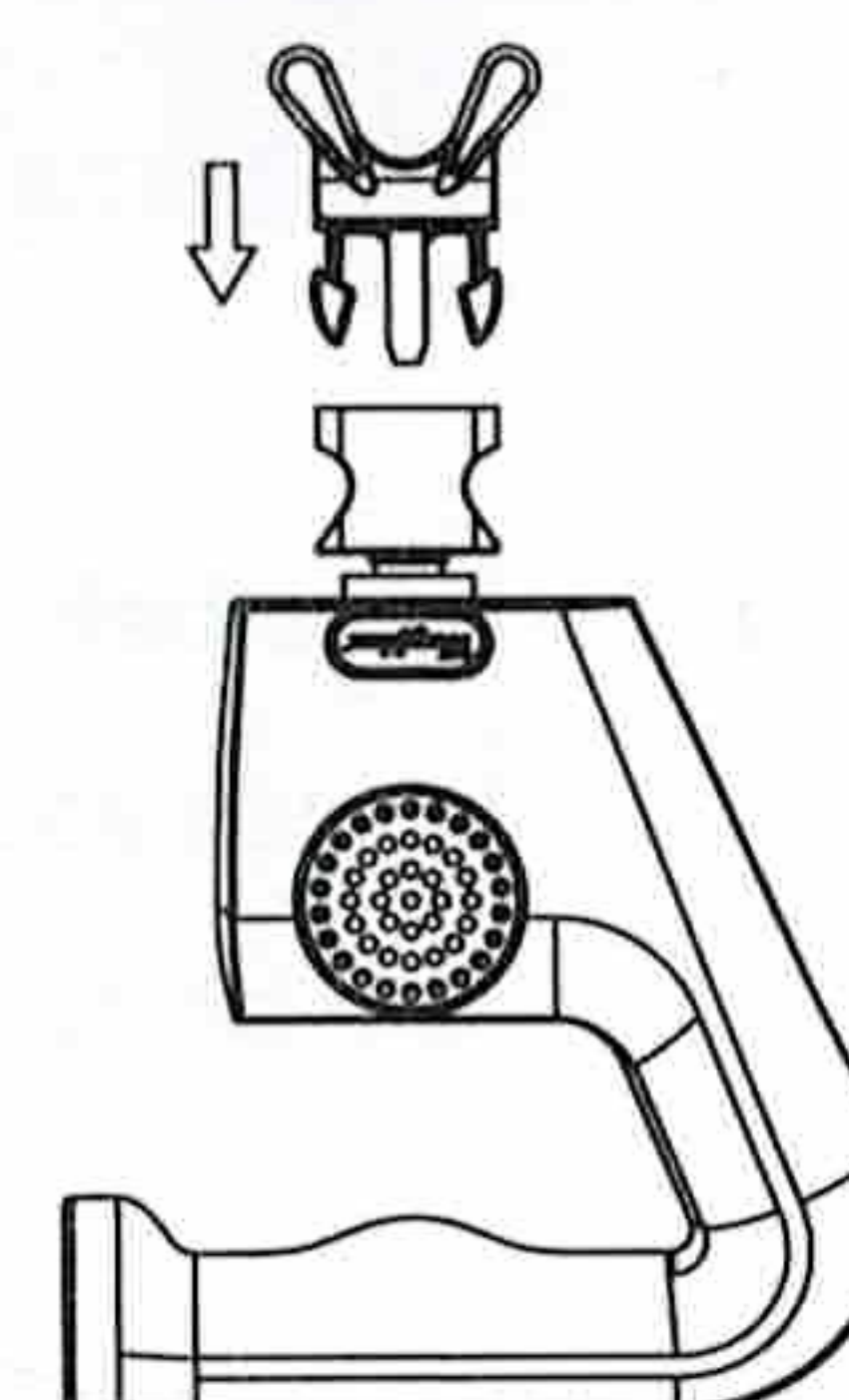

Ball head attachment

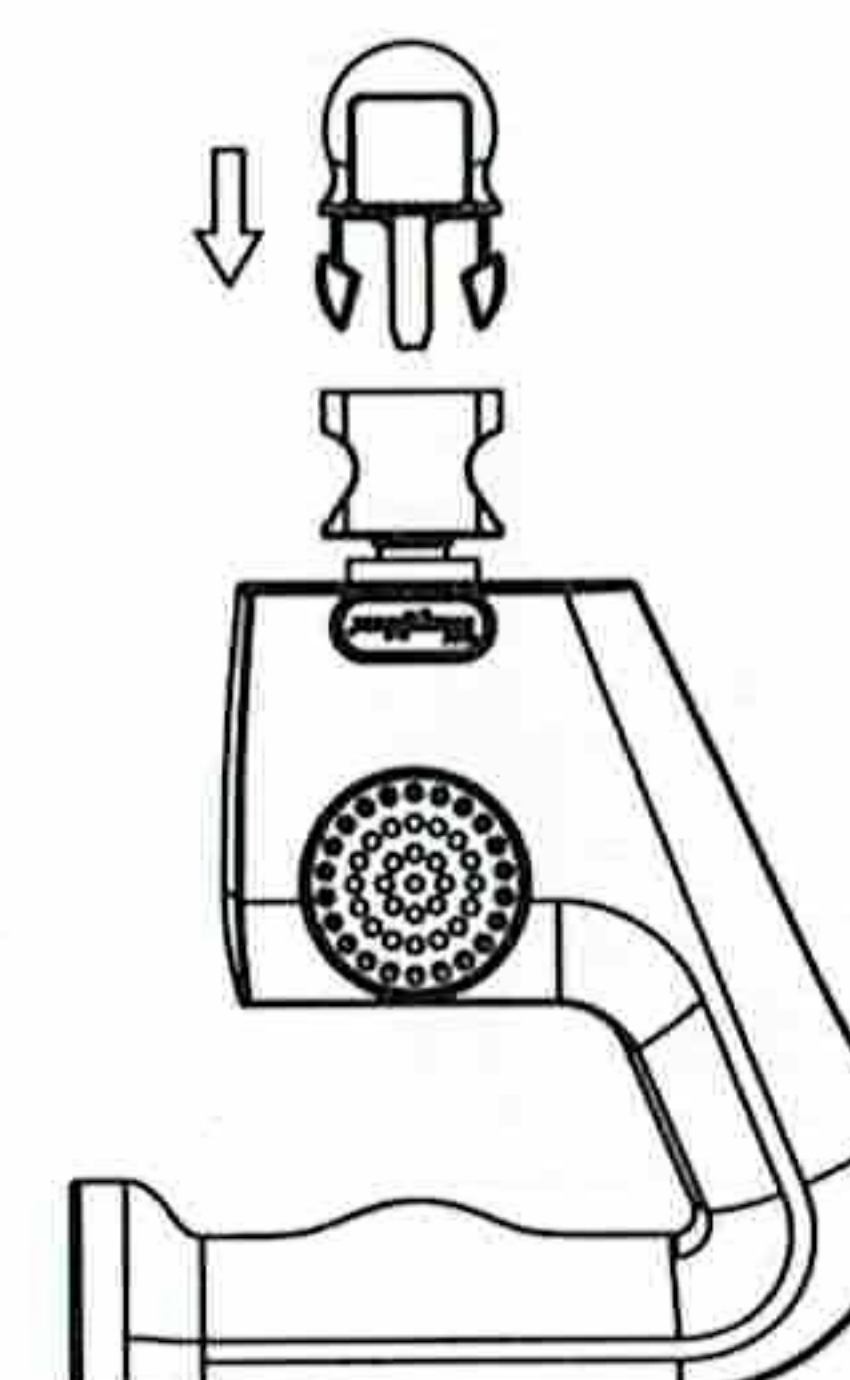

OVERHEAT PROTECTION SYSTEM

The device will shut down upon operating for 15 mins due to overheating protection. We suggest restarting the device about 5 mins after the device shut down. However, you can restart the device immediately without any issues.

ERROR CODES AND SOLUTIONS

Error codes	Solutions
E1	The device run out of battery. Please charge the battery.
E2	Singal failure between PCB board abd battery. Please contact the service center.
E3	Motor failure. Please contact the sevice center.

CAUTION:
Do not disassemble the product by yourself. For any technical issues, please contact our service center.


WAYFLEX MASSAGER :
Evertone Inc.
Glen Cove NY 11542
Customer Service: cscare2021@gmail.com
www.avibaba.com



Wayflex
Total Wave

Massage Device

ENGLISH

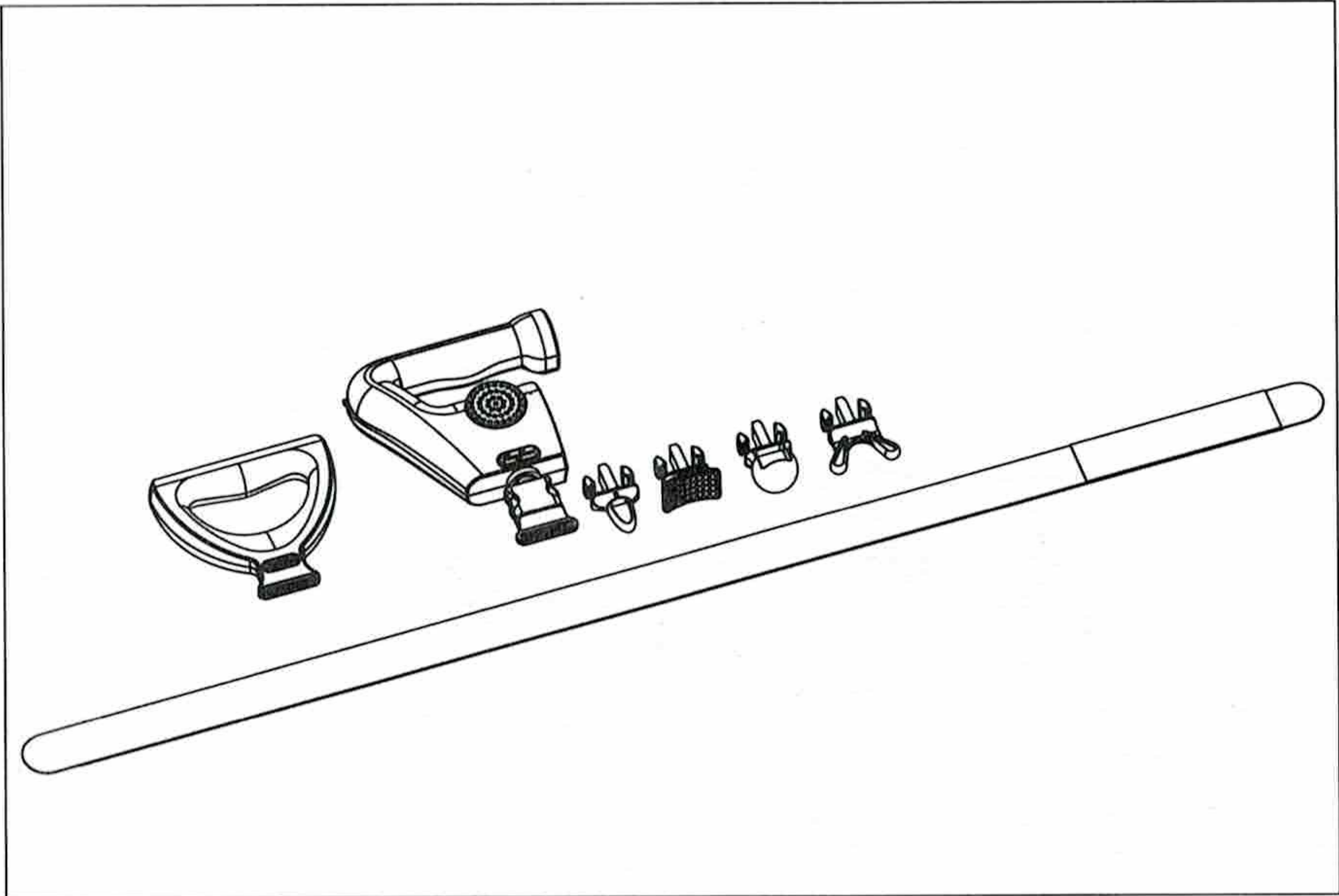

USER'S GUIDE

- ▲ Read all precautions and instructions before using this equipment.
- ▲ After reading, keep this user's guide in a safe place for future reference.