

NOTHING IS MORE IMPORTANT TO PEOPLE THEIR HEALTH.

OPERATION MANUAL SOLAR FLARE FAR-INFRARED CAPSULE

BEFORE OPERATING THE UNIT, PLEASE READ THIS MANUAL
THOROUGHLY AND RETAIN IT FOR FUTURE REFERENCE

BENEFITS

- FAR-INFRARED HEAT THERAPY.
- IMPROVES BLOOD CIRCULATION.
- STRENGTHENS CARDIOVASCULAR SYSTEM.
- EASES JOINT PAIN AND STIFFNESS.
- BURNS HUNDREDS OF CALORIES.
- RELAXES MUSCLES AND INCREASES FLEXIBILITY.
- REDUCES PAIN AND FATIGUE.
- DEEPLY CLEANSSES AND PURIFIES THE SKIN.
- CONTRIBUTES TO HEALTHY WEIGHT LOSS.
- CLEANSSES THE BODY THROUGH PERSPIRATION.
- IMPROVES CELLULAR VITALITY AND PROMOTES OVERALL WELLNESS.
- DIFFERENT TEMPERATURES MAY BE USED FOR DIFFERENT FUNCTIONS.

*30-35 degrees: soothing, relaxing, and stimulating a sense of well-being.

**35-40 degrees: increased soothing effect and muscle relaxation.

***40-50 degrees: as an adjunct to physical therapy, such as for rheumatoid arthritis.

****50-75 degrees : increased sweating for weight loss and deep heat therapy.

- **Please read the instructions carefully before operating this machine, and please pay special attention to safety issues and warning sections.**
- **Please read the safety precautions and proper use of the product, and note the following signs:**



WARNINGS

The machine cannot be modified! Do not break down and repair, as this may cause fire, electric shock, or injury.	Do not damage the wire: Excessive bending, pulling, or putting heavy objects on the wire will damage it and may cause fire and electric shock.
Do not unplug if hands are wet; it may cause injuries or electric shock.	Do not use a voltage other than stipulated. May cause fire and damage.
Do not use if the wire is damaged, it may cause injury or electric shock.	Clean the dust above the plug regularly. If dust accumulates on the plug, it may cause poor contact.
This product is for indoor usage for health care and slimming only. Do not use for other purposes.	Do not place machine in an unstable position. The impact of an earthquake or other external disruption may result in injury or damage to the product.

WARNING	POSSIBILITY OF BODILY INJURY.	
Caution	May injure or dama the product	
Signs meaning	A	Warning or note
	○	Prohibited
	◯	Must be implemented



CAUTION

- Children and elderly people should be supervised when using the capsule.
- Please stop using device if you feel unwell or uncomfortable.
- Stay away from fire or other heat sources when usingit.
- Please turn off and unplug the product after using.





PLEASE NOTE THE SAFETY GUIDELINES.

When you leave the device unattended or unused for a long time, please unplug it.



Please do not spray any volatile substance.



Please do not leave this product in sunshine.



Keep product away from children.



Please do not use this product in a wet environment



Please do not place this product nearby fire or flame.



Please do not dismantle this product.



Warning! Any person experiencing the following symptoms should consult their physician before use.

- Any person with skin diseases, such as peeling or rash.
- Any person with splanchnic diseases.
- Any person with a heart defibrillator or other implantable devices.
- Any person with heart disease.
- Any person who has skeletal fragility or damage to bones, especially the backbone.
- Any person with a malignancy.
- Any person who does not feel well or is being treated for sickness or disease.
- Any person recovering from surgery.
- Any person with osteoporosis.
- Any person whose blood pressure is abnormal.
- Pregnant women.
- Any person experiencing a headache.

OPERATION

- You can use essential oils, slimming creams, and skin care creams with the device to enhance beneficial effects. IF TOPICALS INCLUDE CAMPHOR, MENTHOL, OR OTHER MEDICATIONS, USE ONLY AT LOW TEMPERATURES AS FAR INFRARED CAN ELEVATE EFFECTS AND/OR SENSATIONS AND POSSIBLY CAUSE DISCOMFORT.
- Connect the power to the controller and to the dome carefully.
- Press the power button on the controller, press the "M" button to choose the temperature or time mode, and then set the temperature and time by pushing the "+" and "-" buttons.
- Time settings are from 0-9 hours; normally, 1 hour is ideal.
- The temperature of the capsule will be a little lower than you set because you do not contact the heating wire directly; the heat transfers from the air.
- When the device reaches the desired temperature, the internal thermostat will control the temperature safely for the user.

NOTES

1. Always use this product in the presence of others to ensure proper safety.
2. Do not use this product within half an hour before and after meals.
3. Unplug the product while not in use.
4. Stay far away from moisture when using this product.
5. The storage location should be dry and away from direct sunlight.
6. Use a soft cloth to clean; avoid using strong or abrasive cleaners.

PRODUCT CARE AND STORAGE

1. Do not use sharp objects when handling this product.
2. After use, clean the device and store it away from direct sunlight.
3. Keep the product clean; wipe it with a soft cloth.
4. Do not place product in water, near flames, chemical solutions, or in humid places.
5. Do not place heavy objects on this product.

OPERATION INSTRUCTION OF CONTROL PANEL


Controller

Temperature display from 35-75 degrees

M is the toggle button for time and temperature

Increase Button

Decrease Button

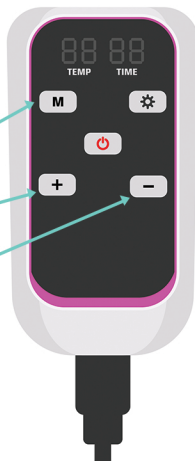

TECHNICAL DETAILS

Product Packaging Specifications Qty: 1 pc / box

Net weight: 5.5 kg **Gross weight:** 6.5 kg

Carton size: 65*45*61cm

Specifications Voltage: 110-220V **Power:** 400W

STORAGE ENVIRONMENT

Temperature: 0-60C **Humidity:** 85

Working Environment: Temperature: 0-40C

Humidity: 85%