




尺寸:135*88mm
材质:105克铜板纸
工艺:拆页(4页)

Mini Portable Oral Irrigator Instruction



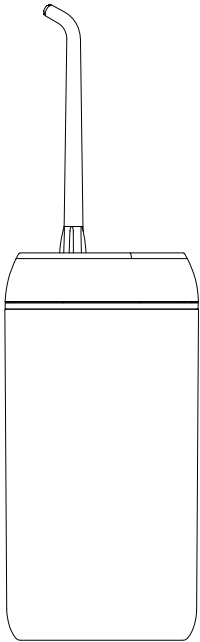


USB charging
Long-lasting use




3 modes
comfortable
flushing teeth



FDA food grade
safety confirmation

FC CE RoHS

SPECIFICATION

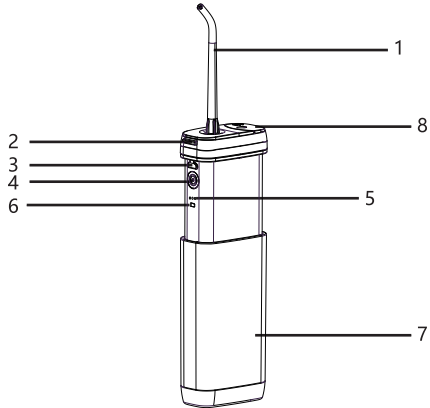

1.Nozzle 2.Nozzle button 3.Charging port
4.Power switch 5.Indicator Light 6.Battery indicator
7.Water tank 8.Water inlet and nozzle storage

Package inside with: Oral irrigator x1 Nozzle x4 Charging cable x1
Instructions(including product certification and warranty card)x1

CHARGING INSTRUCTION

Due to the gums are not get used to the product's stimulation, the gums may bleed at the initial using. (please use low mode first),the bleeding will be stopped in a week.If the situation last for more than two weeks,it may be caused by gum problem(periodontitis etc.).In this case,please stop using the product and consult your dentist.

SAFETY INSTRUCTIONS

DANGER

1. If you do not know the operating process and safety instructions, please do not rush to use the product to avoid hurts caused by mis-operation.
2. Do not disassemble the product or change the parts of the product to avoid damage on the product in case may caused in-valid warranty.
3. Do not use above 50°C hot water into the water tank for cleaning mouth to avoid scald
4. If the power cable or plug is damaged, please replace it before using.
5. Do not put the product into the fire or the heating, and do not charge, use or leave it at high temperature.

WARNING

1. Do not use this product for any other purpose other than oral cleaning.
2. Do not push the nozzle hard on the gums or teeth to avoid hurting on the gums.
3. Do not immerse this product into water for a long time to avoid failure operation or cause fire, electric shock and other dangers.
4. Do not damage, destroy or modify the power cable, and do not stretch, bend or twist the power cable. Do not press heavy objects on the power cable or get it caught between objects, and keep the the product and cable away from hot objects.
5. Do not store the product into a hot or humid places, and avoid direct sunlight.
6. Do not charge the product in hot or humid places, such as the bathroom or toilet.
7. If not in use of the product for a long time, please unplug the power plug and charging plug to avoid fire or danger.

NOTICE

1. There are no replacements or repairs for this product. If repair required, please return it to the seller.
2. Do not use without a nozzle.
3. This product is not suitable for children and infants.unless they are supervised or guided by their safety responsible supervisor on how to use this product.
Keep out of reach for children or infants. Do not use it as a toy.
4. If you have periodontal disease or have had oral surgery in the past 2 months, please consult your dentist before using.
5. Do not add mouthwash into the water tank which will cause corrosion and the tank will be adhered to the device body can result in tank opening failed

Product warranty card

Product Name	
Country of Original	
Warranty period	
Invoice No.	
Date of purchasing	

Product Certification

Product Name: Mini Portable Oral Irrigator
Mode No.: MS6P
Production Time:
QC Inspector::
Inspection Result: Qualified

TECHNICAL PARAMETERS


Product Name	Mini portable oral irrigator	Product model	MS6P
Nozzle diameter	0.6mm	Pulse frequency	1600times/minutes
Water pressure	40~140PSI	Water tank capacity	140ml
Battery capacity	1100mA	Noise	≤72db
Input parameter	5V±1A	Rated power	8W
Charging port	Type-c	Charging time	About 4h
Waterproof Leve	IPX7	Water pressure mode	3 mode
Product size (Height*Width*thickness)	124*62*32mm	Net weight	183g

FAULT DIAGNOSIS


Problem	Possible reasons	Solution
No Working	Newly purchased or not used for more than three months	Please charge it for more than 4 hours
Low water pressure when working	Something inside the nozzle and blocked the nozzle into deformation.	Replace the nozzle
Battery Only can be used for a short time when fully charged.	The battery is at the end of its useful life	Contact your seller
	Insufficient charging times	It refers to recharging for more than 4 hours before using
The charging indicator light does not light up when charging	The power plug is not properly inserted.	Please connect the power correctly.
	Product charging circuit is faulty	Please contact your seller

PRODUCT INSTRUCTION


1.Install the nozzle

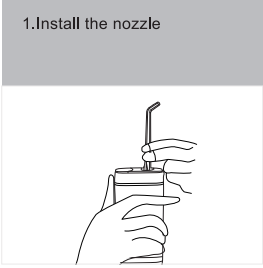


2.Open the tank and the knob to fill water

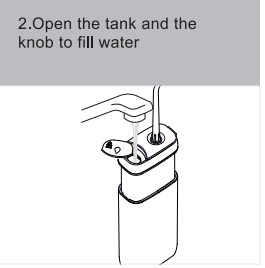


3.Put the nozzle into the mouth and press the power button to start flossing

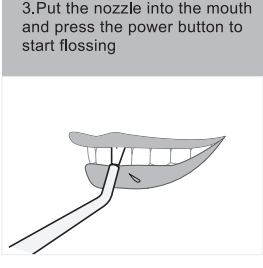


4.Long press the power switch button to switch the mode

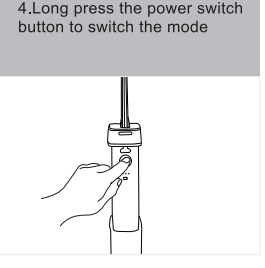

5.Storage after using: remove nozzle, open water tank knob, and store the nozzle

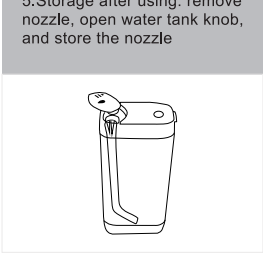


6.Regular cleaning : unattach the tank box and clean the tank

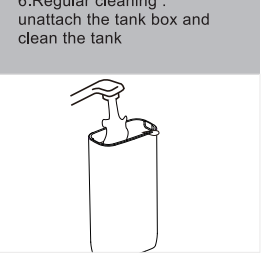

MODE INDICATION

Mode	Mode switch	Mode indicato color	Tar get user	Water pressure
Gentle mode	Press the switch to start flossing	Green	Gentle water pressure,gentle cleaning,suitable for sensitive teeth and first time users	40-60PSI
Pulse mode	Press the switch and hold for 0.5' to switch to pulse mode	Blue	Variable frequency pulse pressure alternate the strong and the weak mode can deeply clean and massage the gums	60-110PSI
Strong mode	Press the switch and hold for 0.5' to switch to strong mode	White	Strong water pressure, efficient cleaning, meeting the daily cleaning for the most of the people.	120-140PSI
Battery indicator	1.The indicator turns red when run out of power 2. The indicator flashes red when charging 3.when fully charged, the indicator turns green.			

PRODUCT MAINTENANCE

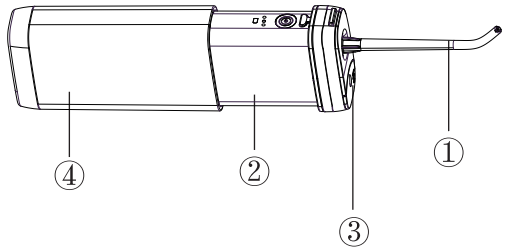

1.Nozzle

Press the nozzle button and remove the nozzle, wash it with water and then wipe dry with a cloth.
The dentist recommends to change the nozzle every 6 months.
If the top of the nozzle is deformed,replace it with a new one to avoid malfunction.

2.Product body

Remove the dirt with a cloth dampened with liquid soap, then wipe clean with a dry cloth. If the product body is dirty, wash it with water.
Do not wipe the product with thinner, gasoline, alcohol etc. Do not wash with hot water above 50°C.Please wipe dry after cleaning.

3.Water inlet

Wipe off the dirt, otherwise it will cause a leakage.

4.Water tank

Open the tank and clean it regularly.
Wash the internal with a brush or other object and wipe with a dry cloth.
When not in use of the product for a long time,please keep the water tank dry.Thank you for purchasing our product.
If you have any questions,please contact your seller.