




EMS Slipper

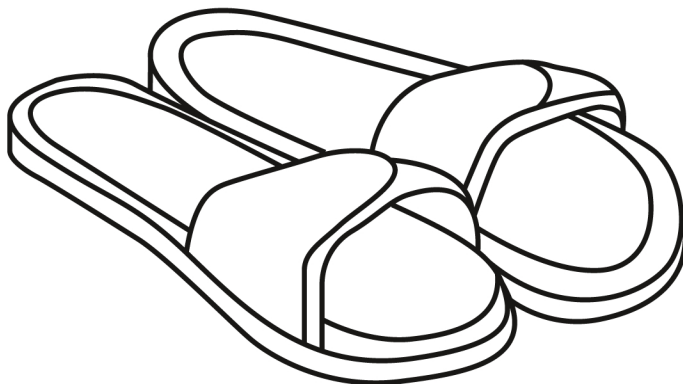

Please read the instructions carefully before using the product

Warranty Card

Product number: _____

Product number: _____ - After sales

Date: _____

Date of purchase: _____

Warranty: _____

Warranty

1. The product will be guaranteed free of charge within one year from the date of sale, and will be repaired at cost after the warranty period.
2. Products damaged by man-made or force majeure factors are not covered by the warranty.
3. If the product is not operated according to the instruction manual or disassembled by itself, it will not be replaced or guaranteed.
4. For after-sales questions, please contact the dealer or local agent.

Product Instructions

1. It is recommended to use within 15 minutes, and the product will automatically shut down after 15 minutes by default
2. It is recommended to use it before going to bed/lunch break for better results, especially in winter when the temperature of the soles of the feet is low and the body's cold delays entering the sleep state.
3. The product is divided into left and right, each with a separate switch, and the left switch and the right switch need to be operated at the same time. 4. The EMS function is only available on the right host, if necessary, adjust the EMS operation on the right.

Little Knowledge About Foot Care

It is still necessary to massage and care the acupoints on the soles of the feet. Because the blood circulation on the soles of the feet is poor and the temperature of the skin is low, the ancients knew how to wash their feet with hot water every day and massage their soles with their hands. This is regarded as health care. Tips for keeping healthy. In fact, the reason is very simple, if the blood circulation of both feet is good, and the transportation line from the heart to the remote parts remains unimpeded it will also drive the blood circulation of the whole body.



Body blood and lymphatic circulation. This is the origin of foot massage.

Precautions for product use

- 1 This product is only for home use, it is forbidden to use it as daily shoes.
2. The product is equipped with energy batteries and electronic circuits.
It is strictly forbidden to use it in a place with humidity control or in a place with water to avoid damage to the product.
3. The material that fits the skin in the store is made of PU leather.
If you need to clean it, please use a wet wipe to clean it. It is better to clean the bottom of the skin before use.
4. It is recommended not to take power when using the product. The effect is more obvious when sitting down.
5. It is recommended that the elderly or children use the product under the guidance of the buyer/reader.
6. The number of products is one size fits all. It is normal if there is an abnormal size. The acupoints are the same for everyone, so there is no size difference. The product mainly promotes blood circulation through heating and massage. It is a massage and health care product, which is very suitable for family use .

Product storage and maintenancet

1. Charge the product before use,
2. After the product is used, put the product away separately. It is not recommended to put it together with the household program to avoid others from using it as particles.
3. It is widely recommended that the electrode sheet massage the baby's position and clean the part that needs to be massaged with a wet towel before use.
4. The protective film needs to be removed before the electrode sheet is used, and the film needs to be pasted back to the electrode sheet after use to prevent the electrode sheet from being damaged. 5. If the electrode sheet is used too many times and the viscosity is visibly affected, you can use full water to clean the anti-skinned water glue.


Features

Press and hold the touch switch to switch the machine, and press it to change the mode, mod

1. vibration + fever (EMS low-level), mode
2. non-motion + heat picture (EMS mid-range), mode
- 3 non-motion + hot drink (EMS document).

All 3 modes have EMS function. If necessary, please install the electrode sheet on the wire first, stick the electrode sheet on the position to be massaged, and then twist the wire in the EMS socket in the product.

Introduction to product parts

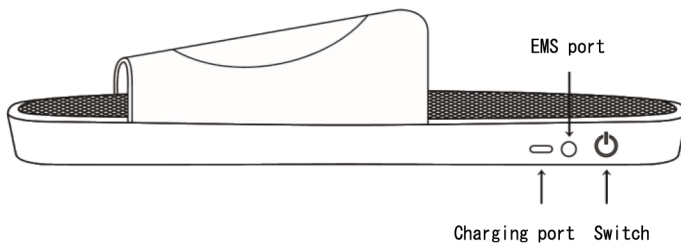

For Any Customer Related Inquiry, Phone: 516-653-0666, Email: cscare2021@gmail.com

Family health starts from the bottom of the feet!

Caring for family

Smart health

