



Turbo-logy Shiatsu Foot Massager with Heat Therapy



Model: KH263



3D Rotating Foot Massage

3D rotating massage, combining the shiatsu and kneading technology, can stimulate the acupuncture points and acupressure regions of the foot and bring you great relaxation and soothe your tired, aching feet and relieve tightness and pain on ankles.



Soothing Heat Therapy

The soothing heated therapy can improve blood circulation, which will deliver the nutrients to all parts of body and bring you optimum state of good health. Moreover, it can alleviate accumulation of lactic acid and relieve muscle pain and stiffness.

IMPORTANT SAFETY INSTRUCTIONS

When using the product, all instructions and basic precautions should always be followed, including the following:

DANGER - To reduce the risk of electric shock:

Always unplug this product from the electrical outlet immediately after use and before cleaning.

WARNING — To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) The product should never be left unattended when plugged in. Unplug from outlet when it is not in use and before putting on or taking off parts.
- 2) NEVER put your feet on massager if your feet are wet.
- 3) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 4) Close supervision is necessary when this product is used by, or near children, invalids, or disabled persons.
- 5) Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.

- 6) Never operate this product if it is dropped, has a damaged cord or plug, a malfunction, or it is dropped into water. Return the product to service centre for examination and repair.
- 7) Do not carry this product by its supply cord or use cord as a handle.
- 8) Keep the cord away from heated surfaces.
- 9) Never operate the product with the air openings blocked. Keep the air openings free of lint, hair and the like.
- 10) Never drop or insert any object into any openings.
- 11) Do not use outdoors.
- 12) Do not operate where aerosol (spray) products are being used or where oxygen is administered.
- 13) To disconnect, switch off all controls and remove plug from outlet.
- 14) The product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or with the lack of experience and knowledge - unless they have been given supervision or instruction concerning use of the product by a person responsible for their safety.
- 15) Children should be supervised to ensure that they do not misuse the product.
- 16) Use heated surfaces carefully. May cause serious burns. Do not use over insensitive

skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.

- 17) Any individual who is pregnant, has a pacemaker, suffers from diabetes, or suffers from malignant tumors, excessive bleeding or tuberculosis should consult with a physician before using this device. This device is designed for home use only.

SAVE THESE INSTRUCTIONS

NOTE: This product is intended for household use only

Name and function of the machine parts



Operation Guide

There is a main power switch in rear of the massager. Make sure the main power switch is in OFF position before you plug in the power adaptor.



1. Plug in the power adaptor to the massager and to the electrical outlet.
2. Switch the **main power switch** in rear to ON position.
3. Press the **ON/OFF** once (front of massager).
The product will automatically start its kneading massage and heated therapy at speed level 3 for 15 minutes.
4. Press the **ON/OFF** again, the heated therapy function will stop and the light on the rotating rollers is off, while the kneading massage continues.
5. Press the **Speed button (UP & DOWN)** to increase or decrease massaging intensity.
(there are 5 intensities).
6. Press the **ON/OFF** third times, the massaging

will stop and the power is OFF.


7. If you do not turn off the power during massaging, the device will stop automatically after 15 minutes. You'll have to press the ON/OFF again in order to turn it on if you want to massage again.
8. This foot massager is equipped with thermal protector. The massager will be turned off automatically when its temperature reaches 185 Fahrenheit degrees. Let it cool down and you can start it again.

After use: When the product is not in use, switch the main power switch in rear to OFF position and the power adaptor should be unplugged from both massager and electrical outlet.

Overheated protection:

This foot massager is equipped with thermal protector. The massager will be turned off automatically when its temperature reaches 185 Fahrenheit degrees.

Massage position:

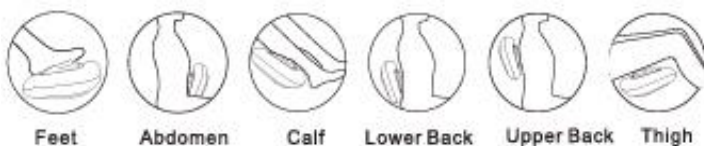

To clean and maintain the product:

1. Cleaning

- Before you clean it, please make sure you switch the main power switch in rear to OFF position and the power adaptor should be unplugged from both massager and electrical outlet.
- Please use dry cloth to clean the product.



Attention

Please make sure that do not use benzene or thinner and other aggressive detergents to clean the product. Otherwise, it may cause damage to the product.

To clean the fabric:

- To open the fabric cover, use your fingers to press on the Red Color plastic ring around the fabric cover and rotate it in opposite direction as indicated on the red color ring. And the red color ring with the fabric will pop up. Take it out slightly.

- Wash it with running water carefully.
- Don't brush it and **let it dry completely by cool air before putting it back.**



2. Maintenance

- When the product is not in use, switch the main power switch in rear to OFF position and the power adaptor should be unplugged from both massager and electrical outlet.
- Keep the product clean and away from the wet environment.

3. More Safety Instruction:

Read the instruction manual thoroughly before using this product. When using electrical products, specifically when children are present, basic safety precaution should always be followed, including the following:

- Always unplug the appliance's power cord from the electrical outlet immediately after using it and before cleaning it.
- Close supervision is necessary when this appliance is used by, close or next to children, invalids, or disabled person. This product should not be operated by child under age 12.
- DO NOT operate this appliance if it has a damaged power cord or damaged electric plug; if it is not working properly; and if it has been dropped or damage or dropped into water. Return the product to service center for examination and repair. Never try to repair this device by yourselves.
- Keep power cord away from heated surfaces.
- The appliance should never be left unattended when it is connected to power source. Unplug the power cord from power outlet when it is not in use or before putting on or taking off parts.
- This device is NOT for outdoor use.
- This device is NOT for use with wet feet.
- DO NOT immerse this device in water or with any other liquid.
- To disconnect the power source, turn all controls to the off position on control panel, and then remove the plug from power outlet.
- Use this appliance only for its intended use as described in this manual.
- DO NOT use attachments that are not

recommended by us with this device.

- DO NOT attempt to stand up on this device. Use it only when you are being seated on a chair or sofa.
- Connect this appliance to properly founded electric outlet only.
- This appliance is not intended for use by person (including children) with reduces physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction in regards to the use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- DO NOT use the device in the bathing room or where you taking show. Use it only in a dry area and with dry feet. Keep it away from humid area, wet floor. Do not spray water onto the massager of the rotating disc or the outside case areas around the appliance.
- DO NOT attempt to plug or unplug unit with power source while your feet are on the machine.
- If you have any concerns regarding your health, consult your doctor before using this appliance.
- DO NOT use this device continuously for more than 45 minutes at a time. Extensive use could cause overheated on the device. Should this occur, discontinue use and allow the unit to cool down before using it again.

- If you experience pain in a muscle or joint for a prolonged period of time, discontinue use and consult your doctor.
- Never use the appliance while sleeping by sitting down or lying down.

Specifications:

Rated Voltage: DC 12V, 2.5A

Rated Power: 48W



Correct Disposal of this product
This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled

waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.