

Instructions

1

Quick Press 1x = Power on
Quick Press 2x = Power off

Power Button



1. Powering on the device

A) Quick press the power button to power on the device. The device will power on with low vibration mode running.


B) Quick press the power button a second time to power off the device.

2

Quick Press 1x = Vibration high
Quick Press 2x = Red light +Vibration low
Quick Press 3x = Red light +Vibration high
Quick Press 4x = Light+vibration off

Red light

Vibration Button

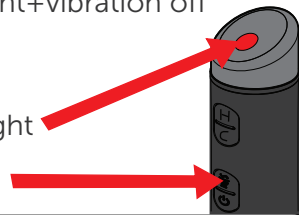

2. Vibration Red light therapy modes

A) Quick press power button to power on device.

B) Quick press vibration button 1x to set the vibration to high.

C) Quick press vibration button 2x to power on the red light and set the vibration to low.

D) Quick press vibration button 3x to power on the red light and set the vibration to high.

E) Quick press vibration button 4x to power off the red light and vibration mode.

3

Quick Press C 1x = Cool Mode
Quick Press H 1x = Warm Mode
Quick Press power button 2x = Reset



3. Warm/Cool modes


A) Quick press power button to power on device.

B) Quick press the C button to set the device to cool setting. The LED will light solid blue indicating the cool feature is active.

C) Quick press the H button to set the device to warm setting. The LED will be solid red indicating the warm feature is active.

D) To power off the warm and cool settings, quick press the power button 2 times.

4

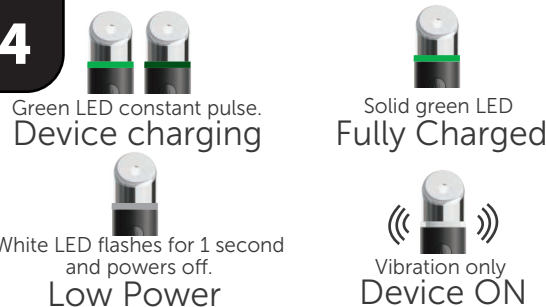

4. Charging the Unit

Charge unit using USB-C port and included charging cable as shown in image. Charge using any 5V USB power source. During charging, the LED light will pulse green to indicate charging. Once the LED turns solid green, this indicates that the device is fully charged.

Note: If the white LED is blinking while powering on, it's time to recharge the device. The device will not run when power is low.

 **BlackWolf**®

revitalizing eye pen



INSTRUCTIONS

©2024 Allstar Innovations

Black Wolf is a registered trademark of Black Wolf Opco, Inc.

Made in China | Patent Pending
Product Color and Appearance May Vary

Elevate your shower experience.
Visit [Blackwolfnation.com](https://blackwolfnation.com) for all
your skin and haircare needs.

IM2

BC

Components



1x Massaging eye pen



1x USB-C Charging Cable



1x Cap

Important safety instructions.
Please read all instructions before use.
Use this product only as indicated in these instructions.

WARNING:

- The device is only intended for domestic/private use, and not for commercial use.
- The device is only intended for the purpose described in this instruction manual.
- To Reduce risk of electric shock, DO NOT get the product, or charging cable wet or immerse in water while charging.
- If you are under a doctors care for any skin condition or if you think you may have a skin condition, please consult your physician before use.
- Do not clean device with harsh chemical agents.
- Do not use the device on broken skin, or open wounds.
- Do not use on animals.
- If the product or its components are damaged, Please discard immediately. Do not attempt to repair, or disassemble the device.
- Avoid contact with eyes.
- Do not stare directly into the red light LED.
- Discontinue use if irritation or rash develops.

CHOKING HAZARD:

Keep out of reach of children. Do not allow children to play with this device or its components. This is not a toy for children.

FCC Certified

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

NOTE: This product has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

FAQ's

How should I clean the device?

The device can be cleaned using an antibacterial wipe. make sure to put the cap back on after use

The power light started to blink red.

The white LED blinking indicates that its time to recharge your device. Plug the provided USB-C charging cable and charge until LED turns solid green.

How do i turn on the red light LED?

Quick press the power button to power on device. Quick press the vibration speed two times to activate the normal vibration and red light LED. Quick press three times to activate fast vibration and red light LED.

How do i make the device turn warm/cool?

Quick press the power button to power on the device. Press the H button to power on the warm feature, and the C button to power on the cool feature. The LED will light up Blue for cool and red for warm when the feature is running.

Getting Started

Read all instructions before first use.

• Charge unit using USB-C cable. When the light on the device is solid green, the device is fully charged and ready to use.

• Apply eye gel or moisturizer to the problem area. For Optimal results we recommend using BlackWolf Eye gel and BlackWolf Moisuizer.

Blackwolfnation.com

• Power on the device by quick pressing the power button. Quick press the vibration speed button once for higher vibration speed. quick press again for regular speed with red light LED. Quick press the H or C button for warm or cool treatment. Quick press the power button to power the device off.

• After use, turn off the device and wipe down the device with a cleaning wipe and reapply the cap.



Product name: Revitalizing Eye Pen

Model No.:

Brand name: BlackWolf

Power input: USB DC 5V 1A