




PL037
EMS Red Blue Light 4
Probes Beauty Massage Device
Operating Manual

Introduction

This is a comprehensive personal skin care EMS Eye Device that combines EMS micro-current and led skin care.

When used, it can be combined with skin care products such as cream or essence to achieve a scientific and efficient skin-beautifying effect.

Part Identification

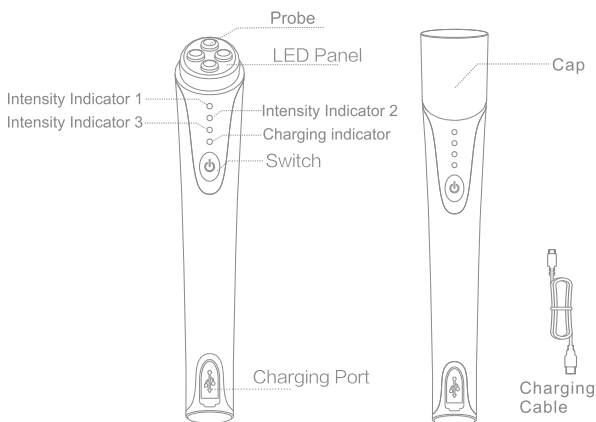

- ※ Turn on/off: Press and turn on the power button for 1 second .
- ※ The machine has a 10-minute timer shutdown function.
- ※ Intensity selection: The default is 1 for the power on. Each time the Intensity selection button is pressed, the 1 position will be changed accordingly, for a total of 3 Intensity.
- ※ For the first time, it is recommended to use 1 gear.

Thank you for purchasing this product! In order to use this product correctly, please read this manual before use; in order to read the manual at any time when you need it, please keep this manual in a safe place.

About this manual:

The product pictures and related product features shown in this manual, as well as some nouns, can only be used as a reference. The user will use the actual product as the standard, and the referenced data comes from the company's laboratory.

Safety Precautions

Be sure to read the safety precautions before use.
Please use this product correctly according to the contents and recommendations of the manual.
The product should be cleaned in time after use to avoid accumulation of grease and its residues.

	Related content: Incorrect actions can lead to personal injury
	Related content: Incorrect operation can result in personal safety or property damage

Danger

- Please use the charging cable supplied with this product or the charging cable that meets the product specifications, and connect the charger plug to the household socket correctly to avoid malfunction or fire caused by heat.
- Since this product has a built-in rechargeable battery, do not put it in a fire or use it in a high temperature environment, which may cause malfunction or fire due to heat.
- If the product is unusually hot or deformed, damaged, or smelly, do not use it to avoid the risk of electric shock and fire.

Caution


- To prevent the occurrence of danger, do not use it with medical electronic equipment.
- It is strictly forbidden to use a person who has undergone cardiac surgery or a medical electronic device such as a pacemaker or an embedded defibrillator placed inside the body.
- For first-time use or sensitive skin, please try a small area of skin and confirm the condition of the skin before using.

How to Use

【 Charging 】

For the first time, please fully charge the machine before using it.

- Connect the machine to the charging cable and charger, and plug into the household power socket;
- During charging, the charging indicator flashes red;
- When fully charged, the charging indicator lights up in green;
- Charging time takes about 2.5 hours;

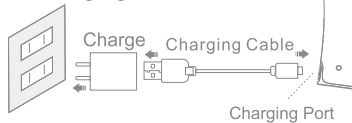

【 Preparation before Use 】

- Please remove the metal objects (rings, necklaces, earrings, bracelets, etc.) from your body.
- Please remove makeup and cleanse before use to ensure that your skin is in good condition.
- When taking care around the eyes, remove the contact lenses first.



Precious metal removal



Remove makeup and clean face before use

※ Please do not use this product in the larynx, eyeball, mucous membranes, etc.

◆ Lightening eye lines ◆

Apply a proper amount of eye cream or beauty gel. Turn on the device. Use the probe to fully touch the skin. The force is moderate. The probe is pulled along the vertical direction of the eyelids and slowly pulled about 1 cm, and stays here for about 3-5 seconds, until the eye lines are stretched and fastened. Repeat 5 times and do the next part.



- If the machine is damaged, in order to avoid danger, it must be repaired by the manufacturer, its maintenance department or a specialist in a similar department.
- Do not use this product (such as alcohol, benzine, spray, etc.) near flammable or explosive materials, which may result in fire or explosion.
- Place packaging materials such as plastic bags where children can't reach them.
- Please turn off the power before cleaning and maintenance.
- Do not use this product for purposes other than those specified in this manual to avoid malfunction or fire due to heat.
- When using, if the product is abnormally hot, please stop using it immediately and shut it down.
- Do not use the machine if it is faulty. Using this product in a faulty state may have adverse effects on the body or other hazards.
- Please do not use this product continuously in the same area, otherwise it may cause a burning sensation.
- This product is not intended for use by persons with physical, sensory or mental problems, or those who lack experience and knowledge (including children), unless supervised by a guardian and instructed in the use of the product, and responsible for the safety of the user.
- Do not immerse the appliance in water or rinse it directly with a water tap.
- If you feel any discomfort in your skin or body, stop using it immediately.
- Please do not store this product in a humid environment.
- When using, remove precious metal objects and watches.
- Do not use this product in the following situations:
 - During pregnancy, lactating or menstruating women.
 - Allergic people have skin diseases, skin inflammation, skin irritation, other skin diseases and severe sunburn.
 - People with skin damage, wetness, swelling, etc.
 - The person who has just finished surgery and is recovering.
 - Due to the upper fracture or the plastic face, the metal has been received on the face or in the body, and the silicon polymer is embedded in the surgery.
 - People with heart disease, malignant tumors, epilepsy, diabetes, infectious diseases, acute conditions.
 - Long-term use of steroid hormone drugs or liver dysfunction caused by capillary vasodilatation.
 - In the past 3 months, surgery has been performed at the site of use, with long-lasting laser or deep chemical stripping.

◆ Lifting and shaping of eye corners ◆

Apply a proper amount of eye cream or beauty gel. Turn on the device. Hold the instrument in one hand and allow the probe to fully contact and press the outer corner of the eye 1 cm below the center. The probe is lifted from bottom to top, slowly pulling about 2 - 3 cm, and staying here for 2-3 seconds. Repeat 5-6 times to help tighten the corner of the eyes. Keep your movements comfortable during use.



◆ Improving of eye circles and pouches ◆

Apply a proper amount of eye cream or beauty gel. Turn on the device. Holding the machine with your hand. Gently massage the dark circles and eye bags with the probe. The force is moderate, the movement is soothing, and the principle of moving from the inner corner to the outer corner of the eye is followed. Repeat 8-10 times as appropriate.



◆ Introduction of eye essence ◆

Apply a proper amount of eye cream or beauty gel. Turn on the device, gently massage the skin around the eyes. The action is soothing, following the principle of moving from the inside out and from the bottom to the top. Repeat this 5-6 times.



【 Power Off 】

After the care is completed, press the power button to turn off the phone.

- ※Tips
- When the care was first started, the skin did not feel the obvious heat. After a few minutes, began to feeling it from the inside out slowly.
 - During the nursing process, the skin may feel hot or slight pain. Some users may feel a slight redness after using the skin. This phenomenon usually disappears within a few minutes after use. If the skin feels abnormally hot during use, it is recommended to stop using and let the skin cool down first.

Fault analysis and elimination

In case of any malfunction during the operation, please carry out following inspections. If failures cannot be eliminated or the product cannot operate properly, please contact Kingdom's after-sales service center. Do not disassemble and attempt to repair it by yourself.

Situation	Possible Reason	Measure
Sudden shutdown or machine does not start	The battery is not fully charged	Fully charge the host and use it again
Charging light is not lit	The main unit, charging stand, and charger are not connected	Reconnect the main unit and charger, or plug the charger into another household socket.
	Battery is dormant	The battery needs 6-30 minutes back to working.It is a normal phenomenon and needs to be fully charged before use.
Massage head contact with skin without irritation or feeling of heat.	The skin is too dry	Please apply the lotion to your skin before use

Maintenance

- Protect the appliance from moisture and do not place the appliance in a dishwasher or microwave oven.
- This product is not waterproof and should not be immersed in water.
- Do not use sharp or corrosive, volatile thinner cleaners.
- It can be cleaned with a soft damp cloth and wiped with a soft, dry cloth.

After-sales

- This appliance is intended for personal use only and is not suitable for commercial use.
- We offer a one-year warranty on the product, which is calculated from the date of purchase. During the warranty period, we will replace parts or the entire product due to the material or the production process.
- To apply for a warranty, the product must be shipped to an authorized Prospera repair station or factory along with a certificate (purchase invoice) that proves the date of purchase.
- The company does not guarantee the following conditions:
 - Damage caused by improper use, such as using the wrong power supply, inserting the product into an improper socket, man-made damage, etc.
 - The product caused by normal use is aging and wear, but it has little effect on the function of the product.
 - Product problems caused by other uses, or improper handling or non-professional repairs or the use of accessories other than the company .
 - Consumables or bonus accessories.
 - Damage caused by force majeure factors (such as lightning, earthquakes, fires, floods, etc.).

Environmental Protection

The symbol on this product or in the instruction manual or on the packaging indicates that this product cannot be disposed of as general household waste at the end of its useful life. Please consult your local government for recycling according to the classification of waste.




Waste Treatment

The marking on this product or packaging material or accessory indicates that the material is made from a recyclable material. Separate waste into different recycling levels to facilitate the recycling of useful materials.



- There is a lithium battery inside the appliance, therefore:
- The battery must be removed from the appliance before it;
 - When removing the battery, the appliance must be turned off;
 - The battery should be disposed of safely;
 - Users cannot replace the battery by themselves.

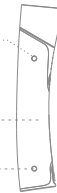

- Steps to remove the battery:
- Unscrew the screw plug of the back cover of the product and unscrew the screw;
 - Use the tool to open the back cover of the machine;
 - Pull out the battery lead and remove the battery.

Specification

Product Name: EMS Red Blue Light 4 Probes Beauty Massage Device
Model: PL037 Rated Voltage: 3.7V ---
Rated Power: 5.0W Charging Voltage: 5.0V ---
Charging Current: 1A Net Weight: about 70g
Product Dimensions: 30x41x172mm
Packing List: 1 main unit, 1 charging cable,1 operation manual.



Warranty exclusions:

This warranty does not apply to any costs, repairs, orservices for the following:

- 1.Repairs necessitated by use other than normal home use.
- 2.Damaged resulting from misuse, abuse, accidents, alternation or improper installation.
- 3.Corrective work necessitated by repairs made by any-one other than Prospera authorized service technicians.

Non-warranty service is provided on a“per incident” basis. Consumer will verify that the product has failed and provide information for servicing the unit. All app-licable repairs, parts, shipping, handling, local tax and“per incident” fees will be charged for non-warranty rep-airs and support calls.

For warranty service, please call 925-292-5978,M-F8:00am to 4:00pm, Pacific Standard Time. Or send email to sales@prosperacorp.com.

Consumer may also mail service request to this address:
392 West Larch Road # 31Tracy, CA 95304-1644

About Prospera

Prospera develops high quality, innovative, and affordable massage products, ranging anywherefrom electronic pulse massagers to jade massagers. Prospera's products are available through leading department stores. consumer electronics stores, and major online retailersL earn more about Prospera at www.prosperacorp.com

Prospera Corporation

392 West Larch Road # 31Tracy, CA 95304-1644
Tel:(925) 292 5978
Email:sales@prosperacorp .com



Warranty card



Warranty Registration Form

Product Name: _____

Product Model: _____ Purchase Price: _____

Date of Purchase: _____ Retailer: _____

First Name: _____ Last Name: _____

Address: _____

Cit: _____State: _____ Zip: _____

Phone: _____ Email: _____

Signature and date: _____

Please mail to:
392 West Larch Road # 31Tracy, CA 95304-1644
Tel:(925) 292 5978
Email:sales@prosperacorp.com

