



Hands Free EMS Neck Massager PL025A

Operating manual



(Real product image)

Introduction

Prospera Hands Free EMS Neck Massager Heated 4 modes 15 levels is with stainless steel electrode sheets and PI heating sheets, micro-computer chip control technology, and helps relax the neck muscles and relieve neck fatigue through micro-current stimulation and hot compress. Hot compress can effectively raise the temperature of the environment around the neck, promote blood circulation, and cooperate with the micro current to drive the muscles to exercise together, giving people a full range of comfortable experience. And it is foldable and mini size for easy carrying.

Product Specifications

- Mini size and foldable
- Cordless
- 4 Massage modes
- 15 intensities options
- 3 temp options 38-43°C/ 100.4 -109.4 °F
- 15 min auto shut off
- Item weight: 55g /0.12 lbs
- Control unit: 1.4” x 1.8” x 0.6”
- Rope length: 450mm
- Box: 0.6” x 0.5” x 0.12”
- USB Charging
- Input: 5V 1A
- Rated Power: 5W
- Voltage: 3.7V, 300mAh
- UPC 850011335772

Content

CONTENTS


Safety Warnings.....	4
Product Illustration.....	4
Operating the Product.....	5
Consumer Warranty.....	6
Warranty Card.....	6

Safety Warnings

Consult your doctor before using this product if you have any lesions, skin conditions, heart issues, pacemakers, or other health-related concerns.

1. This product is intended for adult use only.
2. Do not let babies and children play with the product.
3. Do not use the product on an infant.
4. If you experience discomfort, stop using the product immediately.
5. Do not wash the unit using organic solvent such as thinner to prevent damage.

Product Illustration

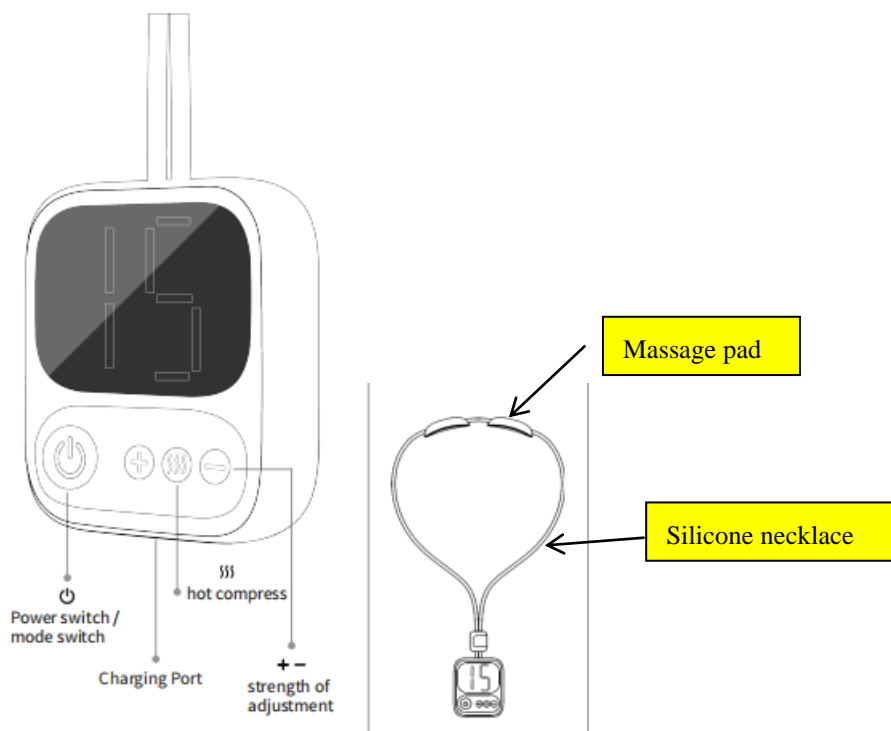

Operating the Product

Step 1: Preparation

Use a wet towel, wet paper towel or clean water to wipe and moisten your neck skin.

Step 2: Wear the massager

Wear the massager to your neck with the conducting plate facing your skin, and adjust the position.

Step 3: Starting up & Shutdown

Press and hold the host button of the product for two seconds, the host will turn on (hear a "di" sound, the blue light is on), enter the default mod (intensity 0) in the power-on state, press and hold the host button for two seconds, you will hear "di" two beeps, the product shuts down.

Step 4: Massage care

Press the switch button to choose preferred mode (four modes optional), press the heating button turn the heating on/off (Three temp. 38-43°C optional), and +/- to adjust the strength (15 intensities options).

Step 5: Device maintenance

Clean the conductive pads with a wet wipe or cloth.

Recommend usage: daily

To order a new health and beauty product, please visit www.prosperacorp.com.

Use coupon code: PROS\$5 for \$5 off any product above \$30.

Consumer Warranty

Prospera offers a 90-day standard warranty on massage products. Labor and parts are included during the warranty period. Please provide the model number and dated proof of purchase (sales receipt) when contacting Prospera for warranty service.

Non-warranty service is provided on a “per incident” basis. The consumer will verify that the product has failed and provide information for servicing the unit. All applicable repairs, parts, shipping, handling, local tax and “per incident” fees will be charged for non-warranty repairs and support calls.

Warranty Card



Warranty Registration Form

Product Name: _____

Product Model: _____ Purchase Price: _____

Date of Purchase: _____ Retailer: _____

First Name: _____ Last Name: _____

Address: _____

City: _____ State: _____ Zip: _____

Phone: _____ Email: _____

Signature and date: _____

Please mail to:

Prospera Corporation

392 West Larch Rd Suite 31

Tracy, CA 95304-1644

Tel:(925)292-5978

Email: sales@prosperacorp.com

