

DOUBLE CHIN REDUCING MASSAGER

Item No. 210199

User Guide

Thank you for purchasing the Double Chin Reducing Massager. Please read this manual thoroughly before operating this device.

INTRODUCTION

The Double Chin Reducing Massager uses electrical muscle stimulation (EMS) and LED light therapy to tone and sculpt your face for a healthier, more radiant look. Simply wear it while you work or relax and use the wireless remote to adjust the modes and intensity.

FEATURES


- Easy-to-use facial toning device
- Cordless and rechargeable
- Electrical muscle stimulation (EMS)
- Blue LED light therapy
- Red LED light therapy
- Rechargeable wireless remote
- Adjustable head strap
- 15-minute shutoff timer

PACKAGE CONTENTS

- Double Chin Reducing Massager
- Reusable gel pads (2)
- Wireless remote control
- USB charging cable
- Adjustable head strap
- Instruction Manual

IDENTIFICATION OF PARTS

1. USB charging port
2. Light emitting area
3. Wireless remote
4. Power On/Off
5. Mode button
6. Charging indicator
7. Intensity Setting
8. Remove Control Display

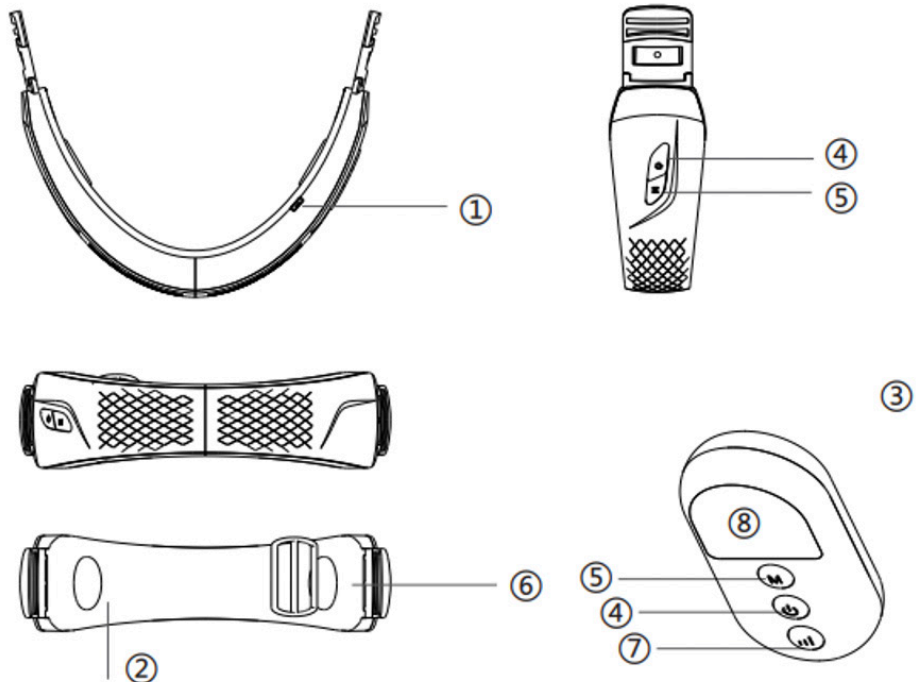

QUICK START

1. Remove Protective Covers:

- Remove the clear protective cover sheets from the round metal electrodes.

2. Prepare Gel Pads:

- Unwrap and place the reusable gel pads over the electrodes.
- After placing the gel pads, remove the blue protective film from each gel pad.
- **Note:** *Do not discard the blue protective films, as they can be used to cover the gel pads when not in use. This helps extend the life and conductivity of the gel pads.*

3. Attach the Head Strap:

- Attach the adjustable head strap. Each buckle should click easily into place.
- If it does not click into place, flip the buckle over and try the other side. Do not force it into place.

4. Position the Massager:

- Place the Double Chin Reducing Massager over your head and adjust the strap for a snug but comfortable fit. The "V" hinge should be on your neck, just behind your chin.

5. Turn On the Device:

- Long press (3 seconds) the power button on the Double Chin Reducing Massager to turn the device on.

6. Turn On the Remote:

- Long press (3 seconds) the power button on the remote to turn it on.

CHARGING THE DEVICE

1. Connect to Power:

- Connect the device to a USB power source using the included USB cable.

2. Charging Indicator:

- The charge indicator will light up red while the device is charging.

3. Charge Time:

- The charge time is approximately 3 hours.

4. Charge Completion:

- When charging is complete, the red LED will turn off. Disconnect the device from the charging cable.

CHARGING THE WIRELESS REMOTE

1. Connect to Power:


- Connect the remote to a USB power source using the included USB cable.

2. Battery Icon:

- The battery icon will flash while the built-in battery is charging.

3. Charge Completion:

- When the icon stops flashing (when the battery icon is a solid 4 bars), charging is complete. Disconnect the remote from the charging cable.

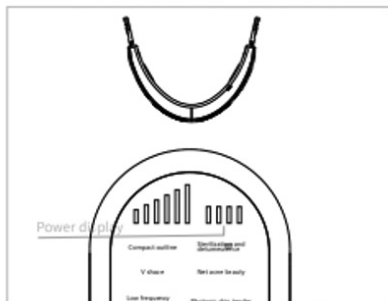

CLEANING AND MAINTENANCE

1. Wipe Clean:

- Wipe clean with a soft cloth if needed.

2. Avoid Water:

- Do not immerse in water or any other liquid.

3. Avoid Harsh Chemicals:

- Do not use any harsh chemicals or cleaning solutions on this device.

4. Storage:

- When not in use, store this device in its original carton in a cool, dry place.

INSTRUCTIONS FOR USE

Step 1: Using the Double Chin Reducing Massager

1. Preparation:

- Ensure the device is charged and the wireless remote-control is functioning properly.

2. Application:

- Place the device on your face.
- To adjust the strap, use the buckle located at the bottom of the device [near the chin joint] for optimal contact with the skin.




Step 2: Operating the Device

1. Power On:

- Press and hold the power key on the remote control for 1 second to turn on the device.

2. Mode Selection:

- The beauty meter can be used independently or in conjunction with the remote control.
- For independent use, press the power key for 1 second.
- To use with the remote control, follow these steps:
 - Press the remote-control key.
 - The device will emit a beep indicating it has turned on.
 - To shut down, press the power key for 2 seconds.

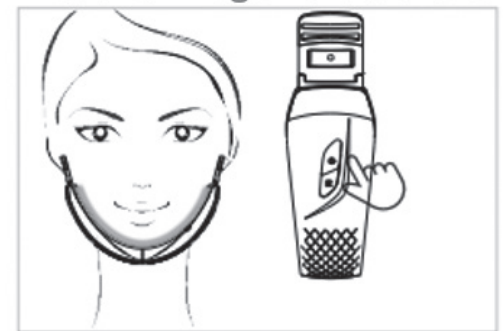

3. Mode Cycle:

- Use the remote-control button labeled "M" to select and cycle through different modes.

4. Available Modes: The Instrument has six modes:

To switch to the next mode, press the "M" button on the remote control

1. Compact Contour: For a firm and defined facial contour.

2. V Plastic Face: To sculpt and enhance the V-line of the face.

3. Low-Frequency Magnetic Therapy: For relaxation and improved skin elasticity.

4. Sterilization and Deswelling: To reduce swelling and sterilize the skin.

5. Clean Acne Beauty: For acne treatment and prevention.

6. Photon Skin: Uses light therapy for skin rejuvenation.



Step 3: Adjusting the Strength:

1. Adjusting the Strength:

- Use the gear strength adjustment button (I/II).
- The default strength is none. Each press increases the strength by one level.
- When you reach the maximum strength, press the button again to return to none.
- To cancel all current strengths, press and hold the strength key.
- Press the button again to return to low strength mode.



FUNCTION DESCRIPTION

Mode Functions

Mode	Name on Remote	Description
1	Compact outline	EMS + Blue LED for contour tightening
2	V shape	EMS for face shaping
3	Low frequency magnetic therapy	EMS + Red LED for lifting and brightening
4	Sterilization and detumescence	EMS + Red LED + gentle heat for swelling reduction
5	Net acne beauty	Red LED for acne cleansing and repair
6	Photonic skin tender	EMS + Blue LED + micro vibration for photo rejuvenation

NOTE: There is a 15-minute automatic shutoff timer. You may restart the device 2 or 3 times for additional treatment if needed. We recommend using this device 2 or 3 times per week for best results.

DETAILED MODE DESCRIPTIONS

1. Daily Maintenance - Red and Blue Intermittent Alternation and Vibration Mode:

- The device will irradiate 3 times with Red Light from left to right and right to left.
- Then, it will switch to Blue Light irradiation 3 times.
- This cyclical performance of Red Light and Blue Light irradiation, combined with vibration therapy, is used for daily skin care.

2. Emergency Skin Repair - Red and Blue Continuous Alternation and Vibration Mode:

- The Red Light and Blue Light continuously circulate from the middle to both sides to switch the irradiation and vibration therapy.
- This mode is used for emergency repair care, providing quick and synchronized skin repair with Red Light and Blue Light.

3. Intelligent Skin Blue Light Flashing Mode:

- The Blue Light flashes slowly on the skin without vibration.
- This mode is used for relieving over-fatigued skin, lightening the skin burden, and providing skin relief.

INTRODUCTION OF COLOR LIGHT THERAPY AND VIBRATION THERAPY

The Double Chin Reducing Massager has four major therapies: red light irradiation, blue light irradiation, EMS micro-current, and vibration. These can be applied alone or in combination to enhance skin elasticity and shape a V face.

THERAPIES

1. Redray [Red Light Irradiation]:

- **Benefits:** Improves the activity of skin cells, promotes metabolism, and stimulates collagen secretion and fibrous tissue formation.
- **Effects:** Enhances blood circulation, skin elasticity, anti-aging, anti-oxidation, and skin repair.



2. Blueray [Blue Light Irradiation]:

- **Benefits:** Suppresses common skin inflammation, calms skin, and regulates keratinocyte secretion.
- **Effects:** Sterilizes, protects skin, strengthens skin absorption, repairs, balances skin grease, whitens skin, reduces wrinkles, and pores.



3. EMS [Electro Muscle Stimulation]:

- **Benefits:** Releases intermittent pulse currents, stimulates ATP production, and activates subcutaneous tissue.
- **Effects:** Improves facial muscle fiber elasticity, activates collagen activity, tightens skin, and enhances overall skin appearance.



4. Vibration:

- **Benefits:** Relieves facial muscles, improves fine muscle strength, and increases facial elasticity.
- **Effects:** Fully pulls and tightens facial muscles, shaping the V face.



SAFETY WARNINGS

- Do not use this device for anything other than its intended purpose.
- Discontinue use if you experience discomfort or adverse reactions.
- This device is not a toy. Keep out of reach of children at all times.
- Do not attempt to open or repair this device. There are no user-serviceable parts.
- Do not use this device near water [sink, shower, tub, pool, etc.].
- Do not use this device around high heat or flames.
- Dispose of this device in accordance with local laws.

WARNING

Consult your physician before using any at-home massage device. If you have sensitive skin, we recommend you try it for 1 minute. If there are no adverse reactions, continue the regimen as usual.

BEFORE USING THIS DEVICE, CONSULT YOUR PHYSICIAN IF YOU:

- Have facial trauma
- Have severe skin inflammation
- Are recovering from surgery
- Have facial fractures
- Are sensitive to metal
- Have high/abnormal blood pressure
- Are under 12 years old
- Are pregnant
- Are recovering from cosmetic surgery
- Wear hearing aids
- Have eardrum inflammation

SPECIFICATIONS

Charging Voltage DC 3.7V~5.0V

Rated Power: 3.5W

Charging Time 3 hours

WARRANTY AND CUSTOMER SERVICE

This item from SharperImage.com includes a 1-year limited replacement warranty. If you have any questions not covered in this guide, please call our Customer Service department at 1 (877) 210-3449. Customer Service agents are available Monday through Friday, 9:00 a.m. to 6:00 p.m. ET.

SHARPER IMAGE®

The SHARPER IMAGE® name and logo are registered trademarks.

©Sharper Image. All Rights Reserved.