
RECHARGEABLE CALLUS VACUUM

Item No. 210203

User Guide

Thank you for purchasing the Rechargeable Callus Vacuum. Please read this manual thoroughly before operating this device.

DESCRIPTION OF PARTS

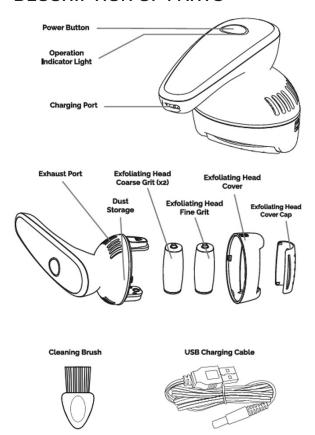

IMPORTANT SAFEGUARDS

Read All Instructions Before Using. When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

Danger

To reduce the risk of electrocution:

- · Do not use while bathing.
- Do not place or store the product where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquids.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Do not put the unit or power cord in water or other liquids.
- · This unit is not waterproof or water-resistant. Do not use near water.
- · Unless charging, always unplug from the outlet when not in use.
- Do not abuse or damage the power cord.
- · Unplug the unit before cleaning.
- Do not attempt to alter the plug or cord in any way.

Warning

To reduce the risk of burns, electrocution, fire, or injury to persons:

- The product should never be left unattended when plugged in.
- · Close supervision is necessary when this product is used by, on, or near children or invalids.
- · Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- · Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
- · Keep the cord away from heated surfaces.
- · Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
- · Never drop or insert any object into any opening.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- · Only use on clean, dry, lotion-free skin.
- · Ensure the plug and outlet are a safe distance from any water source.
- To clean the exterior of the unit, unplug it and use a damp cloth with optional mild detergent.
- If the device appears to be damaged, do not use it.
- · Any attempt to open the unit will void the warranty.
- · Always attach the plug to the unit, then to the USB port. To disconnect, turn off the unit and unplug from the outlet.
- · When unplugging, pull on the plug and not on the cord.
- Do not wrap the cord around the unit. Never handle the unit or power cord with wet hands.
- · Ensure there are no foreign objects inside the unit.
- · Make sure the unit is properly assembled and all moving parts are locked into place before use.
- Do not use attachments not recommended by the manufacturer.
- · For hygienic reasons, this unit should only be used by one person.
- · Only use the device for its intended purpose.
- · Household use only. Do not use outdoors.
- · Keep out of reach of children.
- · Store in a cool, dry place. Do not leave in direct sunlight.

CAUTION: To reduce the risk of fire or electric shock, use indoors only.

CAUTION: High voltage. Keep out of reach of children.

OTHER PRECAUTIONS

When the unit is on, keep it away from hair on your scalp, eyebrows, eyelashes, clothes, brushes, wires, cords, shoelaces, or anything else that might become caught in the unit.

Consult a doctor before use if you have any of the following conditions:

- · Allergic or sensitive skin that reacts to cosmetics, metals, clothing, etc.
- · Acne or body acne
- · Warts
- Moles or loose skin in the area to be treated
- Niahetes
- · Raynaud's disease
- · Eczema
- Psoriasis
- · Hemophilia
- · Immune deficiency
- · Rosacea
- · Couperose
- · Varicose veins
- · Wounds
- Rash, sunburn, inflammation, or irritation
- · If pregnant or recently pregnant
- · If you have a pacemaker

OPERATION

The exfoliating head is designed for removing coarse, callused skin from feet and heels. Do not use anywhere else on the body. Do not use on broken or bleeding skin.

1. Preparation:

- Wash and completely dry the skin that is to be treated.
- Remove the cap and press the power button to turn on to low speed (1). The exfoliating head will start spinning, and the operation light will light up blue.

2. Exfoliation:

- Touch the exfoliating head to the area on your foot that you would like to treat.
- Press gently at first, then apply more pressure and/or increase to high speed (2) by pressing the power button again (the operation light will be solid red on high-speed mode).
- If you press too hard, the exfoliating head will stop spinning.

3. Movement:

- Gently move the unit back and forth or side to side for no more than 2-3 seconds at a time until the callused skin is removed.
- For a smoother appearance, switch to the fine grit exfoliating head and repeat these steps.


4. Turn Off:

- Press the power button a third time while the operation light is solid red to turn the device off.

NOTE: As soon as you feel the exfoliating head touch new, healthy skin, or when the callus has been removed, stop use.

CHANGING EXFOLIATING HEADS

- 1. Turn the product off and remove the exfoliating head cover by pulling off the exfoliating head.
- 2. Press the release button as shown in the diagram and pull the exfoliating head upward to remove.
- 3. Insert one side of the replacement exfoliating head into the square groove on the side of the device with the release button.
- 4. Press the opposite side of the exfoliating head into place in the opposite square groove. You will hear the new head click into place when installed.
- 5. Replace the exfoliating head cover.

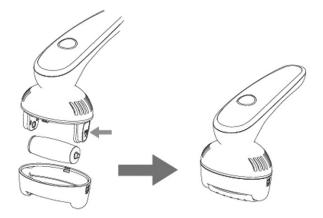

CHARGING

- Low Battery: When the battery is low, the operation light will blink red. If the battery is below operating values, the operation light will blink three times and the device will shut down.
- · Charging:
 - To charge the device, plug the charger into the charging port.
 - While charging, the operation light will glow red when the battery is low and turn green when fully charged.
 - Low power is indicated by a blinking red light in all modes.

CLEANING

- · Cleaning the Device:
 - Remove the exfoliating head cover.
 - Rinse the device under running water to clean.

By following these charging and cleaning instructions, you can maintain the proper functioning and hygiene of your device.

SPECIFICATIONS

Charging Voltage: 5V Charging Current: 1A

WARRANTY AND CUSTOMER SERVICE

This item from SharperImage.com includes a 1-year limited replacement warranty. If you have any questions not covered in this guide, please call our Customer Service department at 1 [877] 210-3449. Customer Service agents are available Monday through Friday, 9:00 a.m. to 6:00 p.m. ET.

SHARPER IMAGE®