

REJUVENATING HEATED COMPRESSION LEG WRAPS


Item No. 210337

User Guide

Thank you for purchasing the Rejuvenating Heated Compression Leg Wraps. Please read this manual thoroughly before operating this device.

DESCRIPTION OF PARTS

1. Leg Massagers
2. Air Pressure Hose
3. LCD Screen
4. Control Panel

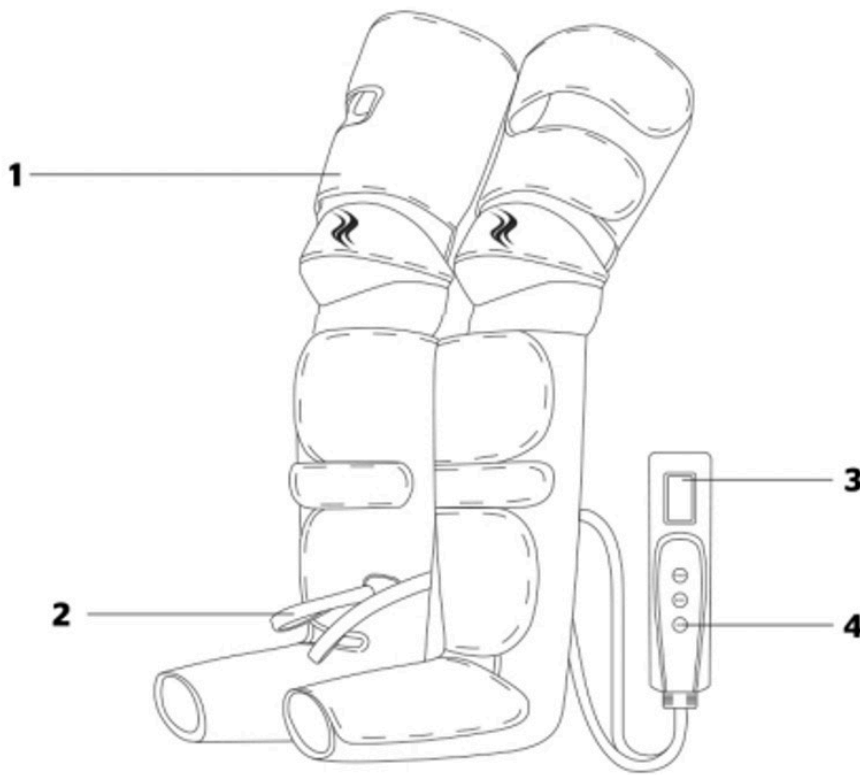

SAFETY PRECAUTIONS

- Do not use voltage beyond the requirement of this product.
- Do not attempt to repair or dismantle the product. There are no user-serviceable parts.
- Do not use the product in a place where it may become wet or damp. If the unit becomes wet, do not use it.
- Do not use this product in a high-temperature environment.
- Stop using the product immediately if you feel uncomfortable or if any abnormalities occur.
- Do not place the unit in a place where it may be in contact with flammable items, explosive items, or chemicals.
- Do not touch the unit with any sharp object, to avoid damage.

Do not use this product if any of the following apply to you:

- Your skin is infected, swollen, burnt, or otherwise damaged.
- You have a pacemaker or other medical implant.
- You are pregnant.
- You are suffering from heart disease.
- Your bones are fragile or injured.
- You are suffering from cancer.
- You are receiving treatment, or your blood pressure is not normal.
- You are suffering from osteoporosis.

By adhering to these safety precautions, you can ensure the safe and proper use of this product.

CONTROLS

1. Power Button [POWER]


- **Turn On:** Press and hold for 3 seconds to turn on the machine.
- **Adjust Intensity:** Press a second time to adjust the intensity to LOW, MIDDLE, or HIGH.
- **Turn Off:** Press and hold for 3 seconds again to turn off the product.

2. Mode Button [MODE]

- **Select Mode:** There are 3 modes to choose from:
 - **Mode 1:** Calf + Foot Airbag
 - **Mode 2:** Only Calf
 - **Mode 3:** Only Foot

3. Heat Button [HEAT]

- **Heat On/Off:** Press to turn the heat on or off. The heat is applied to the bottom of the knee.

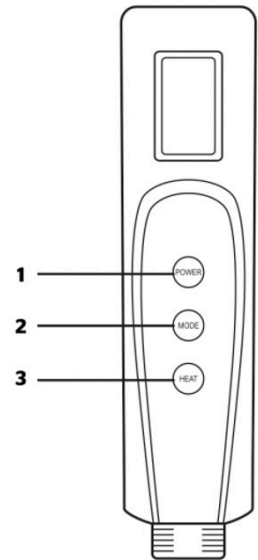

OPERATION

1. Prepare the Massager:

- Slip on each massager, keeping your feet flat on the bottom. Ensure each is marked L for left foot and R for right foot.

2. Secure Your Feet:

- Use the strap to tightly secure your foot into the massager.

3. Secure Your Legs:

- Use the strap to tightly secure your leg into the massager.

4. Check Fit:

- Move your foot to ensure the massager is securely fastened.

5. Connect Air Pressure:

- Insert the two air pressure connections into the controller.

6. Power On:

- Press the power button. The LED will stay on when the unit is powered.

7. Adjust Settings:

- Adjust the air pressure strength and mode by hand. You can see your operation on the LCD screen.

Note: The device will automatically shut off after 20 minutes. For best results do not use the leg massager for more than an hour.

CARE AND MAINTENANCE

- **Storage:** Store the product in a dry place, away from moisture and direct sunlight.
- **Cleaning:** Wipe clean with a dry towel.

TROUBLESHOOTING

Problem	Cause	Solution
The device does not work.	The unit is not connected to a power source.	Make sure the unit is fully connected to the power.
	The power button was not pressed.	Press the power button.
The device suddenly stopped working.	The unit was running more than 15 minutes.	Restart the unit to resume use.

SPECIFICATIONS

- **Rated Voltage:** 12V 2A
- **Power Consumption:** 24W
- **Default Time:** 20 minutes

WARRANTY AND CUSTOMER SERVICE

This item from SharperImage.com includes a 1-year limited replacement warranty. If you have any questions not covered in this guide, please call our Customer Service department at 1 (877) 210-3449. Customer Service agents are available Monday through Friday, 9:00 a.m. to 6:00 p.m. ET.

SHARPER IMAGE®

The SHARPER IMAGE® name and logo are registered trademarks.
©Sharper Image. All Rights Reserved.