

CONGRATULATIONS ON YOUR REALTOUCH™ EYE MASSAGER.

Please read all instructions and warnings prior to use. Keep this sheet as it contains important information for future reference.



THE BLUETOOTH® WORD MARK AND LOGOS ARE REGISTERED TRADEMARKS OWNED BY BLUETOOTH SIG, INC.



SHARPER IMAGE® NAME AND LOGO (INCLUDING THE 77 LOGO, TOMORROW'S TOMORROW, AND REALTOUCH) ARE TRADEMARKS OR REGISTERED TRADEMARKS, DISTRIBUTED BY AND © 2022 MERCHSOURCE, LLC. IRVINE, CA 92618. ALL RIGHTS RESERVED.

FOR CUSTOMER SUPPORT, PLEASE CALL 1-800-374-2744.

Model # 1016957

PRINTED IN CHINA.



SHARPER IMAGE®

REALTOUCH™ EYE MASSAGER HOT + COLD THERAPY

SAFETY WARNINGS:

PLEASE READ BEFORE USE

- Not a children's item. Keep out of the reach of children.
- Do not use this item if you are pregnant or have an injury, infection, or illness.
- Consult your physician if you have a medical condition that may prevent safe use of this item.
- Test your comfort level and check massaged facial areas often.
- Individual reactions to massage, skin contact and temperature therapy may vary. Stop use immediately if skin becomes irritated or inflamed. If you feel any increase in physical discomfort or pain, or if you experience any loss of sensation or numbness. Consult your physician immediately if symptoms persist.
- Never wear the eye mask during charging.
- Do not fall asleep while using the eye mask.
- Always allow eye mask to cool if it becomes uncomfortably hot during use.
- Do not expose to extreme temperatures or moisture.
- Never touch the charging cable with wet hands.
- Do not carry or pull the eye mask by the charging cable.
- Do not use with attachments of any kind other than those that may be included.
- Thoroughly inspect products before every use. If the product appears worn, frayed, splintered, cracked, or broken in any manner, discontinue use and discard immediately.
- This product is not intended to treat, cure, or prevent medical conditions.

RECHARGEABLE BATTERY WARNINGS:

- This product contains batteries that are non-replaceable.
- There is a risk of fire and personal injury if battery is punctured, damaged, or misused.
- Do not dismantle, crush, incinerate or short-circuit battery.
- Never expose battery to extreme temperatures or moisture.
- Keep away from flammable materials and direct sunlight.
- Only charge using the provided charger, cable, and/or adapters.
- Allow battery to cool to room temperature before charging.
- Charge on a non-flammable surface and keep away from flammable products.
- Discontinue charging if battery becomes hot, swells, smokes, or gives off a strong odor.
- Battery chargers used with the product are to be regularly examined for damage to the cord, plug, enclosure, and other parts. In the event of such damage, chargers must not be used until the damage has been repaired.
- Recycle or dispose of battery according to federal, state, and local laws.
- Keep these instructions for future reference.

FCC STATEMENT


- This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:
1. This device may not cause harmful interference, and
 2. This device must accept any interference received, including interference that may cause undesired operation.
- NOTE:** This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
- Reorient or relocate the receiving antenna.
 - Increase the separation between the equipment and receiver.
 - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 - Consult the dealer or an experienced radio/TV technician for help.
- WARNING:** Changes or modifications to this massager not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

IMPORTANT: MAINTAINING BATTERY HEALTH

- Your massager is designed with a lithium-ion battery pack for cordless convenience, increased massage duration per charge, and frequent use and recharging.
- To maintain charge capacity and prevent battery damage, all lithium-ion batteries must be regularly used and recharged.
- For optimal battery health, keep the battery at approximately 50% state of charge when not in use.
- All batteries naturally lose charge when not in use. If you will not be using the massager for an extended period, check the massager at least every 30 days to ensure the battery is at approximately 50% state of charge. Recharge as needed to maintain a 50% state of charge but do not overcharge.
- **IMPORTANT: DO NOT** store the battery for an extended period with above 60% or below 40% state of charge. An incorrect state of charge during storage could irreparably damage the battery and may result in reduced charge capacity, inability to maintain charge, failure to recharge, or complete inoperability of the massager.
- For the health of the battery, always store the massager in a cool, dry place. The storage location should have a consistent temperature between 40°F and 70°F (5°C and 20°C) as well as a relative humidity at 70% or lower.

WHAT YOU GET

- a. RealTouch™ Eye Massager
- b. USB-C to USB-A charging cable

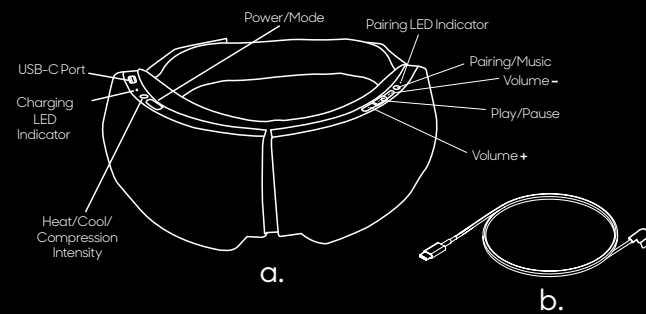


GETTING STARTED

CHARGING

1. Plug the USB-C side of the charging cable into the port located on the top of the eye mask.
2. Plug the USB-A side of the charging cable into a power adapter (not included).
3. Plug the power adapter into a wall outlet. The white LED on the eye mask will flash during charging.
4. Unplug the eye mask when it is fully charged. The LED will stop flashing and remain lit. It may take up to 3 hours to fully charge the eye mask.


Note: The eye mask cannot be used while charging. When the battery is low, a voice prompt will state "low battery" and a red LED will begin flashing. Recharge the eye mask to continue use.

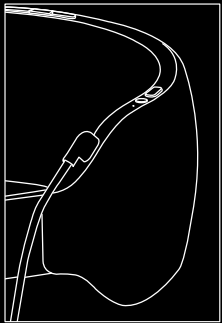

HOW TO USE

COMFORT FIT

1. Adjust the strap at both ends until you can easily and comfortably wear the eye mask. The strap features hook-and-loop fasteners that attach to the strap itself.
2. Make sure the eye mask sits comfortably on your face and that the fit is not too tight. Loosen the strap if necessary.

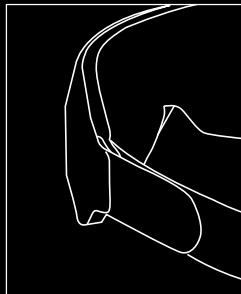


POWER

1. Press and hold the power/mode button for 2 seconds to power on the eye mask. The "power on" voice prompt will play.
2. Press and hold the power/mode button again to manually power off the eye mask. If left on, the eye mask will automatically power off after 15 minutes.

Note: For safety, the eye mask will automatically power off if the internal temperature rises too high. If the eye mask beeps as it powers off, allow it to return to room temperature before next use. Limit each massage or temperature therapy session to 15 minutes per use.


PAIRING

1. To connect your device via Bluetooth® wireless technology, press and hold the pairing button for 2 seconds. The eye mask will enter pairing mode and the LED will start flashing blue.
2. Select "SI MASSAGING EYE MASK" from your device's Bluetooth menu to connect the eye mask. The "Bluetooth connected" voice prompt will play and the LED will remain lit for 5 seconds.
3. To disconnect the eye mask and return to built-in soundscapes, press and hold the pairing button for 1 second.

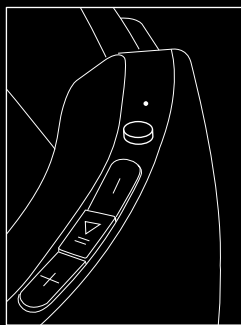


MUSIC SETTINGS

1. Press the play/pause button to play or pause audio.
2. Press the plus + button to increase the volume or the minus - button to decrease the volume. The eye mask powers on with audio at max volume.
3. To cycle through the built-in soundscapes, quick press the pairing/music button while the eye mask is not paired with another device.

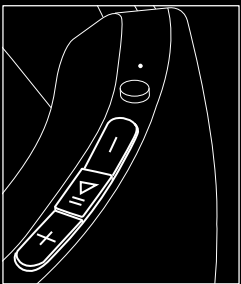

MASSAGE SETTINGS

1. Power on the eye mask to automatically activate air compression and vibration massage.
2. Quick press the power/mode button to cycle through the 3 massage modes or turn off the massage functions.
 - Mode 1: air compression only
 - Mode 2: air compression and vibration
 - Mode 3: vibration only
 - Off: disables air compression and vibration
3. To switch between low and high massage intensities, press and hold the hot + cold/intensity button for 1 second.

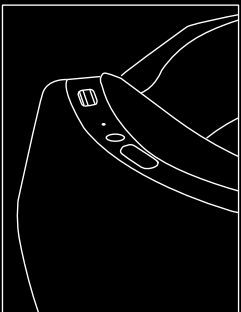

TEMPERATURE SETTINGS

1. Power on the eye mask to automatically activate alternating hot and cold therapy.
2. Quick press the hot + cold/intensity button to cycle through the 3 hot + cold modes or turn off the hot + cold functions.
 - Hot + Cold: alternates hot (3 minutes) and cold (2 minutes)
 - Cold: remains cold
 - Hot: remains hot
 - Off: disables hot and cold