



## **Prospera Penguin Red Light Body Massager PL032**

### **Operating manual**



### **Introduction**

---

The Prospera Penguin Red Light Body Massager is a compact, powerful device designed to relieve tension in the back, shoulders, legs, and feet. It features six adjustable speeds and comes with three interchangeable heads:

- **Roller head** for deep massage.
- **Brush head** for thorough body cleansing.
- **Soft brush** for gentle exfoliation.

With an IPX7 waterproof rating, this massager is suitable for a variety of environments, including the shower or bath.

## **Product Specifications**

- **Applicable Areas:** Suitable for the entire body.
- **Functions:** Body massaging and cleansing.
- **Attachments:** Three interchangeable massage heads.
- **Power:** Cordless design for portability.
- **Modes:** Two rotation massage modes.
- **Intensity Levels:** Six adjustable speed settings.
- **Safety Features:** Auto shut-off after 15 minutes.
- **Charging Time:** Approximately 4 hours with a 1.5V charger.
- **Usage Time:** Up to 2 hours per charge.
- **Rotation Speed:** 2500 RPM  $\pm$  200.
- **Weight:** 503 g / 1.11 lbs.
- **Dimensions:**
  - Main unit: 6.18" x 4.06" x 3.31"
  - Box: 8.23" x 5.71" x 3.9"
- **Charging:** USB powered (Input: 5V 1.5A).
- **Battery Capacity:** 1800 mAh.
- **Input:** 5V 1.5A
- **Rated Power:** 16W
- **Voltage:** 7.4V
- **UPC** 850011335796

## Content

---

### CONTENTS

Safety Warnings.....	4
Product Illustration.....	4
Operating the Product.....	5
Consumer Warranty.....	6
Warranty Card.....	6

## Safety Warnings


---

1. Consult your doctor before use if you have skin conditions, heart issues, pacemakers, or other health concerns.
2. This product is for adult use only.
3. Keep out of reach of children and infants.
4. Stop use immediately if discomfort occurs.
5. Do not clean the device with solvents such as thinner, as this may cause damage.

## Product Illustration


---

(下面的刀线图需要用高清可印刷的版本，下面只是截图示意)  
(刀线是原说明书里面的，标红部分的文字有改动)

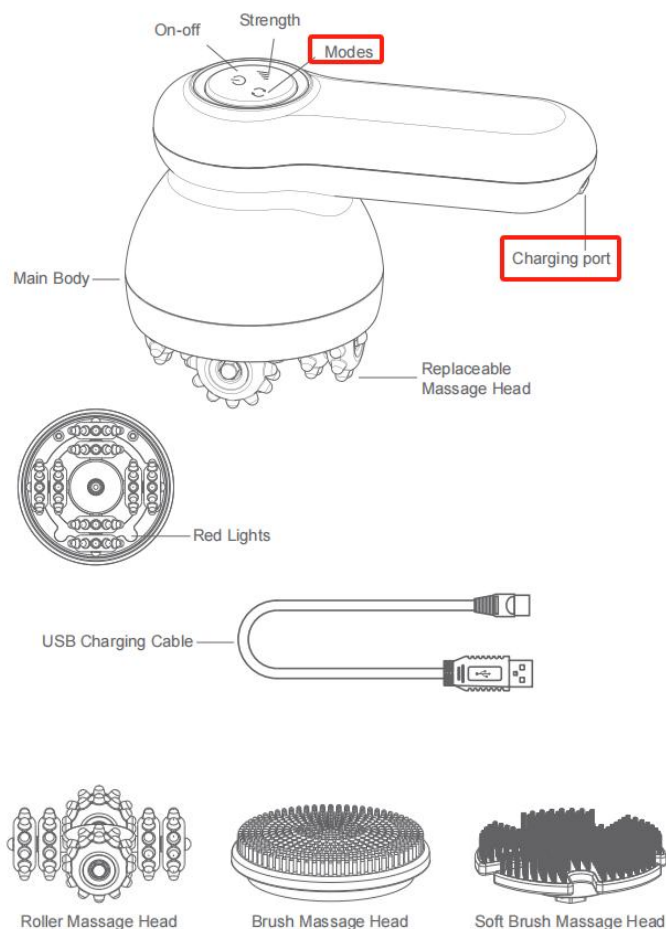

## Operating the Product

---

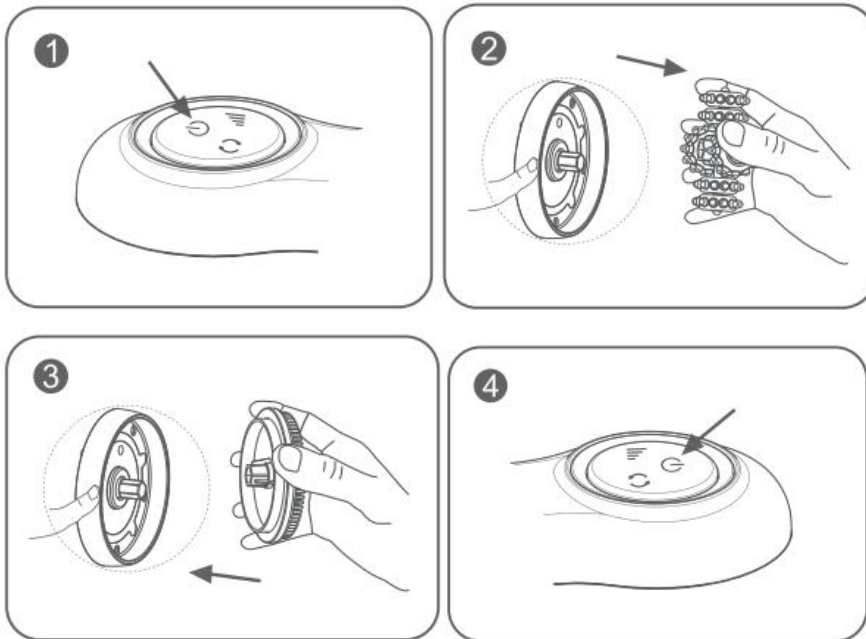

### Step 1: Preparation

Ensure the device is powered off before changing the massage head.

### Step 2: Changing the Massage Head

1. Hold the device securely.
2. Pinch the edge of the massage head and pull outward with firm pressure.
3. Attach the desired massage head and ensure it is securely connected.

### Step 3: Turning on the Device

1. Press the power button to turn the device on.
2. Select your preferred rotation direction.
3. Adjust the intensity to your comfort level.

### Step 4: After Care

1. Press the power button to turn the device off.
2. Clean the massage head thoroughly to remove residue.
3. **Recommended Usage:** Daily

To order a new health and beauty product, please visit [www.prosperacorp.com](http://www.prosperacorp.com).

Use coupon code: PROS\$5 for \$5 off any product above \$30.

## Consumer Warranty

---

### Standard Warranty:

Prospera offers a 90-day warranty covering labor and parts. Please provide the model number and proof of purchase when requesting service.

### Non-Warranty Service:

For repairs outside the warranty period, fees for parts, labor, shipping, and taxes will apply.

## Warranty Card

---



### Warranty Registration Form

Product Name: \_\_\_\_\_

Product Model: \_\_\_\_\_ Purchase Price: \_\_\_\_\_

Date of Purchase: \_\_\_\_\_ Retailer: \_\_\_\_\_

First Name: \_\_\_\_\_ Last Name: \_\_\_\_\_

Address: \_\_\_\_\_

City: \_\_\_\_\_ State: \_\_\_\_\_ Zip Code: \_\_\_\_\_

Phone: \_\_\_\_\_ Email: \_\_\_\_\_

Signature and date: \_\_\_\_\_

### Mail To:

Prospera Corporation  
392 West Larch Rd, Suite 31  
Tracy, CA 95304-1644

### Contact Information:

Phone: (925) 292-5978  
Email: [sales@prosperacorp.com](mailto:sales@prosperacorp.com)

