LIFE BALANCE
SMARTSCALE

PÈSE-PERSONNE INTELLIGENT

#535220, 535221


**READ THIS MANUAL CAREFULLY
BEFORE YOU BEGIN USING THIS PRODUCT**

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Scale will not be operational until these steps are completed. Smart phone required for operation.

STEP 1

Slide battery compartment tab (bottom of the scale) in the direction of the arrow and lift to release. Remove the insulating strip to allow installed batteries - power to connect to the scale battery terminals.

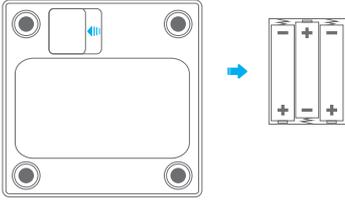Future battery replacement: Remove expired batteries and install 3 x AAA new, high quality batteries, carefully note polarity when installing. Replace cover.

STEP 2

Download the APP and complete Set Up.

Download the app "**MovingLife**" on iOS (iPhone®) or Google Play™ store (free to use, no in app additional costs).


App is available in over 25 languages
– Click on the "Me" icon, then go to settings and select "language".
Integrates with the Apple Health app.

Important: Keep updates turned ON for the APP for the most up to date version and optimal use.

STEP 3

Register and Connect

- Open the Life Balance Smart Scale APP and create your personal profile.
- Before connecting, ensure your Bluetooth® is turned ON – on your smartphone and the scale LED light is ON.

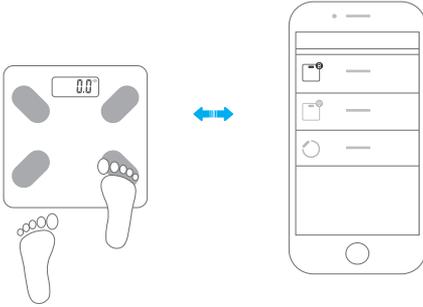Important: The first time you log in to the app you will need to click on your platform app (iOS (iPhone®) or Google Play™ store) and add your personal information. Occasionally check your personal information on the Me page to ensure accuracy.


- Complete the connection from the App to the scale, choose Smart Scale and click Bind when prompted.
- The APP itself is user friendly, if you need help, when in the home page – click on the "... " icon (top right) and then click on Help.

MOVINGLIFE APPLE HEALTH ACTIVATION

INSTRUCTIONS


1. Start with the movinglife App first launch wizard

1. Launch the movinglife App for the first time and log in.
2. After logging in, the App will automatically pop up a reminder to activate App Health. Click "Turn On All" to start data synchronization, and then click "Allow" to complete.

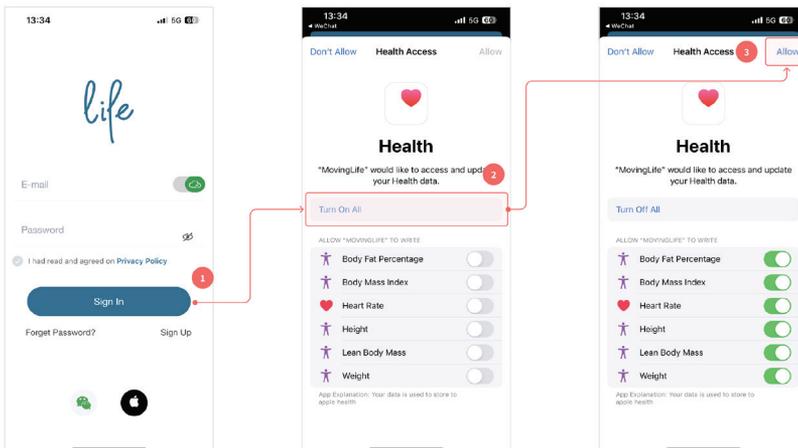II. Enter Apple Health and enable movinglife App Synchronization

1. Find the Apple Health app in the system and launch it.
2. After entering Apple Health, select the "Sharing" tab in the navigation bar. Then slide down the page, find the "Apps" option, and click to enter.
3. Find movinglife in the Apps list and click to enter.
4. Click "Turn On All" to start data synchronization, done.

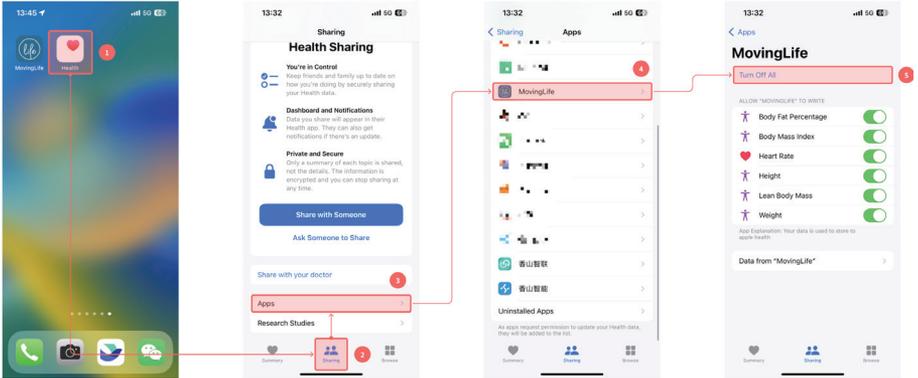STEP 4

Start Weighing


Place the scale on a hard, flat, dry surface like a floor. Stand on the scale with BARE FEET and step on the electrodes. Center your body and balance in the middle of the scale to avoid tipping. The weight reading is stable and flashes. The scale shuts OFF after several seconds. The APP will collect, calculate, show and save the data.

Under Body Composition, Weight, BMI etc, you can click on each measurement for additional information.

MULTIPLE USERS AND MULTIPLE DEVICES ARE SUPPORTED

On the Me page, click on User Management to add users to the App or click on My Devices to add a separate device.

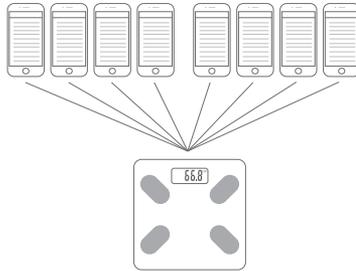WARNINGS, INDICATIONS ON LCD

ALWAYS ENSURE THE SCALE IS ON A FLAT, HARD, DRY SURFACE.

CAUTION: GLASS PLATFORM CAN BE SLIPPERY WHEN WET – ENSURE YOUR FEET ARE DRY BEFORE STEPPING ON THE SCALE.

Calibration – Error exists. and re-test.

Err – Overload indication or unstable measurement.

Lo – Low Battery Indicator – change the batteries.

IMPERIAL / METRIC MEASUREMENTS

LCD can be set to display in pounds (lb) or kilograms (kg). You can press the UNIT button on the bottom of the scale to set the default or switch.

ADVICE, USE AND CARE

This product is not a medical device and is not for use for medical treatment or intended for use by medical institutions. If you have any medical concerns about the use of this product, consult with a medical professional prior to use.

ABOUT THIS PRODUCT

This product is a smart body fat scale with a high-precision strain gauge sensor which can help you store and track weight and other measurements. The scale transmits information via Bluetooth® and charts your results. Device measurements and information provided are intended as a guide to help individuals measure, track and improve their health. Feedback provided is recommended advice for general health, but those with significant health concerns should always consult with a medical professional.

LOW BATTERY

After some time, the LED will stop working and you will not be able to connect to the App. Replace low batteries promptly to avoid battery leakage. Dispose of old batteries according to local regulations.

TEMPERATURE

Do not use the scale in humid or hot conditions. Operational temperature is 5°C – 35°C.

CARE

This is a precision instrument, do not drop and do handle with care. This product uses tempered safety glass to avoid shattered glass.

Keep the scale clean. Do not use acidic cleaning agents.

TROUBLESHOOTING

1. The APP can't obtain weight measurements.

- Confirm connection and signed into the APP and Bluetooth® is ON and within range.
- AAA Batteries may need to be replaced
- Relaunch the APP or Update the APP
- Turn Blue Bluetooth® OFF and then ON again, then bridge the connection again.

2. Fails to connect to your device

- Please ensure ONLY ONE DEVICE is setting up a connection
- Ensure your device is using Bluetooth® 4.0 or above

3. LCD does not light up

- Ensure you are on a flat hard dry surface like a floor
- Check Batteries

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1 YEAR WARRANTY

Relaxus Products Ltd. warrants this product to be free from manufacturer's defect for one year from the date of purchase. Original purchase receipt is required for any warranty claim.

Warranty excludes:

1. Damage caused by improper or incorrect installation or use.
2. Breakage or damage from mishandling or accidents.
3. Use of product for purposes for which it was not intended, or commercial use.

Relaxus Products Ltd. retains discretion to repair, replace parts or supply a new product of the same or similar model. Warranty is not transferable. Any associated shipping costs are not included in warranty coverage.

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