


140mm

LUNA'MOUR
5-Mode V-LIFT Ionic
Anti-Aging RF EMS & LED Light Therapy Facial Device



User Manual ▶

Please read the instruction carefully before use and store for future reference.

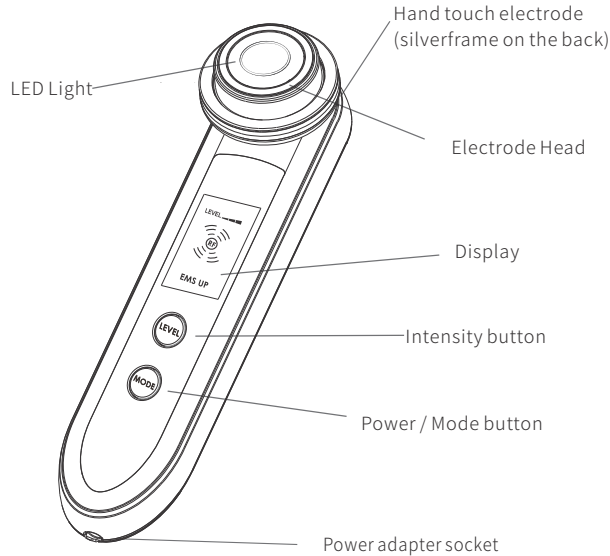

5-Mode V-LIFT Ionic Function Introduction▶

FEATURES ▶

Looking for a quick and easy way to improve the appearance of wrinkles and signs of aging at home?

The 5-Mode V-LIFT Ionic Facial Device with Anti-Aging RF, EMS & LED Light Therapy uses a combination of technologies including Radio Frequency (RF), Electrical Muscle Stimulation (EMS) and Light Therapy technology, as well as Heat Therapy, Cold Therapy, and Ionic Therapy.


These different types of technologies will leave your skin feeling refreshed and radiant!

- Radio Frequency (RF) helps to tighten skin, stimulate production of collagen, and brighten complexion.
- Electrical Muscle Stimulation (EMS) helps to stimulate muscle contractions in the face, which in turn, may help to strengthen and lift muscles for facial contouring.
- Light Therapy (Red | Blue) helps reduce wrinkles, fine lines, age spots, and acne-causing bacteria.
- Heat Therapy helps to increase the skin permeability, which results in better absorption of the skincare products to reach deeper into the layers of the skin.
- Cold Therapy helps to reduce puffiness and redness.
- Ionic Therapy – Positive ions helps to attract the encapsulated impurities from the skin (such as dirt, makeup, etc) while negative ions help to improve the skin's moisture barrier, allowing for a more hydrated complexion.

P1

OPERATION GUIDE ▶

● Display screen / Key

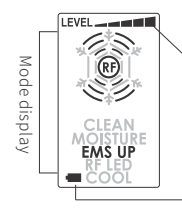


- Hold the power/mode button for 3seconds and turn on the power supply. (automatically enter RF cleaning mode)
- Each mode can be switched to other modes by pressing the power mode button.
- Adjustable intensity levels can be selected for the RF cleanse mode, RF moisturizing mode, RF EMS UP. Each mode has 5 intensities.
- RF temperature can not be adjusted.

Please note there is a 2-minute auto-shutoff for each mode. After that time period, the power/mode button will flicker, and device will be in standby state. Press the power / mode button to start the device again. If the device is not in use for two minutes, the power will turn off automatically (Auto-Off) after "buzzing" two times.

P2

● Mode Display

Each time you press the power / mode button, you will switch in the order of RF cleaning mode-RF moisturizing mode-RF EMS UP-RF LED-cool pack mode-RF cleaning mode (repeat).

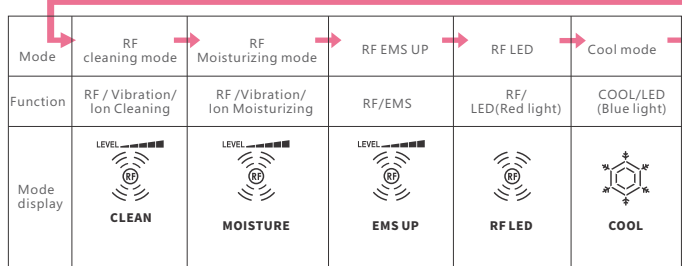

- The RF on the model display is shown or hidden when performing RF cleaning mode, RF moisturizing mode, or RF EMS mode care. When a short beep is heard, RF therapy has started (the mode display shows the "RF" word). When two "beep" sounds, it indicates that the ion or EMS has started (the RF display of the mode display disappears).

● Intensity display

After selecting RF clean mode or RF moisturizing mode, press the strength button to select the ion output strength. After selecting RF EMS UP mode, press the strength button to select EMS strength. Strength is divided into five levels, please choose the individual comfortable strength. Each time the strength button is pressed, the strength level will increase gradually



For the first use, it is recommended to start at Level 1 and increase intensity according to comfort levels.

P4

Modes Overview

Mode	Treatment Type	Display / Screen Function	Can be used directly on skin VIN	Suggested Treatment	Recommended Treatment Time	Frequency
RF Cleaning Mode CLEAN	Thermal therapy, positive ion attracts impurities from the skin	RF Ion Cleaning Vibration	NO Use with a cotton pad	To remove makeup	2-3 minutes per session	
RF Moisturizing Mode MOISTURE	Thermal therapy, negative ion to neutralize free radicals and boost absorption	RF Ion Moisturizing Vibration	NO Use with a cotton pad (A sheet mask can be used as is)	To apply skincare products, or sheet mask	2-3 minutes per session	2-3 sessions per week
RF EMS Mode: EMS UP	Thermal therapy, EMS stimulates muscle contractions to firm and lift	RF EMS	NO Use with a cotton pad	To apply skincare products (lotion, creams, serums)	2-3 minutes per session	
RF LED Mode: RF LED	Thermal therapy, LED phototherapy for fine lines reduction and skin radiance	RF LED	YES		2-3 minutes per session	
Cold Therapy Mode: COOL	Cold mode Shrink pores and depuffs skin	COOL	YES		2-3 minutes per session	Can be used daily

* Results may vary per individual.

On clean skin, apply your favorite skincare product so the device glides over skin (except in COOL mode). Do not use oils, scrubs with this device. Water-based skincare is preferred.

After using the COOL mode, it is normal for the device to feel warm. Please turn off the device after use of COOL mode and wait for the device to cool down.

Charging Mode >

This product operates with a built-in rechargeable battery and cannot be used while charging.

When using for the first time (or after a long time), make sure to fully charge the device before use. When charged, the lower left side of the screen will display the battery symbol and flicker. When the charge is complete, the battery symbol will remain static and will disappear when in use.

Insert the power adapter plug into the charging socket of the device and plug the power adapter into a power outlet.

Device can only be used with provided cable and power adapter.

After charging, please unplug the power adapter and the charging plug connected to the device.

The use time will vary according to the mode, intensity, and frequency of treatments.

Other Helpful Information >


- To install a cotton pad, tissue, or soft cloth (for CLEAN, MOISTURE, EMS UP Modes),
1. Remove the magnetic ring from the electrode surface.
 2. Place a piece of cotton on the electrode surface
 3. Place the magnetic ring back in and make sure the cotton pad, tissue, or soft cloth is securely in place. (Use a thinner type of cloth if the magnetic ring does not feel secure).
 4. Apply lotion or moisturizer on the cotton / cloth that is on the electrode surface. After use, wipe the remaining lotion and moisturizer on the electrode surface.

While using the device:

1. Make sure to constantly move the device on the skin. Doing otherwise may cause burns.
2. RF radiofrequency temperatures can feel different depending on the state of the skin.
3. Do not use on eyelid and periorcular care (except in COOL mode).
4. RF LED mode LED lights can shine on the eyelids and around the eyes. When using, close your eyes and the electrode surface should be more than 5 cm away from the skin.
5. Do not look directly at the LED lights in RF LED mode.
6. The ambient temperature is suggested to be 15°C ~ 30°C.
7. When using, please make sure that the electrode surface is in full contact with the skin.

After using the device:

Apply your favorite skincare regimen after your treatment to further nourish skin. Remove the cotton pad and wipe any residue of skincare products remaining on the electrode face or magnetic ring. This product is not waterproof, please wipe the residue on the electrode surface thoroughly after use with a dry cloth; do not wash with water!

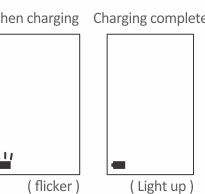


● How to Use Diagram

- <RF CLEAN, RF MOISTURE, RF LED modes> Please glide on the skin in the direction indicated by the arrow. RF LED mode LED lights can be directed toward the eyelids and around the eyes for care. when using, please close your eyes first, and the electrode surface should be above the skin 5cm

- Unusable position

- A position that can only be used under LED mode LED lights
- When using, please close your eyes first and the electrode surface should be at least 5cm away from the skin/eye area.

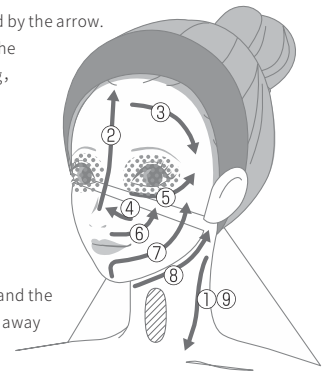

● <RF EMS UP mode>


Please glide on the skin in the direction indicated by the arrow

- Unusable position

- To slide in the direction indicated by the arrow

● <COOL Mode>

- Repeat [close-release-move] Steps to care on the face.
- Can be used on eyelids and around the eyes. <on eyelids and around the eye, please close your eyes first, and the electrode surface should be at least 5cm away from the skin /eye area.

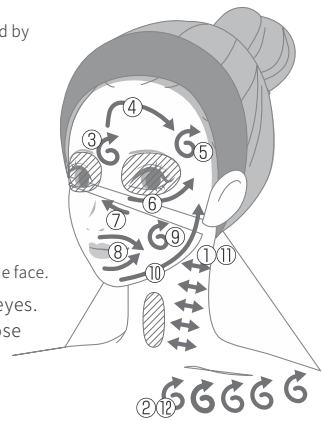

CAUTIONS ▶

Please read this manual carefully before use and use this product according to directions of use.

If there is an abnormal reaction, please stop using device immediately.

When using this product do not remain in the same area of the face for a long period of time.

Do not expose/approach to strong magnetic fields, do not directly expose to the sun, do not contact with corrosive liquid, or the device may work abnormally and get damage.

Do not put this device next to high temperature, humidity and large power appliances. This may result in fire, electric shock and other safety accidents.

This product is only for beauty care, please do not use it outside its specifications.

Please use the original charger to charge device, otherwise it may lead charging failure or burn out the device which may result in fire, electric shock and other safety accidents.

Do not disassemble, repair or modify this product by yourself, or it may cause short circuit, electric shock, or human injury.

This product is not waterproof, do not wash or immerse in water.

P7

SPECIFICATION ▶

Model Number: LL-VLIFT	Cold Model: ≤ 15°C
Vibration Frequency: 10000 rpm	Material: ABS/PC/Alloy
Product Size: 191*52*46mm	Charging voltage: 5V
Net Weight: 205g	Power: ≤8W
Charging Time: ≈3H	Country of origin: China

Warning:

1. Please stop using the product immediately when it is unusually hot or has an unusual odor, so as not to cause burns or fire.

2. Please conduct a skin patch test prior to first use. If you find any abnormal feeling/reaction during use, stop using it immediately. If abnormal skin occurs, stop using it and consult a doctor.

3. Use intensity 1 for initial use.

4. When using RF LED mode and RF COOL mode, do not look directly at LED lights, or you may damage your eyes or skin.

5. The following should not be using this device: Infectious patients, heart disease patients, high fever patients, malignant tumors, acute patients, medical appliance, skin disease patient, pregnant/lactating person, excessive sunburns, muscle particularly sensitive, vasodilator, Being hospitalized.

P10

WARRANTY INFORMATION :

Please use this as your reference for any customer service or warranty needs you may have when the device is purchased within the United States. When purchasing outside the United States, please contact your place of purchase for local warranty information. Our warranty is limited to the terms set out below:

What is Covered in the Warranty:
This lifetime warranty covers only defects in materials and workmanship in your V-LIFT Ionic Device when purchased through an authorized partner or directly from LUNA'MOUR. Your warranty is valid for you alone, the original purchaser, and cannot be transferred. To validate your warranty, please keep your original purchase receipt together with these warranty conditions for the duration of the warranty period. If the product was a gift, obtain a copy of the original receipt or have the original purchaser process the warranty claim on your behalf.

What is Not Covered by this Warranty?
This warranty does not cover: wear, tear, or deterioration, damage, failure caused by neglect, abuse or misuse, improper care, impact or insertion of foreign objects, use with a voltage connector, repair or alteration carried out by parties other than LUNA'MOUR or its authorized agents, or loss of device.

What Will LUNA'MOUR do if your device is defective?
If you discover a defect in your device covered by the warranty, contact Customer Care at hello@lunamour.co for assistance. If we confirm a defect covered by this warranty, LUNA'MOUR will at its sole discretion decide whether to replace your device with either a new or refurbished device. If we find that your device is not defective, as determined by LUNA'MOUR in its sole discretion, we will notify you and return the device to you. All warranty services must be authorized and performed by LUNA'MOUR or a service provider authorized by LUNA'MOUR.

WARRANTY LIMITATION & DISCLAIMERS
This warranty gives you specific legal rights and is in addition to your statutory rights. Your remedies for the breach of this limited warranty are limited solely to the repair and replacement of your device as set forth in this warranty statement.

To the maximum extent permitted by applicable law, LUNA'MOUR disclaims all implied warranties. In no event shall LUNA'MOUR be liable for any special, incidental, consequential, punitive, or indirect damages based upon breach of warranty, breach of contract, tort (including negligence or strict product liability), or otherwise, even if LUNA'MOUR is advised or aware that such damages may arise. Such damages include but are not limited to: loss of device, cost of any substitute equipment, the claims of third parties, or damage to property.

If the applicable law does not allow limits on warranties, or on the remedies for breach in certain transactions, the limits in this paragraph and the disclaimer or exceptions from warranty may not apply. This limitation does not apply to claims for death or personal injury caused by a defect in the device or any other liability that cannot be excluded or limited by law.

These warranty terms, conditions, and limitations allocate the risks of product failure between LUNA 'MOUR and you, the purchaser of a LUNA'MOUR branded device. This allocation is recognized by both parties and is reflected in the price of LUNA'MOUR branded device. By using the LUNA'MOUR branded device, you acknowledge that you have read these terms, conditions, and limitations, understand them, and are bound by them.

P11

Important Safety Information:

- If the product is faulty, please contact the customer service center, please do not disassemble or modify the machine by yourself. Do not operate device where it would come in contact with water. Do not use or charge device with power adapter that is damaged and/or faulty. Do not harshly bend or abuse the wire of USB charging cord. Do not attempt to take the device apart. Device is not serviceable by user. This will void any warranty.

Stop using the device immediately if you think it is defective. Please contact Client Care at: hello@lunamour.co

Product Representation :

LUNA'MOUR is constantly improving its products with the latest technology. As a result, the actual product may differ slightly from the product shown in the images on the packaging and in the user guide.

Skin Patch Test:

Prior to using the device for the first time, it is recommended that you conduct the following skin test on your forearm:

- 1) Keep your forearm dry and clean.
- 2) Turn on the device which will activate level one of the EMS mode, and conduct a three-minute treatment. Gently press device over your skin so it is in full contact. If rash, itch or irritation occurs on the tested skin, please stop the test immediately. If skin feels abnormal, also stop use immediately.
- 3) If the area is not red and no negative condition or fever appears after 24 hours, then you may keep using the device.
- 4) It is recommended that you conduct the foregoing skin test prior to using each of modes prior to usage.

5-Mode Facial V-LIFT Device

Client Care :
For information regarding device usage, warranty, and additional information, please contact Client Care at hello@lunamour.co

Disclaimer :
This product is not intended to diagnose, cure or prevent any disease. LUNA'MOUR makes no claims, representations, or warranties regarding the ability of this product to cure any mental, physical, or skin conditions. A qualified health professional should be consulted with regards to any condition requiring medical attention.

www.lunamour.co
hello@lunamour.co

Manufacture SKU / Model #: D919
D991

P13