

LUNA'MOUR

LUNALIGHT
Ultimate Radiance LED Light Therapy Device



USER MANUAL

LUNALIGHT from LUNA'MOUR - The Ultimate Radiance LED Light Therapy Device!
Illuminate your path to timeless beauty and holistic well-being with the remarkable LUNALIGHT.
Our cutting-edge handheld device is meticulously engineered to redefine your skincare and wellness regimen, combining 41 potent Light Emitting Diodes (LEDs) with innovative interchangeable heads.
With the power to seamlessly switch between Red, Yellow, Green, or Blue Light Therapies, paired with the soothing option of incorporating sonic vibration into each 5-minute treatment, the LUNALIGHT was created to help transform your self-care routine from head to toe!



Product Specifications:
Model: LL-PHOTO
Dimensions: 57X65mm
Net weight: 98g
Rated power: 2.5W
Power specification: input: AC100-240V 50/60Hz
Output: 5V---500mA
Built-in lithium battery: 3.7V---700mAh

Product List:



Charge:
● Built-in lithium battery which can stand for about 2-3 months after being fully charged.
● Full charge will take about 2.5 hours.
● When fully charged, a blue light indicator will illuminate and turn off once charging is complete.
● When the device is low on battery, a red light indicator will illuminate signaling that the device needs to be plugged into a power source using the provided Micro USB cable.
● The device cannot be used while it is charging.

Turning On/Off:
● Press and hold the "Power" button for 2 seconds to turn the device on and press once more for 2 seconds to turn the device off.
● Please note the device has an auto-shutoff feature which goes into effect after 5 minutes provided no additional buttons are selected.

- To turn on the Light Emitting Diode Attachment, press the middle button indicated by a small icon of a light bulb (See page 2 for Device Table)
- To turn on the Sonic Vibration, press the top button indicated by a small icon of 3 wavy lines (See page 2 for Device Table)
- The Sonic Vibration mode features 3 intensity levels which can be adjusted by a short press on the Sonic Vibration button located at the top of device and indicated by a small icon of 3 wavy lines.

To eject or replace the Light Emitting Diode Attachment, press the Light Emitting Diode Ejection Button located on the side of the device (See page 2 for Device Table)

To replace the Light Emitting Diode Attachment, align the 2 prongs and 3 pins and press to attach to the head of the device.

- Use only on clean, dry skin.
- Be sure to apply the provided protective goggles prior to starting a treatment as the light source is strong.


Sonic Vibration; Specs & Benefits:
Sonic Vibration Therapy assists as a highly effective modality that helps to exercise muscles, stimulate tissue, and enhance blood flow which, in turn, may help reduce muscle and tendon soreness. Additionally, gentle vibration helps to open pores allowing for deeper and better absorption of various Light Therapy wavelengths.

Light Therapy; Specs & Benefits:
Red Light Therapy may be effective in helping stimulate collagen production and reduce the appearance of fine lines and wrinkles. The LUNALIGHT features Red Light Therapy at 625nm(±5nm) wavelength which provides deeper penetrating properties that can help stimulate skin cells and improve cell vitality.

Yellow Light Therapy is best known for skin cell rejuvenation and its soothing anti-inflammatory benefits. Studies have shown that yellow light therapy can effectively reduce inflammation and help promote healing.
The LUNALIGHT features Yellow Light Therapy at 587nm(±5nm) wavelength which helps boost lymphatic flow, remove toxins and waste from the skin, and promotes recovery.

Green Light Therapy is best known for helping to correct uneven skin tone, minimize hyperpigmentation, discoloration, and redness. Studies have shown the efficacy of green light therapy in reducing melanin production and enhancing skin tone uniformity. The LUNALIGHT features Green Light Therapy at 525nm(±5nm) wavelength.

Blue Light Therapy harnesses antibacterial powers which have been shown to effectively reduce inflammatory acne, redness, and irritation.
The LUNALIGHT features Blue Light Therapy at 465nm(±5nm) wavelength which can help reduce the production of oil secretion and help alleviate skin discomfort caused by environmental stresses.

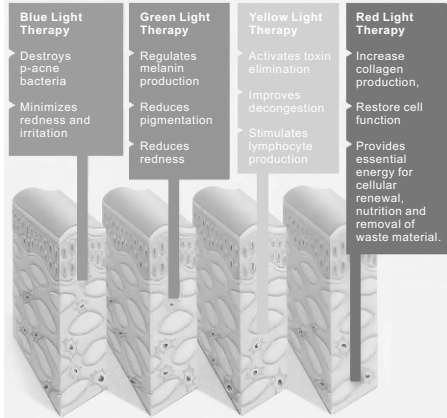

*Benefits and results may vary per individual.

Warnings:
CONTRAINDICATIONS:
Ask your doctor prior to using if you:
Are taking prescription or over the counter drugs known to cause light sensitivity.
Are, or suspect you are, pregnant.
Have any suspicious or cancerous lesions.

Please read the instructions carefully before use and use this product in the correct order.

Please keep out of reach of children. This device is designed for use on adults over 18 years of age and not recommended for use on children under the age of 12.

Please do not use if you are experiencing the following:
Heart disease, high fever, infectious diseases, malignant tumors, skin wounds, acute diseases, using medical appliances, skin diseases, pregnant women and/or breastfeeding, excessive sunburn, particularly sensitive skin, and those who are in the hospital.

Please consult your doctor prior to use if you are, or suspect you are, pregnant. Please consult your doctor prior to use if you have, or suspect you have, cancer. Do not use device on open wounds.
Do not use this device on skin that is painful to touch, raised (higher than the surrounding skin), crusting, inflamed, or bleeding. This may be a sign of a more serious condition and should be treated by your doctor.
No modification of this device is allowed.

- Do not place magnets, metal or conductive objects on the body of the device as it may cause malfunction.
- Do not rinse the instrument with water, use protective cloth provided to wipe clean after use.
- Do not wipe with alkaline detergent, alcohol, thinner, gasoline, nail polish remover, etc.
- If not in use for a long time, please make sure that the instrument is fully charged.
- Please stop using it immediately when there is abnormal heat, so as not to cause burns and fires and other safety issues.
- When the power cord is damaged, please do not use it to avoid short circuit and fire.

-DO NOT submerge the device into water as it will malfunction and/or cause injury.
-Do not place near flammable materials (spray, oil, natural gas, etc.)
-If product has abnormal noise, odor, or smoke, do not disassemble it at will, please contact the customer service email: hello@lunamour.co

Important Safety Information:

- If the product is faulty, please contact the customer service center, please do not disassemble or modify the machine by yourself.
Do not operate device where it would come in contact with water.
Do not use or charge device with power adapter that is damaged and/or faulty.
Do not harshly bend or abuse the wire of USB charging cord.
Do not attempt to take the device apart. Device is not serviceable by user. This will void any warranty.

Stop using the device immediately if you think it is defective.
Please contact Client Care at: hello@lunamour.co

Product Representation :
LUNA'MOUR is constantly improving its products with the latest technology. As a result, the actual product may differ slightly from the product shown in the images on the packaging and in the user guide.

Skin Patch Test:
Prior to using the device for the first time, it is recommended that you conduct the following skin test on your forearm:
1) Keep your forearm dry and clean.
2) Turn on the device.
3) Secure the protective glasses over your eyes and turn on the LED function & turn on the Sonic Vibration
4) Conduct a 2-minute treatment; gently press device over your skin so it's in full contact. If rash, itch, or irritation occurs on the tested skin, please stop the test immediately. If skin feels abnormal, also stop use immediately.
5) If the area is not red and there are no negative conditions after 24 hours, you may keep using the device.
6) It is recommended that you conduct the foregoing skin test prior to using each of the 4 Light Emitting Diode Attachments provided.

WARRANTY INFORMATION :

Please use this as your reference for any customer service or warranty needs you may have when the device is purchased within the United States. When purchasing outside the United States, please contact your place of purchase for local warranty information. Our warranty is limited to the terms set out below:

What is Covered in the Warranty:

This lifetime warranty covers only defects in materials and workmanship in your LUNALIGHT Device when purchased through an authorized partner or directly from LUNA'MOUR. Your warranty is valid for you alone, the original purchaser, and cannot be transferred. To validate your warranty, please keep your original purchase receipt together with these warranty conditions for the duration of the warranty period.
If the product was a gift, obtain a copy of the original receipt or have the original purchaser process the warranty claim on your behalf.

What is Not Covered by this Warranty?

This warranty does not cover: wear, tear, or deterioration, damage, failure caused by neglect, abuse or misuse, improper care, impact or insertion of foreign objects, use with a voltage connector, repair or alteration carried out by parties other than LUNA'MOUR or its authorized agents, or loss of device.

What Will LUNA'MOUR do if your device is defective?

If you discover a defect in your device covered by the warranty, contact Customer Care at hello@lunamour.co for assistance. If we confirm a defect covered by this warranty, LUNA'MOUR will at its sole discretion decide whether to replace your device with either a new or refurbished device. If we find that your device is not defective, as determined by LUNA'MOUR in its sole discretion, we will notify you and return the device to you. All warranty services must be authorized and performed by LUNA'MOUR or a service provider authorized by LUNA'MOUR .

WARRANTY LIMITATION & DISCLAIMERS

This warranty gives you specific legal rights and is in addition to your statutory rights. Your remedies for the breach of this limited warranty are limited solely to the repair and replacement of your device as set forth in this warranty statement.

To the maximum extent permitted by applicable law, LUNA'MOUR disclaims all implied warranties. In no event shall LUNA'MOUR be liable for any special, incidental, consequential, punitive, or indirect damages based upon breach of warranty, breach of contract, tort (including negligence or strict product liability), or otherwise, even if LUNA'MOUR is advised or aware that such damages may arise. Such damages include but are not limited to: loss of device, cost of any substitute equipment, the claims of third parties, or damage to property.


If the applicable law does not allow limits on warranties, or on the remedies for breach in certain transactions, the limits in this paragraph and the disclaimer or exceptions from warranty may not apply.

This limitation does not apply to claims for death or personal injury caused by a defect in the device or any other liability that cannot be excluded or limited by law.

These warranty terms, conditions, and limitations allocate the risks of product failure between LUNA'MOUR and you, the purchaser of a LUNA'MOUR branded d evice. This allocation is recognized by both parties and is reflected in the price of LUNA'MOUR branded device. By using the LUNA'MOUR branded device, you acknowledge that you have read these terms, conditions, and limitations, understand them, and are bound by them.

LUNALIGHT Ultimate Radiance LED Light Therapy Device

HOW TO USE VIDEOS:

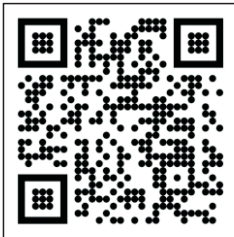

Client Care :

For information regarding device usage, warranty, and additional information, please contact Client Care at hello@lunamour.co

Disclaimer :

This product is not intended to diagnose, cure or prevent any disease. LUNA'MOUR makes no claims, representations, or warranties regarding the ability of this product to cure any mental, physical, or skin conditions. A qualified health professional should be consulted with regards to any condition requiring medical attention.

www.lunamour.co
hello@lunamour.co

Item Number: BRM-5700

