


# ULTRA Sculpt System

## 4-in-1 CV RF LED Face Body Device

### User Manual

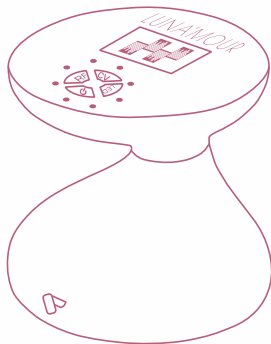

[www.lunamour.co](http://www.lunamour.co)  
[hello@lunamour.co](mailto:hello@lunamour.co)

Please read this user manual carefully before use as it includes important safety information.



Read all instructions prior to use.

**SAVE THESE INSTRUCTIONS**


# Table of Contents

•Meet Your Ultra Sculpt .....	2
•Device Features: Buttons .....	3
•Device Features: Intensities .....	4
•Featured Technologies .....	4
•Auto Power OFF function .....	5
•How to Charge .....	5
•Skin test before use .....	6
•How to Use .....	7
•Treatment Areas and Recommendations .....	8
•Cleaning Your Device .....	10
•Device maintenance .....	10
•Troubleshooting .....	10
•Helpful Tips for Best Results .....	12
•Product specification .....	13
•Client Care .....	13



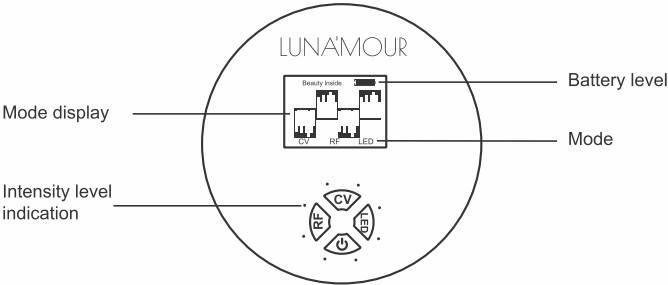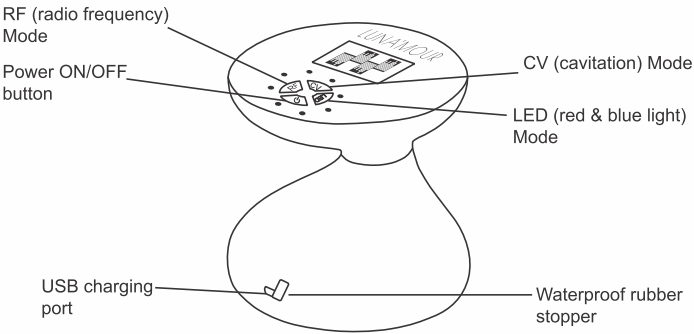

# Meet Your Ultra Sculpt

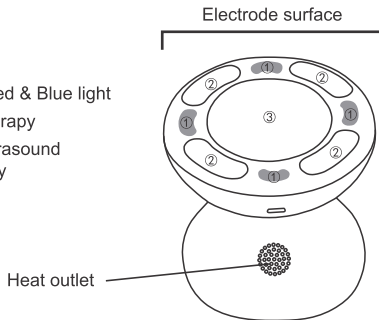

- ① Led Red & Blue light
- ② RF therapy
- ③ CV ultrasound therapy




## Device Features: Buttons


Button	Function	Indication lights
	Power on Power off	Blue light lights up Blue light lights off
	Two levels	RF mode blue light lights up
	Two levels	CV mode blue light lights up
	Press full mode on	LED mode light lights up

## Device Features: Intensity Levels

Two levels for RF therapy: Press once for level 1, press twice for level 2, 2 indication lights light up.

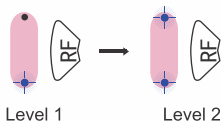

Two levels for CV ultrasound therapy: Press once for level 1, press twice for level 2, 2 indication lights light up.

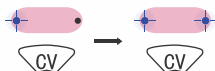

One level for LED red&blue light therapy. Press once 2 indication lights light up.



\*For each mode, higher level means a higher intensity.

## Featured Technologies

By combining ultrasound and RF Technology with red + blue light therapy, the Ultra Sculpt System helps reduce the appearance of cellulite as well as slim and sculpt various areas of the body such as your buttocks, thighs, waist, hips, abdomen and arms.

There are three main modes for this device, which LED therapy can be combined with the other two modes.

## Instruction for 3 Modes:

### **RF therapy mode:**

100Khz/ Min RF( radiofrequency) Therapy reaches deep into the skin layers creating a deep thermal effect that activates collagen production and decreases the volume of adipocytes(fat cells) , which in turn helps tighten skin, reduce cellulite and fine lines.

### **CV ultrasound mode:**

The 1Mhz ultrasound therapy which is often used by professional slimming machines. High frequency vibration caused by ultrasound therapy can remove cellulite deep under skin, reduce stubborn fat, achieve the slimming and body shaping effect.

### **LED red blue light therapy mode:**

Red light therapy and blue light therapy are often used together and complement one another. Together, red and blue light therapy can treat signs of aging, including age spots and wrinkles as well as sun damage.

## Auto power off function

The default time setting for each mode is 10 minutes. The device will enter Standby mode after 10 minutes.


Standby mode is 2 minutes. Device turns off automatically after 2 minutes.

## How to Charge


Charge the Ultra Sculpt System for at least 2 hours before you use it or anytime the battery indicator is "empty".

To charge :

1. Plug the USB cord into the Ultra socket and into a power outlet.
2. The battery indicator will show it is charging while the device is plugged
3. Once fully charged, unplug the charger

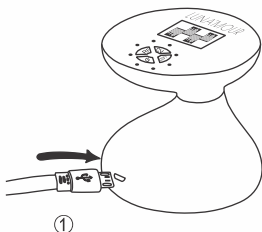

①

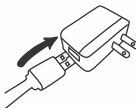

②



③




**Attention:**

- Do not charge the device in a damp environment.
- If charge indication light does not light up, take off the usb cable and charge again.
- Remove any gel, essence or dirt on the device before charging.

## Skin test before use

It is recommended to test if the device is suitable for your skin before first use.

- ① Keep your inside arm dry and clean.
- ② Apply gel or lotion on your inside arm, turn on RF mode, CV ultrasound mode and use for about 3 minutes. Gently glide the device over your skin. If no erythema or fever appear after 24 hours, then you can keep on use the device. Remember to use Level 1 during test.

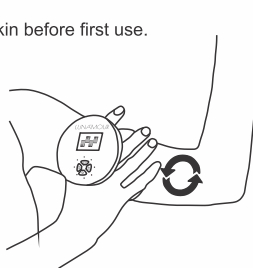






**Attention:**

- If rash, itch or irritation occurs on the tested skin, please stop test immediately.
- If skin feels abnormal, stop use immediately.

## How to use each mode

Apply gel, essence or lotion on your skin evenly before use. During the treatment, keep your skin moist for a successful treatment.

Button	Function	Treatment area	Frequency
	Reduce cellulite	body area need treatment	2-3 times each area, about 10 minutes for each area.
	Body slimming	body area need treatment	
	Body shaping & skin toning	Body or face area need treatment	

**Attention:**

Do not stay at same skin area for too long when using RF & CV mode to avoid over use. Glide the device in circles over your skin.

## Do not use on the below areas

Oral, private parts, mucosa etc



Neck, adam's apple, right below heart (near heart), finger, wrist, ankle etc area with thin skin.

Scratched skin, following plastic surgery, itch, fever, skin purulent inflammation areas.

RF & CV mode can not be used on face, head.

### How to contact your skin

Please make sure the device electrode surface fully contact your skin as following pic on the left.

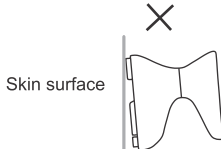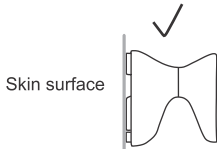

## Treatment Areas and Recommendations

Make sure device is in contact with skin at all times during session. Repeat below steps as needed.




: Parts for RF therapy



: Moving direction


### Arms( 5 minutes each)

- ① Hands face up, stretch out your arm.
- ② Glide the device from elbow towards armpit.
- ③ Same step for the other arm.

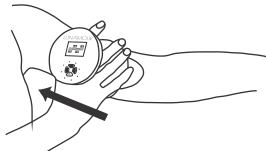



## Thighs (5 minutes each)

- ① Inside part of thigh, move from knee towards thigh root.
- ② Outside part of thigh, move from knee towards thigh root.
- ③ Same step for the other thigh.

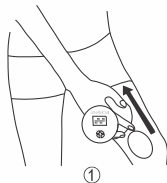



## Belly (Step ① for 5 minutes, step ② + ③ for 5 minutes.)

- ① Move from outside parts of belly toward inside parts in a Z route.  
Glide from Belly (waist) towards belly button slowly.
- ② Glide from belly towards thigh root.
- ③ Circle movement around belly button.

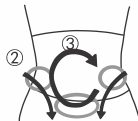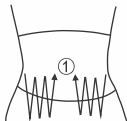

## Buttocks( Step ① for 5 minutes, Step ② + ③ 5 minutes)

- ① Glide from inner parts of buttocks toward outer parts in a Z route.
- ② Glide from inner parts toward outer parts.
- ③ Glide from top to bottom, from center towards thigh root.

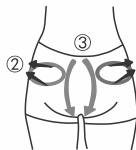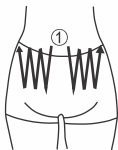

## Cleaning Your Device

Turn device off and unplug from any power source.

Wipe the tops of the electrodes with damp cloth to remove any residue from skin oils after each use.

Clean device body by using damp or microfiber cloth. Do not use chemical or alcohol based cleansers on exterior as harsh chemicals may discolor the metallic panels.

Allow device to dry completely prior to use.

## Device Maintenance

- Charge device only with the supplied USB power/ charging cord and a universal USB power adapter (not included)
- Turn off the device and disconnect from power supply when not in use.
- Do not place the device under high temperature or damp environment, or with direct sunlight.
- Keep it away from children.

## Troubleshooting

Please refer following chart for trouble shooting.

Issue		Possible causes		Solution
Unable to turn on the device	→	Do not press Power button	→	Press power button
	→	No battery	→	Charge first and then try
Can use only a few minutes	→	Low battery	→	Charge first and then use
	→	Default timer setting	→	power on the device again

Can not charge	→	USB charging cable does not insert in completely	→	Insert usb cable into charging port correctly
	→	Already fully charged	→	When fully charged, battery icon is full
Can not feel RF	→	Treatment area is not moist or too dry	→	Apply skin care products or gel before use
	→	Level too low	→	Adjust to higher level
	→	Move speed too fast	→	Recommend move at speed about 2-5cm/s
RF too hot	→	Electrode surface does not contact skin completely	→	Contact electrode surface completely with skin
	→	Stay at an area for too long	→	Move at a speed about 2-5 cm/s
	→	Electrode surface does not contact skin completely	→	Contact electrode surface completely with skin
	→	Intensity level too high	→	Adjust to lower level
Can not feel ultrasound vibration	→	Ultrasound vibration is not obvious, it's normal		
Device make noise during use	→	It's noise caused by ultrasound vibration Normal		

## Helpful Tips for Best Results

Problem	Solution
When is the best time to use?	It is recommend to use after bath and your skin is warm.
Is it suitable for sensitive skin?	Please refer to “skin test”, if no abnormalities occurred, then ok to use.
What type of lotions should I use?	Normal lotion, essence or gel is ok.

## Product specification

Power supply: DC 5V -1A USB

Measurement: 85.8(L)x85.8(W)x85mm(H)

N.W.: 270g. (Device only)

Ultrasonic frequency: 1Mhz

Power: Max 10W

Accessories: Main device, USB cable, Manual, Pouch

Charging time: Approximately 2H

## Client Care

You can reach Client Care by sending an e-mail to [hello@lunamour.co](mailto:hello@lunamour.co) for help with your device and warranty-related questions.

## Disclaimer:

This product is not intended to diagnose, cure or prevent any disease. LUNA'MOUR™ makes no claims, representations, or warranties regarding the ability of this product to cure any physical, skin, or mental condition. A qualified health professional should be consulted with regards to any condition requiring medical attention.

# Warnings:

## CONTRAINDICATIONS:

Ask your doctor prior to using if you:

Are taking prescription or over the counter drugs known to cause light sensitivity.

Are, or suspect you are, pregnant.

Have any suspicious or cancerous lesions.

Please read the instructions carefully before use and use this product in the correct order.

Please keep out of reach of children. This device is designed for use on adults over 18 years of age and not recommended for use on children under the age of 12.

Please do not use if you are experiencing the following:

Heart disease, high fever, infectious diseases, malignant tumors, skin wounds, acute diseases, using medical appliances, skin diseases, pregnant women and/or breastfeeding, excessive sunburn, particularly sensitive skin, and those who are in the hospital.

Please consult your doctor prior to use if you are, or suspect you are, pregnant.

Please consult your doctor prior to use if you have, or suspect you have, cancer.

Do not use device on open wounds.

Do not use this device on skin that is painful to touch, raised (higher than the surrounding skin), crusting, inflamed, or bleeding. This may be a sign of a more serious condition and should be treated by your doctor.

No modification of this device is allowed.

- Do not place magnets, metal or conductive objects on the body of the device as it may cause malfunction.
- Do not rinse the instrument with water, use protective cloth provided to wipe clean after use.
- Do not wipe with alkaline detergent, alcohol, thinner, gasoline, nail polish remover, etc.
- If not in use for a long time, please make sure that the instrument is fully charged.
- Please stop using it immediately when there is abnormal heat, so as not to cause burns and fires and other safety issues.
- When the power cord is damaged, please do not use it to avoid short circuit and fire
- DO NOT submerge the device into water as it will malfunction and/or cause injury.
- Do not place near flammable materials (spray, oil, natural gas, etc.)
- If product has abnormal noise, odor, or smoke, do not disassemble it at will, please contact the customer service email: [hello@lunamour.co](mailto:hello@lunamour.co)



## **Important Safety Information:**

- If the product is faulty, please contact the customer service center, please do not disassemble or modify the machine by yourself.

Do not operate device where it would come in contact with water.

Do not use or charge device with power adapter that is damaged and/or faulty.

Do not harshly bend or abuse the wire of USB charging cord.

Do not attempt to take the device apart. Device is not serviceable by user. This will void any warranty.

Stop using the device immediately if you think it is defective.

Please contact Client Care at: [hello@lunamour.co](mailto:hello@lunamour.co)

### **Product Representation :**

LUNA'MOUR is constantly improving its products with the latest technology. As a result, the actual product may differ slightly from the product shown in the images on the packaging and in the user guide.

### **Skin Patch Test:**

Prior to using the device for the first time, it is recommended that you conduct the following skin test on your forearm:

- 1) Keep your forearm dry and clean.
- 2) Turn on the device which will activate level one of the EMS mode, and conduct a three-minute treatment. Gently press device over your skin so it is in full contact. If rash, itch or irritation occurs on the tested skin, please stop the test immediately. If skin feels abnormal, also stop use immediately.
- 3) If the area is not red and no negative condition or fever appears after 24 hours, then you may keep using the device.
- 4) It is recommended that you conduct the foregoing skin test prior to using each of three modes prior to usage.

## WARRANTY INFORMATION :

Please use this as your reference for any customer service or warranty needs you may have when the device is purchased within the United States. When purchasing outside the United States, please contact your place of purchase for local warranty information. Our warranty is limited to the terms set out below:

### What is Covered in the Warranty:

This lifetime warranty covers only defects in materials and workmanship in your ULTRA Sculpt System 4-IN-1 CV RF LED Face Body Device when purchased through an authorized partner or directly from LUNA'MOUR. Your warranty is valid for you alone, the original purchaser, and cannot be transferred. To validate your warranty, please keep your original purchase receipt together with these warranty conditions for the duration of the warranty period. If the product was a gift, obtain a copy of the original receipt or have the original purchaser process the warranty claim on your behalf.

### What is Not Covered by this Warranty?

This warranty does not cover: wear, tear, or deterioration, damage, failure caused by neglect, abuse or misuse, improper care, impact or insertion of foreign objects, use with a voltage connector, repair or alteration carried out by parties other than LUNA'MOUR or its authorized agents, or loss of device.

### What Will LUNA'MOUR do if your device is defective?

If you discover a defect in your device covered by the warranty, contact Customer Care at [hello@lunamour.co](mailto:hello@lunamour.co) for assistance. If we confirm a defect covered by this warranty, LUNA'MOUR will at its sole discretion decide whether to replace your device with either a new or refurbished device. If we find that your device is not defective, as determined by LUNA'MOUR in its sole discretion, we will notify you and return the device to you. All warranty services must be authorized and performed by LUNA'MOUR or a service provider authorized by LUNA'MOUR.

### WARRANTY LIMITATION & DISCLAIMERS

This warranty gives you specific legal rights and is in addition to your statutory rights. Your remedies for the breach of this limited warranty are limited solely to the repair and replacement of your device as set forth in this warranty statement.


To the maximum extent permitted by applicable law, LUNA'MOUR disclaims all implied warranties. In no event shall LUNA'MOUR be liable for any special, incidental, consequential, punitive, or indirect damages based upon breach of warranty, breach of contract, tort (including negligence or strict product liability), or otherwise, even if LUNA'MOUR is advised or aware that such damages may arise. Such damages include but are not limited to: loss of device, cost of any substitute equipment, the claims of third parties, or damage to property.

If the applicable law does not allow limits on warranties, or on the remedies for breach in certain transactions, the limits in this paragraph and the disclaimer or exceptions from warranty may not apply. This limitation does not apply to claims for death or personal injury caused by a defect in the device or any other liability that cannot be excluded or limited by law.

These warranty terms, conditions, and limitations allocate the risks of product failure between LUNA 'MOUR and you, the purchaser of a LUNA'MOUR branded device. This allocation is recognized by both parties and is reflected in the price of LUNA'MOUR branded device. By using the LUNA'MOUR branded device, you acknowledge that you have read these terms, conditions, and limitations, understand them, and are bound by them.

# ULTRA Sculpt System 4-IN-1 CV RF LED Face Body

HOW TO USE VIDEOS:

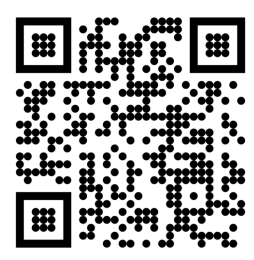

Client Care :

For information regarding device usage, warranty, and additional information, please contact Client Care at [hello@lunamour.co](mailto:hello@lunamour.co)

Disclaimer :

This product is not intended to diagnose, cure or prevent any disease. LUNA'MOUR makes no claims, representations, or warranties regarding the ability of this product to cure any mental, physical, or skin conditions. A qualified health professional should be consulted with regards to any condition requiring medical attention.

[www.lunamour.co](http://www.lunamour.co)  
[hello@lunamour.co](mailto:hello@lunamour.co)

Item Number :BP-1705  
LL-ULTRA

[www.lunamour.co](http://www.lunamour.co)

