




Panther Glow 3X LED Hair Rejuvenate Massage Comb PL045

Operating manual

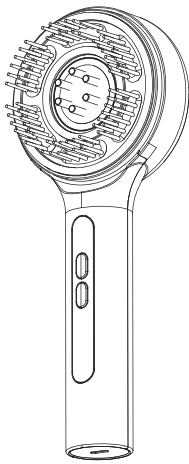

Introduction

The Prospera Panther glow 3X LED hair rejuvenate massage comb enhances scalp care with three-color light therapy—red to soothe, blue to control oil, and purple to boost vitality. Its roller ball liquid dispensing ensures even application of hair growth serum for better absorption, while the kneading massage function relaxes the scalp and promotes circulation. Designed for convenience, it features USB charging, a long-lasting battery, and an IPX6 waterproof massage head, making it easy to use and clean anytime, anywhere.

Product Specifications

• Three-color light therapy: red, blue, purple

• Adjustable 5-speed kneading

• Roller ball liquid dispensing

• Brush head waterproof

• Cordless

• 10 mins auto shut off

• Recharge time: about 2.5hrs

• Working time: ≥4hrs

• Item weight: 281g /0.62 lbs

• Main unit: 8.3” x 3.15” x 2.05”

• Box: 8.78” x 5.04”x 2.24”

• USB Charging

• Input: 5V 1A

• Rated Power: 3W

• Voltage: 3.7V

• Battery capacity: 1500mAh

• UPC 850073782002

Content

CONTENTS

Safety Warnings4

Product Illustration4

Operating the Product5

Consumer Warranty.....6

Warranty Card6

Benefits

Red, blue, and purple light therapy are used for different hair and scalp treatments based on their wavelengths and effects on hair follicles and scalp health. Here’s how each type of light works:

• **Red Light (630-670nm)**

• **Stimulates Hair Growth:** Red light penetrates deep into the scalp, promoting blood circulation and stimulating hair follicles.

• **Boosts Collagen & ATP Production:** Enhances cellular energy (ATP) production, which helps hair cells grow and repair.

• **Reduces Inflammation:** Helps in reducing scalp inflammation, which is linked to hair thinning and loss.

• **Prevents Hair Loss:** Strengthens hair follicles and prolongs the growth (anagen) phase of the hair cycle.

• **Blue Light (400-470nm)**

• **Targets Dandruff & Scalp Conditions:** Blue light has antimicrobial and antibacterial properties, making it effective against dandruff, seborrheic dermatitis, and fungal infections.

• **Reduces Scalp Oiliness:** Helps regulate excessive oil production, which can contribute to clogged follicles and hair thinning.

• **Calms Irritation:** Soothes inflammation and redness on the scalp.

• **Purple Light (Combination of Red & Blue)**

• **Balances Scalp Health:** Combines the benefits of red and blue light to promote hair growth while addressing scalp issues like dandruff and inflammation.

• **Enhances Hair Quality:** Improves hair strength, shine, and thickness.

• **Soothes Sensitive Scalp:** Works well for those with both hair loss and scalp irritation.

Best Uses for Hair & Scalp Treatments

• Red light is best for hair growth and preventing thinning.

• Blue light is ideal for oily scalps, dandruff, and bacterial issues.

• Purple light is great for overall scalp health and hair quality.

Safety Warnings

Consult your doctor before using this product if you have any lesions, skin conditions, heart issues, pacemakers, or other health-related concerns.

1. This product is intended for adult use only.


2. Do not let babies and children play with the product.

3. Do not use the product on an infant.

4. If you experience discomfort, stop using the product immediately.

5. Do not wash the unit using organic solvent such as thinner to prevent damage.

Product Illustration

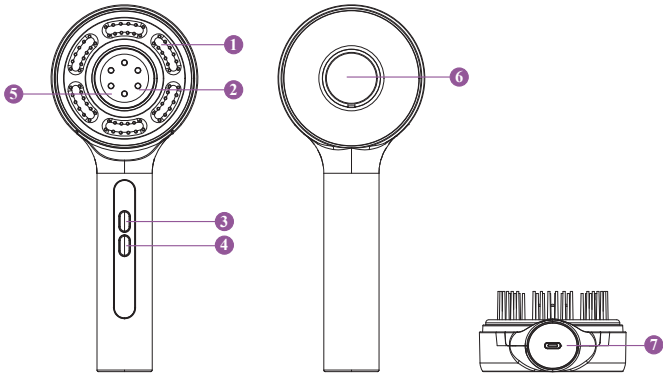

1 Kneading brush

2 Liquid rooler balls

3 Speed

4 ON-OFF-LED

5 LED ring

6 Liquid tank

7 Charging port

Operating the Product



Step 1: Preparation

Open the silicone cover of the liquid tank, add an appropriate amount of scalp care serum, and close the cover. Remove the transparent dust cover.



Step 2: Turn on The Power

Press and hold the power button to turn on the device. It defaults to red light and speed level 0. Short press the power button to adjust the LED light, and press the SPEED button to adjust the massage comb speed.



Step 3: Scalp Massage & Nourishing Care

Keep the liquid roller balls close to the scalp and move slowly for a soothing massage. The scalp care serum will gradually dispense, allowing you to enjoy both scalp massage and nourishing care at the same time.



Step 4: After Care

After the treatment, press and hold the power button to turn off the device. Clean any remaining serum from the liquid tank, wipe the massage head with a damp cloth, or rinse it under running water for 3–5 seconds. **Avoid getting the handle wet.**

Recommend usage: daily

To order a new health and beauty product, please visit www.prosperacorp.com.

Use coupon code: PROS\$5 for \$5 off any product above \$30.

Consumer Warranty

Prospera offers a 90-day standard warranty on massage products. Labor and parts are included during the warranty period. Please provide the model number and dated proof of purchase (sales receipt) when contacting Prospera for warranty service.

Non-warranty service is provided on a “per incident” basis. The consumer will verify that the product has failed and provide information for servicing the unit. All applicable repairs, parts, shipping, handling, local tax and “per incident” fees will be charged for non-warranty repairs and support calls.

Warranty Card



Warranty Registration Form

Product Name: _____

Product Model: _____ Purchase Price: _____

Date of Purchase: _____ Retailer: _____

First Name: _____ Last Name: _____

Address: _____

City: _____ State: _____ Zip: _____

Phone: _____ Email: _____

Signature and date: _____

Please mail to: _____ Tel:(925)292-5978

Prospera Corporation _____ Email: sales@prosperacorp.com

392 West Larch Rd Suite 31

Tracy, CA 95304-1644