



## Hot Compress EMS Deep Massage Belt SL011-A Operating manual



### Introduction

Experience wellness on the go with the Prospera EMS Belt—an elegant fusion of style and functionality. This cordless marvel utilizes medium frequency for a speedy, pain-relieving massage that enhances blood circulation, engages muscles, and boosts metabolism. Tailor your experience with eight modes—beat, knead, massage, press, shape, slap, exercise, and relax. With a remote controller and a visible LED screen, you can enjoy a more convenient user experience.

### Product Specifications

- Suits various sizes
- Comes with an extension
- With remote controller and LED screen
- Cordless
- 8 Massage modes
- 9 intensities options
- 3 temp options 40°C, 45°C, 50°C
- 30 mins auto shut off
- Item weight: 198g /0.44 lbs
- Main unit: 3.2" x 1.2"
- Remote controller: 4.2" x 1.85" x 0.7"
- Belt length: 38.6" x 4.5"
- Extension size: 13.1" x 2.3"
- Box: 8.7" x 4.7" x 2.3"
- USB Charging
- Input: 5V 1A
- Rated Power: 10W
- Voltage: 3.7V, 1200mAh (main device), 200mAh (remote controller)
- UPC 850011335833

### Content

#### CONTENTS


Safety Warnings.....	4
Product Illustration.....	4
Operating the Product.....	5
Consumer Warranty.....	6
Warranty Card.....	6

### Safety Warnings

Consult your doctor before using this product if you have any lesions, skin conditions, heart issues, pacemakers, or other health-related concerns.


1. This product is intended for adult use only.
2. Do not let babies and children play with the product.
3. Do not use the product on an infant.
4. If you experience discomfort, stop using the product immediately.
5. Do not wash the unit using organic solvent such as thinner to prevent damage.

### Product Illustration

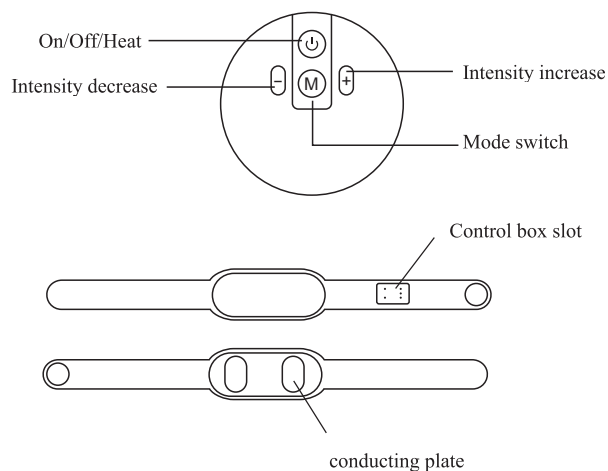

### Preparation before use

Connect the controller 3 and 2 buttons to the belt.

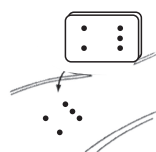

### Operating the Product

1. Press and hold the switch button for 2 seconds to start the machine.
2. Press the heating button briefly to adjust the hot compress temperature: first gear 40°C, second gear 45°C, third gear 50°C.
3. Short press the M mode button to select massage mode (8 modes in total).
4. Press the +/- key briefly to adjust the massage force (1-15 gear force). The display screen displays the current force.
5. The host can be disassembled, and it can work normally within 30 seconds after the function is set. If it is not installed after 30 seconds, the host is in standby mode.
6. After 30 minutes of operation, the massager will automatically shut down and enter into standby mode.

#### Step 1: Apply Gel

Apply cream or gel to the treatment area (abdomen, waist, arm, etc.) for better electric conduction.

#### Step 2: Install the Control Box

Attach the control box to the slot, ensuring it is in the correct attachment.

#### Step 3: Wear the Belt

Place the belt on your body with the conducting plate facing your skin. Adjust the tightness. Use the extension tape if needed for a proper fit.

#### Step 4: Massage Care

**Use either the main device control or remote controller.** Press and hold the switch/heat button on the main device to power it on. If using the remote controller, ensure it is turned on as well.

Short press the switch/heat button to turn the heating on/off (Temperature options: 40°C, 45°C, 50°C).

Short press the "M" button to choose your preferred mode (eight modes available: beat, knead, massage, press, shape, slap, exercise, and relax).

Use the +/- buttons to adjust the strength (9 intensity levels).

Long press the switch/heat button to turn the device and remote controller off and end the care session.

#### Step 5: Maintain the Device

Clean the belt with a wet wipe or cloth.

Recommend usage: daily

To order a new health and beauty product, please visit [www.prosperacorp.com](http://www.prosperacorp.com).

Use coupon code: PROSS5 for \$5 off any product above \$30.

### Consumer Warranty

Prospera offers a 90-day standard warranty on massage products. Labor and parts are included during the warranty period. Please provide the model number and dated proof of purchase (sales receipt) when contacting Prospera for warranty service.

Non-warranty service is provided on a "per incident" basis. The consumer will verify that the product has failed and provide information for servicing the unit. All applicable repairs, parts, shipping, handling, local tax and "per incident" fees will be charged for non-warranty repairs and support calls.

### Warranty Card



#### Warranty Registration Form

Product Name: \_\_\_\_\_  
 Product Model: \_\_\_\_\_ Purchase Price: \_\_\_\_\_  
 Date of Purchase: \_\_\_\_\_ Retailer: \_\_\_\_\_  
 First Name: \_\_\_\_\_ Last Name: \_\_\_\_\_  
 Address: \_\_\_\_\_  
 City: \_\_\_\_\_ State: \_\_\_\_\_ Zip: \_\_\_\_\_  
 Phone: \_\_\_\_\_ Email: \_\_\_\_\_  
 Signature and date: \_\_\_\_\_

Please mail to:  
 Prospera Corporation  
 392 West Larch Rd Suite 31  
 Tracy, CA 95304-1644

Tel: (925)292-5978  
 Email: [sales@prosperacorp.com](mailto:sales@prosperacorp.com)

