

140mm

100mm

100mm

NuvōMed®

ECM-6/1052 6 Gear Electric Cupping Massager



Instruction Manual Please read all instructions carefully and retain for future use

- Excessive fatigue, drunkenness, staying up late (2 hours before and after meals are also prohibited).
- Patients with skin allergy to metals and red light.
- The product has a heated surface, is not sensitive to heat, and children cannot express their wishes freely or manage themselves.

CAUTION (wrong operation may result in serious personal injury)

- Please use this product correctly according to the instructions.
- When using the product, the transparent can body should be close to the skin, with uniform contact, and avoid contact with wounds and scars.
- Long-term use of the same part of this product may cause low temperature burns or skin strains. Please use carefully.
- If you have any discomfort during or after use, please stop using it immediately and consult a doctor.
- Recommended for personal use. For cross-use, close and clean clear jar before use.
- Do not use chipped or broken transparent can.
- It is normal for redness, swelling or ecchymosis to appear after using this product, and it usually subsides naturally within 2 days (varies from person to person and skin type).
- When the product works for a long time, the surface will heat up, which is a normal phenomenon.
- If there is any abnormality, please remove it immediately.
- Do not use this product while driving.
- Do not plug and unplug the USB charging cable with wet hands to prevent electric shock or accident.
- This product contains lithium batteries. Please do not disassemble, replace, immerse in water or throw into fire to prevent product damage and explosion.
- After use, the blood vessels on the skin surface will dilate. You should rest for 45 minutes before bathing.

DESCRIPTION OF PARTS



HOW TO USE

- Before use**
1. Please read the safety precautions carefully before using this product. Please do not use this product if it does not meet the use conditions.
 2. It is recommended to clean and moisturize the skin with a wet towel to improve the sense of use experience.
 3. Tighten and install the transparent can body (clean the instrument with a wet towel before and after use).
 4. Check whether the product is complete and normal.
 5. Put the product with its transparent can body thoroughly cover your desire therapy position.

LIMITED WARRANTY

This manufacturer's product warranty extends to the original consumer purchaser of the product. Neither the retailer nor any other company involved in the sale or promotion of this product is a co-warranty of this manufacturer warranty.

WARRANTY DURATION: All materials and workmanship are warranted to the original consumer purchaser for a period of thirty (30) days from the original purchase date.

WARRANTY COVERAGE: This product is warranted against defective materials or workmanship. This warranty is void if the product has been damaged by accident, in shipment, unreasonable use, misuse, neglect, improper service, commercial use, repairs by unauthorized personnel or other causes not arising out of defects in materials or workmanship. This warranty is effective only if the product is purchased and operated in USA, and does not extend to any units which have been used in violation of written instructions furnished by manufacturer or to units which have been altered or modified or, to damaged products or parts thereof which have had the serial number removed, altered, defaced or rendered illegible.

WARRANTY DISCLAIMERS: This warranty is in lieu of all warranties expressed or implied and no representative or person is authorized to assume for manufacturer any other liability in connection with the sale of our products. There shall be no claims for defects or failure under any theory of tort, contractor commercial law including but not limited to, negligence, gross negligence, strict liability, breach of warranty and breach of contract. Under no circumstances will Manufacturer's / Distributor's maximum liability exceed the retail value of the product.

WARRANTY PERFORMANCE: During the above 30 days warranty period, a product with a defect will be either repaired or replaced with a reconditioned comparable model (at manufacturer's option). The repaired or replacement product will be in warranty for the balance of the 30 days warranty period and an additional one-month period. No charge will be applicable for such repair or replacement.

SERVICE AND REPAIR: If service is required for this product, you should first contact Nuvomed Inc. Customer Service at info@nuvomed.us or by calling 1 (877) 612-5619, Monday to Friday 10am to 6pm EST.

NOTE: Manufacturer cannot assume responsibility for loss or damage during incoming shipment. As a precautionary measure, carefully package the product for shipment, and insure it with the carrier. Be sure to enclose the following details with the product: your full name, return address and daytime phone number, a note describing the problem you experienced, a copy of your sales receipt or other proof of purchase to determine warranty status. C.O.D. shipments cannot be accepted.

GETTING STARTED

- Remove the product from the box.
- Please read all instructions and retain all packaging until you have checked and counted all of the parts.

PACKAGE CONTENTS

- Cupping Massager
- Type-C Charging Cable
- Instruction Manual

KEY FEATURES

- Helps relieve muscle soreness and stiffness
- Drudge the meridians
- Boost up metabolic process
- Promotes blood circulation
- 6 suction levels and 6 temperature levels for selection
- Built-in rechargeable battery

SAFETY PRECAUTIONS

DANGER (wrong operation can result in serious personal injury)

Do not use with the following medical electronic equipment at the same time, it will cause medical equipment error or damage, and seriously endanger life!

- Implantable medical electronic devices such as cardiac pacemakers;
- Life-sustaining medical electronic equipment such as artificial heart and lungs;
- Electrocardiograph and other medical electronic equipment;
- Patients with other medical implants.

Please do not use in the following situations to avoid accidents or discomfort:

- Pregnancy, menstruation, epidermal damage, varicose veins and taking drugs.
- Patients with infectious diseases, sensitive skin, postoperative patients, high fever or critically ill patients.
- Patients with malignant tumor or local injury and ulceration, traumatic fracture, acute suppurative inflammation, and skin loss of elasticity.
- Those with abnormal heart, cranial nerve and blood pressure.
- Serve nervousness, convulsions, spasms, and manic anxiety.
- Frequent bleeding diseases, such as thrombocytopenia, leukemia, hemophilia, allergic purpura and other bleeding diseases.
- Near the human heart, facial features, private parts or large blood vessels on the body surface, etc.
- Patients with skin perception disorders or skin abnormalities.
- Patients who are being treated by a doctor or feel abnormal.

Power on/off

- Press and hold the Power button for 1.5 seconds to turn on the unit, the suction and temperature indicator light will turn on. (The operation has a 20-minute countdown function, and will automatically shut down after 20 minutes.)
- Press and hold the Power button for 1.5 seconds will turn off the unit, the suction and temperature indicator light will turn off.
- The default mode is normal cupping mode, and the heating function is automatically turned on.

Suction control

- The default suction is level 3 when power is on.
- Short press the Suction strength button to adjust the suction level. Each short press will increase 1 level and the blue LED indicator lights will turn on to show the corresponding level.
- When the suction reaches the highest level, short press the button again will return the suction to level 1.
- It is recommended to start from level 1 for the first use.
- The suction level and the corresponding indicator lights (blue LEDs) are as follow:

Suction Strength (Default at 3)



Temperature (Default at 3)



Temperature control

- The default temperature is level 3 when power is on.
- Short press the Temperature control button to adjust the temperature level. Each short press will increase 1 level and the red LED indicator lights will turn on to show the corresponding level.
- When the temperature reaches the highest level, short press the button again will turn off the heating function, all temperature indicator lights will turn off, and the red light at the bottom of the unit will also go out.

Pause control

During cupping, when you feel the suction is too strong or the temperature is too hot, or you want to pause, there are two ways:


1. Short press the Pause button to pause the cupping and enter the standby state. Press the Pause button again or the Power button twice, the suction will start again.
2. Press and hold the Power button to power off the unit. All functions are turned off at that time.

Warm Tips

There are 6 levels of suction strength, it is recommended to start from level 1. And there are 6 temperature levels, it is recommended to use level 1.

Notes:

- When using a liquid such as massage oil, please sure to apply the massage oil evenly on the skin, there is no obvious liquid on the skin so that the liquid will not go inside the unit and cause a damage.
- The unit is designed with a filter device, which can filter out most of the particles or liquid impurities when using, protect the unit and maintain the service life of the product. Please check and replace the filter regularly according to the frequency of personal use. After replacing the filter, just restore the transparent can.

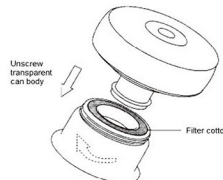

HOW TO RELEASE THE MASSAGER AFTER USING

After finishing the therapy, use a finger to lift the transparent can body at one suitable point on its circumference, so as to release the massager.

HOW TO CHARGE

- Before using it for the first time, charge it fully (about 3-4 hours).
- Connect the Type-C Charging Cable (included) to a suitable USB power outlet rated DC 5V/2A, and then the other end to the DC Charging Port of the device.
- Blue indicator lights will flash, when the battery is being charged.
- Blue indicator lights will turn off, when the battery is fully recharged.
- The battery can be recharged at any time and at any level.

CLEANING AND MAINTENANCE

- Please wipe the product body with dry non-woven fabric. Do not clean with solvents such as acid-base cleaners, alcohol, thinner, gasoline, and nail polish.
- After each use, please turn off the power before cleaning (wipe with a damp towel).
- Please store the product in a safe and dry place.
- Do not leave the product in a state where water remains.
- Do not wipe with a hard brush or stainless-steel ball to avoid damaging the surface.
- Do not directly inhale the instrument with water or other liquids for flushing.
- When the product is not used for a long time, please charge it regularly to ensure its normal operation.
- Do not use or store the product in the following location:
 - Places subject to vibration and shock, places exposed to direct sunlight, hot and humid places; places with water, strong electromagnetic fields; near fireworks; places within the reach of children; etc.

SOLUTION OF ORDINARY MALFUNCTIONS

- It is normal if the motor makes noises during operation.
- If the rated working time is over, the power will automatically go off.
- Do not use this product continuously for over 20 minutes, please allow a 30-minute break before resuming.

PRODUCT SPECIFICATIONS

Product Size: 3.6" (L) x 3.6" (W) x 3.1" (H)
Suction Strength: 60~90kPa
Rated Power: 5W
Charging Voltage: 5V 2A
Working Time: About 2 hours
Built-in Rechargeable Battery: 3.7V 1200mAh

FCC INSTRUCTIONS

Caution: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the product.

NOTE: This product has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This product generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If the product does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the product and receiver.
- Connect the product into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.