

SHARPER IMAGE®

POWERBOOST® PRO+ COMPACT

DUAL-TEMPERATURE
PERCUSSION MASSAGER

CONGRATULATIONS ON YOUR NEW POWERBOOST® PRO+ COMPACT.


Please read all instructions and warnings prior to use.

Keep this manual as it contains important information for future reference.

BC

WHAT YOU GET

- a. PowerBoost® Pro+ Compact
- b. Interchangeable attachments (x4)
- c. USB-C charging cable
- d. Carrying case

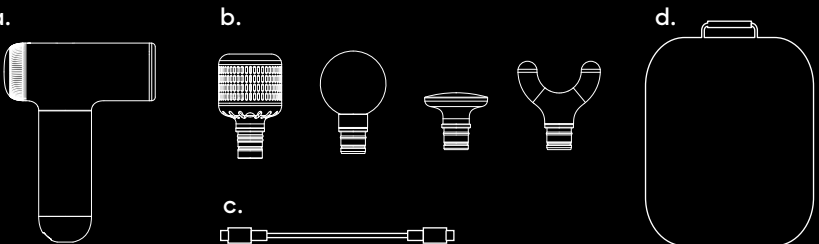


THE FUTURE OF RECOVERY

01

Plug the included USB-C charging cable into the input port on the bottom of the massager handle and then plug the other end of the cable into a USB-C power adapter (not included). Plug the power adapter into a wall outlet to begin charging. The lightning bolt on the display will flash red while the massager is charging and will turn off when charging is complete. The dots under the lightning bolt indicate the charge percentage.

Charging speeds may vary between power adapters. When plugged into a compatible charger, the massager fully charges in approximately 3 hours.

IMPORTANT: If you plug in the massager and the display does not light up, the power adapter may operate at an incompatible voltage. Unplug the massager immediately and use a different power adapter. Only charge the massager with compatible power adapters.

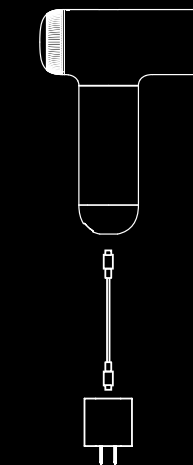


Power adapter not included

02

Select and plug the desired attachment into the front of the massager. Make sure the attachment is securely pressed in before beginning your massage.

Power off the massager before changing attachments. Twist and pull to remove attachments.


IMPORTANT: To use the Hot + Cold attachment, power off the massager and insert the Hot + Cold attachment into the front of the massager. Firmly press the attachment into place and then power on the massager. If the attachment is properly connected, the TEMP icon will light up. If the TEMP icon does not light up, turn off the massager and reconnect the attachment.

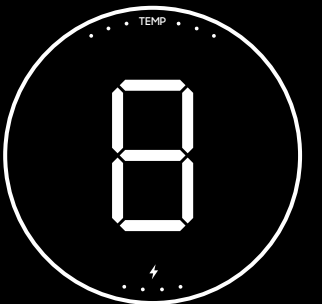

03

To power on the massager, press and hold the power button on the bottom of the massager handle. Press the power button again to toggle through the 6 massage intensity levels. Press and hold the power button to turn off the massager.

Note: The massager is equipped with advanced anti-stall technology. To avoid motor stalling, massage speed will automatically increase when extra pressure is applied.

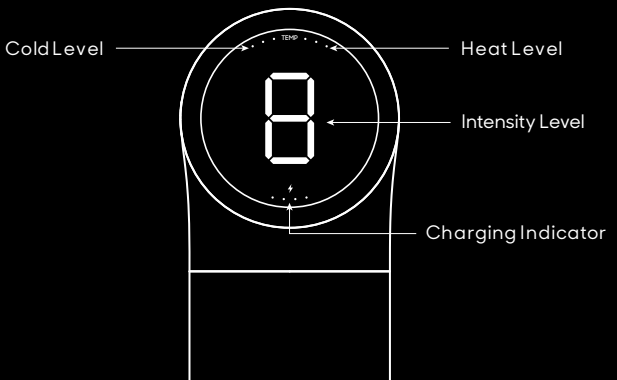

04

When you using the Hot + Cold attachment, twist the display bezel clockwise to cycle through the 3 heating levels or counterclockwise to cycle through the 3 cooling levels.

Illuminated red dots at the top of the display indicate heating level. Illuminated blue dots at the top of the display indicates cooling level. If no red or blue dots are lit, active heating and cooling have been disabled.

Note: Wait a few seconds at each heating or cooling level for the attachment to reach max temperature.


CAUTION: Never continue heat or cold therapy for more than 20 minutes at a time. If skin becomes bright pink or red, stop use immediately.

05

Use the massager for up to 60 seconds per body area.

Enjoy your massage for up to 10 minutes.


Do not use on neck, head or any hard or bony areas of the body.

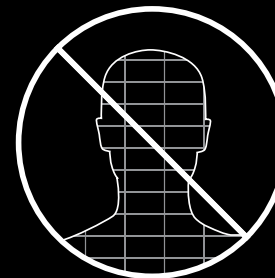

SAFETY WARNINGS: PLEASE READ BEFORE USE


- Not a children's item. Keep out of the reach of children.
- Always wear a shirt/clothing when using. Do not use on bare skin to avoid possible irritation.
- As with any massage, bruising may occur. Test your comfort level and check targeted areas often. Stop massage immediately if you feel any increase in physical discomfort or pain.
- Do not use attachments of any kind other than those that may be included.
- Do not fall asleep while using.
- Do not leave this item unattended when it is plugged in.
- Always allow item to cool if it becomes too hot.
- Do not expose to extreme temperatures or moisture.
- Keep away from sources of water and heated surfaces.
- Never touch the charging cable or ports with wet hands.
- Do not carry this item by the charging cable or attachments.
- For indoor use only.
- Thoroughly inspect products before every use. If the product appears worn, frayed/splintered, cracked, or broken in any manner, discontinue use and discard immediately.
- This product is not intended to treat, cure, or prevent medical conditions.
- Do not use this item if you are pregnant or have an injury or illness. Consult your physician if you have a medical condition that may prevent safe use of this item.

CHOOSE YOUR ATTACHMENT

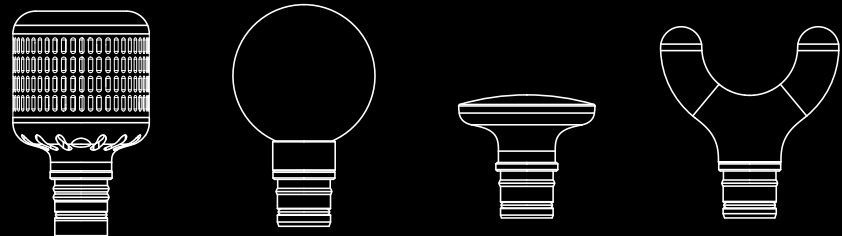


HOT + COLD


Heat or cool therapy for advanced relief.

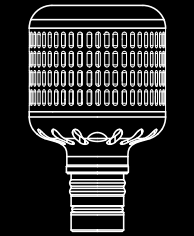

GOOD FOR
All Body Areas

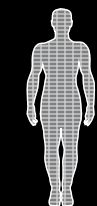


MUSCLE IMPACT
Low

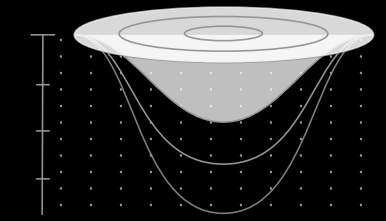


ROUND


Impact large areas for overall comfort.

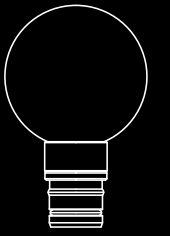

GOOD FOR
Glutes, Arms, Thighs, Back

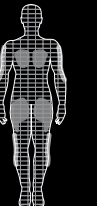

MUSCLE IMPACT
Low

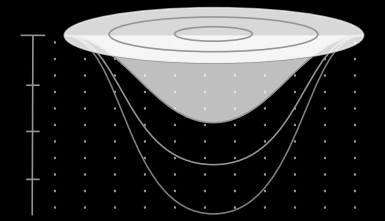


FLAT


Hammer-like impact for deeper recovery.



GOOD FOR
Chest, Shoulders, Back

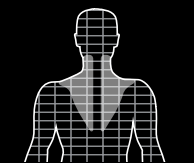


MUSCLE IMPACT
Medium

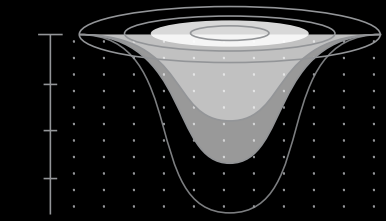


FORK


Dual impact to rock from side to side on muscle knots.

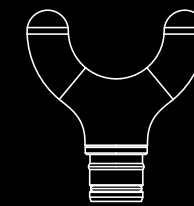

GOOD FOR
Shoulders, Biceps, Calves, Feet

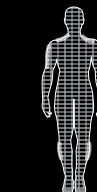

MUSCLE IMPACT
High

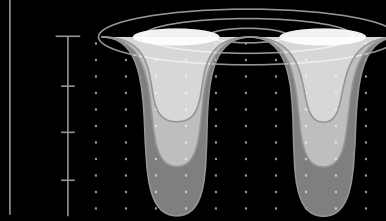

IMPORTANT: MAINTAINING BATTERY HEALTH

- Your massager is designed with a lithium-ion battery pack for cordless convenience, increased massage duration per charge, and frequent use and recharging.
- To maintain charge capacity and prevent battery damage, all lithium-ion batteries must be regularly used and recharged.
- For optimal battery health, keep the battery at approximately 50% state of charge when not in use.
- All batteries naturally lose charge when not in use. If you will not be using the massager for an extended period, check the massager at least every 30 days to ensure the battery is at approximately 50% state of charge. Recharge as needed to maintain a 50% state of charge but do not overcharge.
- **IMPORTANT:** Do NOT store the battery for an extended period with above 60% or below 40% state of charge. An incorrect state of charge during storage could irreparably damage the battery and may result in reduced charge capacity, inability to maintain charge, failure to recharge, or complete inoperability of the massager.
- For the health of the battery, always store the massager in a cool, dry place. The storage location should have a consistent temperature between 40°F and 70°F (5°C and 20°C) as well as a relative humidity at 70% or lower.

RECHARGEABLE BATTERY WARNINGS

- This product contains batteries that are non-replaceable.
- There is a risk of fire and personal injury if battery is punctured, damaged, or misused.
- Do not dismantle, crush, incinerate or short-circuit battery.
- Never expose battery to extreme temperatures or moisture.
- Keep away from flammable materials and direct sunlight.
- Only charge using the provided charger, cable, and/or adapters.
- Rechargeable battery is only to be charged by an adult.
- Adult supervision is required while charging the battery.
- Allow battery to cool to room temperature before charging.
- Charge on a non-flammable surface and keep away from flammable products.
- Discontinue charging if battery becomes hot, swells, smokes, or gives off a strong odor.
- Battery chargers used with the product are to be regularly examined for damage to the cord, plug, enclosure, and other parts. In the event of such damage, chargers must not be used until the damage has been repaired.
- Recycle or dispose of battery according to federal, state, and local laws.

FCC STATEMENT

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: 1. This device may not cause harmful interference, and 2. this device must accept any interference received, including interference that may cause undesired operation.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

WARNING: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

CARE INSTRUCTIONS

- Spot clean only. Unplug and wipe with a damp cloth.
- Air-dry only. Never place in a dryer.
- Allow to dry fully before charging.
- Do not submerge in water.
- Do not use detergent or soap when cleaning.

SHARPER IMAGE® NAME AND LOGOS (INCLUDING THE 77 LOGO, TOMORROW'S TOMORROW, AND POWERBOOST) ARE TRADEMARKS OR REGISTERED TRADEMARKS. DISTRIBUTED BY AND © 2023 MERCHSOURCE, LLC. IRVINE, CA 92618. ALL RIGHTS RESERVED.

FOR CUSTOMER SUPPORT, PLEASE CALL 1-800-374-2744.

PRINTED IN CHINA.