

88mm

150mm



FDA CLEARED K252142

**DL039 Penguin™ LED Sculpt Pro - Red/blue light, FDA Approved,
Whole body Microcurrent EMS System
Operating manual**



Introduction

This handheld EMS Multi-Functional Beauty Device combines Red LED light, microcurrent, EMS, heating, and vibration for professional skincare at home. With four modes—Beauty, Tightening, Skincare, and Refresh—it lifts, firms, and rejuvenates skin. Fully IPX6 waterproof, it's safe, convenient, and perfect for radiant, youthful skin anytime, anywhere.

88mm

150mm

Product Specifications

- FDA approved
- IPX 6 waterproof
- Working Modes: 4
- Strength: 3
- Heating temperature: 104±5.4°F
- Working voltage: 3.7V
- Charging: 5V/1A
- Rechargeable battery: 1,200mA
- Charging time: ≈2hrs
- Using time: ≈1hr
- Item weight: 0.34 lbs
- Main unit: 6.6" x 2" x 2.1"
- Box: 7.68" x 5.67" x 2.44"
- USB Charging
- UPC 850073782156

1

Content

CONTENTS

Safety Warnings.....	3
Product Illustration.....	3
Operating Instructions.....	4 - 8
Consumer Warranty.....	9
Warranty Card.....	9

2

88mm


Safety Warnings

Consult your physician before using this product if you have any skin lesions, skin conditions, heart conditions, a pacemaker, or any other medical concerns.

1. This product is intended for adult use only.
2. Keep out of reach of babies and children.
3. Do not use on infants.
4. If you experience pain or discomfort, discontinue use immediately.
5. Do not clean the unit with organic solvents (such as thinner), as this may damage the product.

Product Illustration

150mm

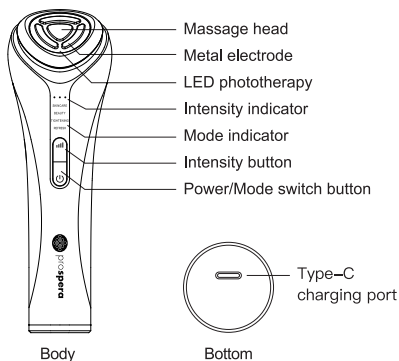

3


Operating Instructions

Getting Started






1. Skin Test

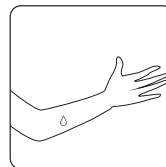For first-time use, please perform a patch test on the inside of your arm to determine whether your skin is suitable for this product.

1. Clean the test area on the inside of your arm.
2. Apply skincare products such as toner, essence, or gel evenly to the test area.
3. Turn on the device. Keep the electrode head in close contact with your skin and move it slowly in a circular motion (1–2 minutes per mode).
4. After the test, remove any residue from your skin and check for any abnormalities.
5. If you experience redness, itching, or a strong irritating sensation in the tested area, stop using the product immediately and rinse with clean water.
6. If you feel any discomfort after use, do not use the product again.
7. If no abnormalities (such as redness, swelling, or itching) occur within 24 hours, you may use the product with confidence.



2. Operation Instructions

- 1. Power On:** Press and hold the Power button  for about 2 seconds to turn on the device. You will hear two beeps. The device will automatically start in Skincare Mode.
- 2. Power Off:** While the device is on, press and hold the Power button  for about 2 seconds to turn it off. You will hear two beeps.
- 3. Mode Selection:** While the device is on, briefly press the Power button  to cycle through the four modes. The indicator light will remain on to show the selected mode.
- 4. Intensity Adjustment:** While the device is on, briefly press the Intensity button  to cycle through the intensity levels (Low–Medium–High). The indicator light will remain on to show the selected level. Changing modes during use does not affect the selected intensity level.
- 5. Vibration Control:** While the device is on, press and hold the Intensity button  for about 2 seconds to start or stop the vibration of the electrode head.
- 6. Timer Function:** The countdown timer starts automatically for each mode. The treatment time varies by mode. During use, you will hear a beep reminder to switch between the left and right sides of your face. When the time is up, the device will beep, the mode and intensity indicators will flash, and the device will automatically shut off after 1 minute of inactivity. The timer resets when you switch modes during use.

4

3. Function Descriptions


Modes	SKINCARE	BEAUTY	TIGHTENING	REFRESH
Functions	Blue LED light EMS (micro-current) Warming Frequency beating	Red LED light EMS (Penetration) Warming Frequency beating	Red LED light EMS (Tightening) Warming Vibration	Red LED light EMS (micro-current) Warming Vibration
Timer	3 minutes	5 minutes	4 minutes	5 minutes

Notes: 1) There are 3 intensity levels for each mode: low, medium, high.
2) The warming temperature: 104±5.4°F.

4. Usage guide

Reminder: Frequency and timing can be adjusted according to your own skin condition.

1) SKINCARE : Blue LED light + EMS (micro-current) + Warming+ Frequency beating

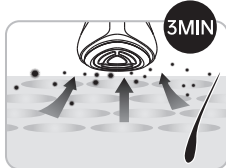


It deep cleanse pores, removes excess oil and dirt, make your skin refresh and healthy.

Recommend use with: facial cleanser/
makeup remover/cleansing water

1. Apply facial cleanser/makeup remover/cleansing water onto face evenly;
2. Select SKINCARE mode and keep electrode head in close contact with skin.
3. Hold the main unit and massage the face in circular motion from inside to outside.


* Frequency of use: 2-3 times a week, or 1-2 times a week for sensitive skin.

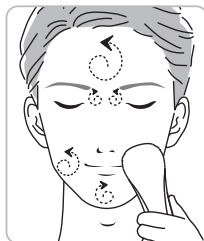

Note: Do not use the product on eyeballs or eyelids.

5

2) BEAUTY: Red LED light + EMS (Penetration) + Warming + Frequency beating

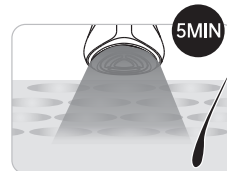

It improves blood circulation through warming and frequency beating, penetrate the nutrition of skincare product into deeper layer of skin through negative ion and promotes skin metabolism with acceleration on microcirculation through red LED light.

Recommend use with: gel/toner/essence or other skin care products


1. Apply skin care product onto face evenly;
2. Select BEAUTY mode and keep electrode head in close contact with skin;
3. Hold the main unit and massage the face in circular motion from inside to outside.

* Frequency of use: Daily or once in every 2 days. 2-3 times a week for sensitive skin.



Note: Do not use the product on eyeballs or eyelids

3) TIGHTENING: Red LED light + EMS (tightening) + Warming + Vibration

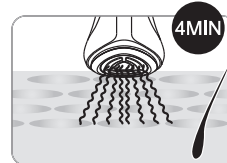

It refines the lines around face, tightens the muscle and skin tissue through EMS and warming functions; it promotes skin metabolism with acceleration on microcirculation through red LED light and warming.

Recommend use with: gel

1. Apply appropriate amount of gel to cheeks and forehead evenly;
2. Select TIGHTENING mode and keep electrode head in close contact with skin;
3. Hold the main unit and massage from center to outside and then pull up from bottom from bottom of the face.
4. If you feel the electrode head does not slide well, you can apply more gel.


* Frequency of use: 2-4 times per week, or 1-2 times per week for sensitive skin.



Note: Do not use the product on eyeballs, eyelids or laryngeal prominence.

6

4) REFRESH : Red LED light + EMS (micro-current) + Warming + Vibration

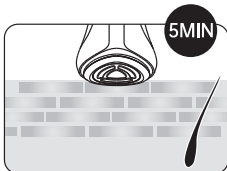


It can be used together with lotion or essence to accelerate absorption and keep skin smooth through EMS, warming and vibration functions.

Recommend use with: Lotion/essence

1. Apply appropriate amount of lotion/essence evenly to cheeks and forehead;
2. Select REFRESH mode keep electrode head in close contact with skin;
3. Hold the main unit and massage the face in circular motion from inside to outside.

* Frequency of use: 1-2 times per week.

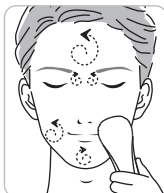


Note: Do not use on eyeballs, eyelids or laryngeal prominence.

5. Recommend Skincare Guide



Morning care


1. Cleanse the face;
2. Apply gel/toner/essence onto face evenly;
3. Operate in BEAUTY mode.

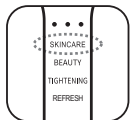

Evening care

1. Apply facial cleanser/makeup remover/cleansing water onto face evenly;
2. Cleanse the face in SKINCARE mode;
3. Operate with skin care products in BEAUTY mode.

6. Charging Instructions

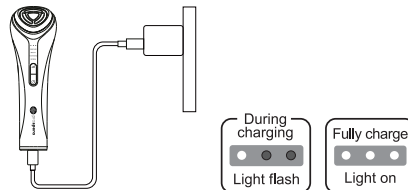

1. A power adapter is not included. Please charge the main unit using a DC 5V/1A power source.
2. Please charge the product at room temperature.
3. Use the provided Type-C cable to connect the device to the charging port located at the bottom of the main unit. The battery indicator will flash while charging. All three lights will remain on once the unit is fully charged.
4. After the unit is fully charged, unplug the power adapter and remove the USB charging cable.
5. The product cannot be operated while charging.

7. Cleaning, Maintenance and Storage

1. Turn off the product before cleaning. Do not clean the device while it is charging.
2. Do not clean the product with hot water, detergent, alcohol, gasoline, or other chemicals.
3. Do not immerse the product in water during cleaning.
4. Wipe the electrodes and main unit with a damp cloth, then dry them with a dry cloth or tissue.
5. Allow the product to air dry in a cool, well-ventilated place, then store it in the storage bag.
6. If the product will not be used for an extended period, charge it once every two months.
7. Do not place the product on uneven surfaces or high locations to prevent injury from falling or damage to the device.
8. Keep the product out of reach of children.

Recommend usage: daily

To order a new health and beauty product, please visit www.prosperacorp.com
Use coupon code: PROS\$5 for \$5 off any product above \$30.

