

NuDerma

High Frequency Hair & Scalp Wand

User Manual



WARNINGS AND PRECAUTIONS

Immediately discontinue use of NuDerma if you experience any rare adverse reactions such as an allergic response, swelling, skin discoloration, extreme redness, hyperpigmentation or if you feel any unusual pain during use. Stop using NuDerma if you experience dizziness, discomfort, or unusual sensations.

Do not use NuDerma if you use a pacemaker, neuro-stimulator, ECG monitor, or any other type of implanted electronic device that is inside your body. If you use a body-worn medical device, such as an insulin pump, consult your physician prior to use. If you wear permanent metal retainers or have braces, consult your dental professional prior to using NuDerma.

Do not use NuDerma if you are or might be pregnant. Consult your physician prior to use if you are breastfeeding.

If you have high or low blood pressure, a heart condition, epilepsy, or seizure disorders, consult your physician prior to use.

If you have any diagnosed skin conditions, consult a Dermatologist prior to using NuDerma. Do not use NuDerma on areas with active skin cancer or suspicious lesions. Do not use NuDerma on any areas of skin that lack sensation or have swelling, bruising, cuts, abrasions or infections.

The NuDerma Wand is not intended for persons with reduced physical, sensory or mental capacities unless supervised by a responsible adult party.

Do not use this device on children or animals. Keep the NuDerma and its attachments out of reach of children and pets.

Prior to use, we recommend removing any metal objects, such as jewelry, smart watches, earbuds or wired headphones, which may come in direct or indirect contact with NuDerma during a treatment session. Do not touch metal objects while NuDerma is in use.

While using NuDerma, do not touch a mobile phone or any electronic device that is plugged-in or charging.

The glass tube that is included with your NuDerma is fragile. Do not handle it by the tip. Only hold and handle it by the base when inserting and removing it from the NuDerma handle.

Do not handle or use your NuDerma with wet hands. Do not immerse any part of the NuDerma in water. Do not use NuDerma near bathtubs, showers, basins, or other vessels containing water.

Never use NuDerma near combustible liquids or gasses, or with oil-based or alcohol-based products.

Always store in a cool, dry place, away from extreme temperatures.

Only use official NuDerma glass tubes sold by Pure Daily Care. Only use the glass tube that is designed for use with the NuDerma wand model you are using: NuDerma Hair and Scalp Health. Using non-approved 3rd party glass tubes will void your warranty and could damage the device. Never power on the NuDerma before inserting the glass tube attachment. Never insert fingers or any objects other than the provided glass tube into the wand even if the device is unplugged.

Do not use any parts of your NuDerma that look visibly damaged. Unplug the device by grasping the plug and not the cable.

Do not plug the NuDerma wand into a power-strip, timer, or connect it via an extension cord. Do not plug NuDerma into an outlet where the electrical supply is unstable or likely to be turned off.

Always unplug the device safely from the socket if the cable becomes frayed or damaged or if the device starts to smoke, emits a loud noise, or appears visibly broken. If the power cord or adaptor is damaged do not attempt to repair, contact support@puredailycare.com

Do not attempt to dismantle or repair your NuDerma for self-service. If the appliance is not working as it should or has been dropped, damaged, left outdoors, or has come in contact with water, do not use and contact our support team at support@puredailycare.com

INTRODUCTION

Welcome to your NuDerma High Frequency Hair & Scalp Wand from Pure Daily Care.

This revolutionary tool is designed to energize the scalp, boost circulation, and support healthy hair vitality using advanced high-frequency technology—an esthetician-trusted modality adapted for at-home use.

Whether you are looking to improve scalp health, reduce buildup, or enhance the effects of a hair serum, you are in the right place. With over 1 million NuDerma users worldwide, we're excited to help you unlock the full potential of your hair and scalp.

Please read this user manual and any additional inserts carefully before use to ensure optimal results and safety.


WHAT IS HIGH FREQUENCY AND HOW DOES IT WORK?

NuDerma uses high frequency energy delivered through a specialized electrode (also called a tube) filled with neon-argon fusion gas, to stimulate the scalp and improve circulation. This advanced, non-invasive technology generates low-level thermal energy and oxygen-enriched molecules, both of which are believed to revitalize the scalp and help create a healthier environment for hair follicles.

Estheticians and dermatologists have relied on this modality for decades. Now, it is available for home use with an attachment specifically designed for the scalp: the comb tube.

When paired with the included scalp brush and Pure Daily Care Ultimate Hair & Scalp Serum, available separately, this system offers a powerful 3-step approach to support hair that looks fuller and healthier.

WHAT'S IN THE BOX

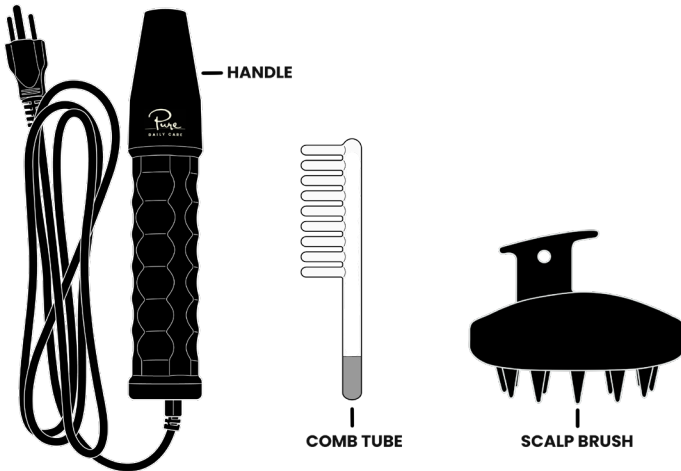


HOW TO INSERT AND REMOVE TUBES INTO HANDLE

1. Remove the NuDerma handle and comb tube from the box.
 - a. Always handle the tube by holding it just above the metal base.
 - b. Only insert or remove the tube while the NuDerma is powered off.
2. Insert the tube by gently pushing the metal base into the opening of the handle.
 - a. Never force or jam the tube – this can break the tube.
 - b. Never handle the tube by the tip.
3. Before removing the tube, ensure the handle is powered off.
 - a. Hold the tube as close to the handle as possible and pull it straight out – never twist it or pull it out at an angle.
 - b. Do not insert your fingers or any other foreign objects into the handle. Only use the included comb tube.

If the tube breaks or becomes stuck in the handle, contact our support team. Do not attempt to remove fragments yourself.

HOW TO USE THE NUDERMA

1. Begin with a clean scalp and clean, dry hair.
2. Insert the comb tube into the wand, then plug your NuDerma into a wall outlet. Do not use extension cords, plug adapters or power adapters.
3. Hold the wand and place one finger on the tube's base before turning the dial on. This helps prevent a harmless static sensation when making contact.
4. Turn the dial to the lowest power setting to start. You may hear a faint buzzing sound as the tube activates.

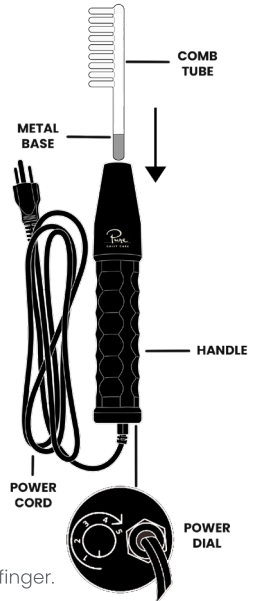

5. With your finger still on the tube, touch the comb attachment to your scalp, then release your finger.


6. Slowly glide the comb tube across your scalp using short, comb-like strokes, or lift and place the comb tube on the scalp and make small circular motions – whichever technique feels the most comfortable for your hair type or hairstyle.

- a. Avoid holding the wand still in one spot for too long.
- b. Never use on broken skin, irritated scalp, or open wounds.
- c. ensure the comb tube is in contact with the scalp using gentle pressure – do not press hard.

7. A full scalp session should last no more than 2 minutes, 2-3 times per week.

8. After your session, apply 1-2 droppers of the Ultimate Hair & Scalp Serum, or enough to create a thin layer across the entire scalp. Gently distribute it with your fingertips or the included scalp brush.

9. The scalp brush can be used daily—on wet or dry hair, with or without serum—for a soothing massage or to help distribute shampoo or treatments.

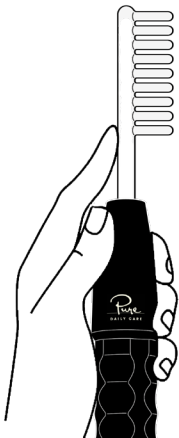

IMPORTANT SAFETY & USAGE GUIDELINES

- Always begin with the lowest intensity power dial setting for the first session. On subsequent sessions you may wish to increase the intensity to your level of comfort, but keep in mind that the system is effective even at the lowest setting. It is not necessary to use a high power level to achieve maximum benefits.
- Do not use the wand on wet hair or while in the shower.
- Do not use on skin conditions such as eczema, psoriasis, or other scalp sensitivities without consulting a healthcare provider.
- Avoid eye area, temples, and direct contact with facial skin.
- This product is intended for use on the scalp only.

HOW TO TAKE CARE OF THE NUDERMA

1. Power off and unplug your wand after each session. Remove the comb tube from the handle and store the handle and tube in their original box.
2. If needed, clean and sanitize the glass comb tube with a soft cloth and a small amount of rubbing alcohol. Dry with a soft, clean cloth. Ensure the comb tube is completely dry before using or storing it.
3. Place the handle and comb tube back into the foam inserts in its original box.
4. To prevent buildup, clean the scalp brush with a gentle cleanser and water. Dry it completely with a soft cloth.

TECHNICAL SPECIFICATIONS

Voltage	AC 110V 220V 230V 240C
Power	10W
Dimension	33 x 21 x 6 cm
G.W.	1 Kg

If there is any problem with your NuDerma device, please contact
support@puredailycare.com



DAILY CARE